

# GeneXpert MTBRIF and Mycobacterium tuberculosis Sputum Culture in Establishing the Diagnosis of Pulmonary Tuberculosis

*by* Berlian Beatrix Rarome

---

**Submission date:** 15-Mar-2021 01:23AM (UTC-0700)

**Submission ID:** 1533444491

**File name:** in\_Establishing\_the\_Diagnosis\_of\_Pulmonary\_Tuberculosis\_a....pdf (1.37M)

**Word count:** 5374

**Character count:** 28637

# Indonesian Journal of Tropical and Infectious Disease

Vol. 8 No. 3 September–December 2020

## Research Report

11

### GeneXpert MTB/RIF and *Mycobacterium tuberculosis* Sputum Culture in Establishing the Diagnosis of Pulmonary Tuberculosis and Rifampicin Resistance in Suspected Childhood Pulmonary Tuberculosis in Soetomo Hospital

Berlian Beatrix Rarome<sup>1</sup>, Nurul Aisah<sup>2</sup>, Retno Asih Setyoningrum<sup>3a</sup>, Ni Made Mertaniasih<sup>4</sup><sup>1,3</sup> Division of Respiriology, Department of Child Health, Faculty of Medicine, Universitas Airlangga / Dr. Soetomo Hospital, Surabaya<sup>2</sup> Medical School, Universitas Airlangga, Surabaya<sup>4</sup> Department of Microbiology, Faculty of Medicine, Universitas Airlangga / Dr. Soetomo Hospital, SurabayaReceived: 2<sup>nd</sup> October 2019; Revised: 8<sup>th</sup> September 2020; Accepted: 9<sup>th</sup> September 2020

#### ABSTRACT

The diagnosis of childhood tuberculosis remains a challenge worldwide. The GeneXpert MTB/RIF test, a rapid *Mycobacteria tuberculosis* diagnostic tool, was recommended for use in children. No pediatric studies of GeneXpert MTB/RIF assessing pulmonary tuberculosis within a hospital setting has been done in Indonesia. We evaluated the performance of the GeneXpert MTB/RIF test compared with sputum culture on Lowenstein-Jensen (LJ) for the diagnosis of childhood pulmonary tuberculosis. This study was conducted in pediatric respirology inpatient and outpatient Dr. Soetomo Hospital, a tertiary care facility in Surabaya between June and August 2015 with a cross-sectional design. We consecutively enrolled 27 children aged 3 months to 14 years who had history of close contact with adult tuberculosis patients and showed symptoms of pulmonary tuberculosis. Sputum collection was performed by induced sputum and three examination methods were performed (microscopic, GeneXpert MTB/RIF and sputum culture) simultaneously followed by a drug sensitivity test for specimens detected with MTB growth. The GeneXpert MTB/RIF test had a sensitivity of 100% (95% CI 100-100) and a specificity of 95% (95% CI 85-100). The positive predictive value for diagnosing pulmonary TB was 89% (95% CI 68-100), the negative predictive value was 100% (95% CI 100-100) and positive likelihood ratio was 20 (95% CI 2.82-128). The GeneXpert MTB/RIF test on one sputum sample rapidly and correctly identified all children with culture-confirmed pulmonary tuberculosis with high specificity. Similar results were obtained between GeneXpert MTB/RIF and sputum culture based on age groups and clinical manifestations. Rifampicin resistance were both detected in GeneXpert MTB/RIF and MTB sensitivity test.

**Keywords:** Childhood pulmonary tuberculosis; Sensitivity, Specificity; GeneXpert MTB/RIF

#### ABSTRAK

Menegakkan diagnosis tuberkulosis (TB) pada anak sampai saat ini masih sulit dikerjakan. GeneXpert MTB/RIF adalah suatu metode diagnostik baru yang dapat mengidentifikasi *Mycobacterium tuberculosis* (MTB) dengan cepat. Walaupun metode ini telah direkomendasikan pada anak-anak, namun penelitian tentang GeneXpert MTB/RIF dalam mendiagnosis TB paru anak di lingkungan Rumah Sakit (RS) belum pernah dikerjakan di Indonesia. Kami membandingkan hasil pemeriksaan GeneXpert MTB/RIF dengan kultur dahak MTB pada media Lowenstein Jensen (LJ) dalam menegakkan diagnosis TB paru pada anak yang diduga TB paru. Penelitian ini dilakukan di poli dan bangsal respirologi anak RSUD Dr. Soetomo antara Juni sampai Agustus 2015 secara cross sectional. Dengan sampling konsekutif mengumpulkan 27 anak usia 3 bulan sampai 14 tahun yang mempunyai kontak erat dengan penderita TB dewasa dan menunjukkan gejala TB paru. Pada setiap anak dilakukan pengambilan dahak dengan cara induksi dahak kemudian dilakukan tiga metode pemeriksaan sekaligus yaitu secara mikroskopis, GeneXpert MTB/RIF dan kultur yang dilanjutkan dengan uji kepekaan MTB bagi spesimen yang terdeteksi ada pertumbuhan MTB. Sensitivitas GeneXpert

\* Corresponding Author:  
retnooedijo@yahoo.co.id

MTB/RIF adalah 100% (95% CI 100-100) dan spesifisitas 95% (95% CI 85-100). Nilai duga positif GeneXpert MTB/RIF adalah 89% (95% CI 68-100), sedangkan nilai duga negatifnya adalah 100% (95% CI 100-100) dan Likelihood positive nya adalah 20 (95% CI 2,82-128). GeneXpert MTB/RIF mampu mendeteksi semua spesimen yang terdeteksi positif MTB oleh kultur dahak MTB namun dalam waktu yang lebih singkat dan dengan spesifisitas yang tinggi. Kesepadanan hasil antara GeneXpert MTB/RIF dan kultur dahak didapatkan berdasarkan kelompok umur dan manifestasi klinis TB. Selain Dalam mendeteksi resistensi Rifampicin, GeneXpert MTB/RIF memberikan hasil yang sama dengan uji kepekaan MTB.

**Kata kunci:** Tuberculosis paru anak, Sensitivitas, Spesifisitas, GeneXpert MTB/RIF

**How to Cite:** GeneXpert MTB/RIF and *Mycobacterium tuberculosis* Sputum Culture in Establishing the Diagnosis of Pulmonary Tuberculosis and Rifampicin Resistance in Suspected Childhood Pulmonary Tuberculosis in Soetomo Hospital. Rarome, BB, Aisah, N, Setyoningrum, RA, Mertaniasih, NM. *Indonesian Journal of Tropical and Infectious Disease*, 8(3), 152–162.

## INTRODUCTION

Difficulty to diagnose tuberculosis in children can lead to under and over diagnosis of TB which can cause higher morbidity and mortality. Confirming the diagnosis of childhood TB is a major challenge. However, research on childhood tuberculosis in relation to better diagnostics is often neglected because of technical difficulties, such as the slow growth in culture, the difficulty of obtaining specimens, and the diverse and relatively nonspecific clinical presentation of TB in this age group. While the classic presentation of childhood TB is prolonged cough and weight loss, HIV infection, with its chronic pulmonary manifestations and wasting, may confound the diagnosis of childhood TB. These difficulties are worsened by the increased incidence of multiple drug resistance.<sup>1,2</sup> Therefore, early diagnosis of TB in children is very important in order to control the incidence of TB disease.

Tuberculosis remains a major problem for the health of mankind. In 2012, the estimated incidence of TB cases were 8.6 million cases / year and 58% of these cases occur in Southeast Asia and the western Pacific. Approximately 50-60% of children living with adult pulmonary tuberculosis (PTB) patients who have positive acid-fast bacilli (AFB) sputum results will be infected with TB as well and about 10% of them will get TB disease.<sup>3</sup> World Health Organization (WHO) in 2012 estimated that there were 530,000 new cases of TB in children with a mortality rate of 74,000.<sup>4</sup> Indonesian TB data in 2012 showed

the proportion of TB cases in children among all TB cases was 8.2%.<sup>5</sup>

The diagnostic approaches that exist today are less sensitive. Although conventional examination with a microscope has a high positive predictive value for detecting *Mycobacterium tuberculosis* (MTB), its sensitivity is low. Examination using media culture with Lowenstein-Jensen (LJ) is still the gold standard for diagnosis but this test is difficult and requires a long time ( $\pm$  6-9 weeks) to get the results with positive results obtained only in 10% - 15% culture examination.<sup>6,7</sup> Polymerase chain reaction (PCR) test provides high sensitivity by multiplying deoxyribonucleic acid (DNA) of bacteria, and has been extensively evaluated in order to detect the DNA of MTB. GeneXpert MTB/RIF is an integrated and automated test with molecular approaches. Sample preparation, amplification and detection is done automatically by PCR. GeneXpert MTB/RIF is able to detect MTB as well as diagnosing resistance to Rifampicin. The results will be obtained in less than 2 hours.<sup>8,9,10,11</sup>

In December 2010, WHO has encouraged the use of GeneXpert MTB/RIF as a tool for the diagnosis of TB due to high sensitivity and specificity but studies on the use of GeneXpert MTB/RIF in children are still rare.<sup>11,12</sup> The aim of this study is to compare the GeneXpert MTB/RIF with MTB sputum culture examination in the diagnosis of PTB and rifampicin resistance in children with suspected PTB in Dr. Soetomo Hospital Surabaya.

## MATERIALS & METHODS

This study was analytical observational to compare the GeneXpert MTB/RIF assay with MTB sputum culture for detection of pulmonary tuberculosis and Rifampicin resistance in new pediatric inpatients and outpatients at the Department of Pediatric and Child Health, Dr. Soetomo Hospital, Surabaya, Indonesia, a tertiary referral center. This study was approved by the research ethics committee of Dr. Soetomo Hospital Surabaya. The parents of all study participants provided written informed consent.

Between June 2015 and August 2015, new outpatients and inpatients children, aged 3 months old - 14 years old, with the diagnosis of suspected tuberculosis were eligible for enrollment in the study. A patient with suspected tuberculosis was defined as having a symptom and risk factor screening (one or more of five factors: tuberculosis contact, cough for more than 3 weeks, weight loss, malnutrition, or fever for more than 2 weeks with unknown origin) according to the Indonesian National TB Program. Patients were excluded if they were deemed to have a poor prognosis, congenital heart defects, severe congenital abnormalities, acute hemodynamic disturbances (hypotension, shock, heart failure, decreased consciousness), critical illness (sepsis, renal failure, impaired liver function severe), have received TB treatment > 1 month and patients with HIV infection or if parents or guardians refused the informed consent.

Procedure of sample collection started with history taking and physical examination taken when administered. The clinical manifestations observed in this study were in accordance with TB scores commonly used in Indonesia, including a history of close contact with adult TB patients, the results of TST, fever  $\geq 2$  weeks are not unexplained, coughing  $\geq 3$  weeks, enlarged neck lymph nodes, inguinal and axillary, nutritional status, swelling of bones/joints and chest X-rays. Sputum samples were collected from children who could expectorate. In children who could not spontaneously expectorate, sputums were collected by induced sputum procedure. Smear microscopy, culture with LJ media and

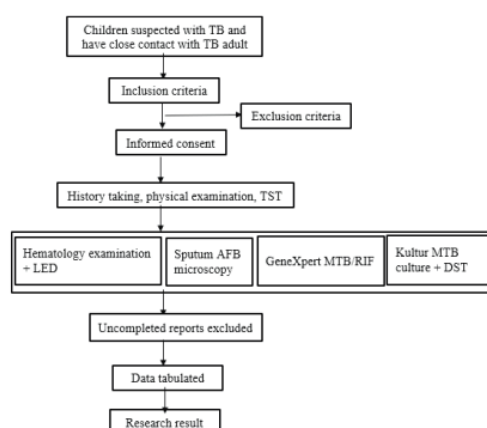


Figure 1. Operational Framework

GeneXpert MTB/RIF were done simultaneously on all samples as described previously. Cultures were classified as negative when no growth were detected after 8 weeks of incubation. Contaminated samples were retreated and re-cultured, and excluded if still contaminated. Drug susceptibility testing was done on LJ media. Sputum was added to the GeneXpert MTB/RIF sample reagent in a 1:1 ratio (1 mL of sputum to 1 mL of the sample reagent). Two mL of this mixture was added to the GeneXpert MTB/RIF cartridge and run in the machine in accordance with manufacturer's instructions.

All clinical and laboratory data were compiled in databases. Selected variables were exported to SPSS (version 21) for analysis. Comparisons of GeneXpert MTB/RIF assay and sputum culture assay were done with Pearson  $\chi^2$  or Fishers exact test. The sensitivity, specificity, and predictive values of the assays with 95% CIs were calculated. The equivalence between GeneXpert MTB/RIF assay and sputum culture were analyzed with McNemar test and Kappa. All statistical tests were two-sided with alpha of 5%.

## RESULTS

Twenty seven children were recruited and had sputum for analysis. Table 1 shows about half (51.9 %) of the children enrolled were younger than 5 years. The most common clinical

**Table 1.** Baseline Characteristics

Characteristics	N(%)
Age < 5 years old	14 (51.9)
≥ 5 years old	13 (48.1)
Gender Male	14 (51.9)
Female	13 (48.1)
Contact Identified	18 (66.7)
Not Identified	9 (33.3)
Contact MDR	5 (18.5)
Non MDR	12 (44.4)
Scar BCG Present	21 (77.8)
None	6 (22.2)
Cough > 3 weeks	18 (66.7)
Fever > 2 weeks	19 (70.4)
Lymph node enlargement	9 (33.3)
Nutritional Status Normal	12 (44.4)
Poor	10 (37.0)
Malnutrition	5 (18.5)
TST Positive	18 (66.7)
Negative	9 (33.3)
Bone destruction	1 (3.7)
Chest X-ray Suggestive TB	16 (59.3)
Not Suggestive TB	11 (40.7)
TB Score ≥ 6	23 (85.1)
< 6	4 (14.8)

manifestation reported was fever with no apparent cause lasting more than 2 weeks. Eighteen subjects (66.7%) showed positive TST and 9 (33.3%) had negative result.

Molecular examination with GeneXpert MTB/RIF provides the highest positive results of 33.3%, followed by microbiological culture (29.7%) and microscopic examination (22.2%) (Table 2). There was no false negative results for GeneXpert MTB/RIF examination.

GeneXpert MTB/RIF correctly identified all 9 rifampicin-sensitive on specimen analysis (Table 3). All 8 LJ culture positive specimens were also analyzed with the LJ drug-sensitivity test and none were Rifampicin resistance. Therefore MTB drug sensitivity cannot be analyzed statistically since the LJ results on solid media were equivalent with the results by the GeneXpert MTB/RIF (Table 3).

Table 4 shows that only 5 (27.8%) of 18 patients with positive TST results were confirmed

**Table 2.** Sputum Results

Sputum Examination	N (%)
MTB culture Positive	8 (29.7)
Negative	19 (70.3)
GeneXpert MTB/RIF Positive	9 (33.3)
Negative	18 (66.7)
Smear microscopic Positive	6 (22.2)
Negative	21 (77.8)

**Table 3.** Drug Sensitivity Test of Positive Results (8 LJ Culture and 9 GeneXpert MTB/RIF)

Rifampicin Sensitivity Test	N (%)
Drug sensitivity test with LJ	
Rifampicin Sensitive (+)	8 (100)
Rifampicin Resistance	0 (0)
INH Sensitive	7 (87.5)
INH Resistance	1 (12.5)
Etambutol Sensitive	8 (100)
Etambutol Resistance	0 (0)
Streptomycin Sensitive	7 (87.5)
Streptomycin Resistance	1 (12.5)
GeneXpert MTB/RIF	
Rifampicin Sensitive (+)	9 (100)
Rifampicin Resistance	0 (0)

**Table 4.** TST Result VS GeneXpert MTB/RIF

		GeneXpert MTB/RIF		Total
		Positive (%)	Negative (%)	
TST	Positive	5 (27,8)	13 (72,2)	18
	Negative	4 (44,4)	5 (55,5)	9
Total		9	18	27

positive by GeneXpert MTB/RIF, whereas 13 (72.2%) others showed negative results by GeneXpert MTB/RIF. There were 4 (44.4%) of 9 children who showed negative TST result, but confirmed positive on GeneXpert MTB/RIF, and 5 (55.5%) of 9 children showed a negative result on both TST and GeneXpert MTB/RIF.

Table 5 shows the equivalence results of GeneXpert MTB/RIF and MTB sputum culture in the age group ≥ 5 years old as many as 12 samples. There is only one sample which showed a positive result on GeneXpert MTB/RIF but confirmed

**Table 5.** GeneXpert MTB/RIF vs MTB Sputum Culture in Age Group

GeneXpert MTB/RIF		MTB culture		Agreement (%)	McNemar	Kappa
		(+)	(-)			
< 5 years old	(+)	2	0	100	p=1,000	1,000
	(-)	0	12			
≥ 5 years old	(+)	6	1	92,4	p=1,000	0,847
	(-)	0	6			

**Table 6.** GeneXpert MTB/RIF vs MTB Sputum Culture in Clinical TB Group

GeneXpert MTB/RIF		MTB culture		Agreement (%)	McNemar	Kappa
		(+)	(-)			
Clinical TB	(+)	8	1	94,2	1,00	0,883
	(-)	0	8			
Non Clinical TB	(+)	0	0	100	-	-
	(-)	0	10			

negative on MTB sputum culture. McNemar test shows no significant difference. Kappa test results shows significant reliability between the results of the GeneXpert MTB/RIF with MTB sputum culture MTB in the age group  $\geq 5$  years. The agreement between GeneXpert MTB/RIF and MTB sputum culture examination was 92.4%.

Table 6 shows the total positive clinical TB assessment in 17 (63.0%) of 27 children. There were 9 (53.0%) clinical TB children who have positive results of GeneXpert MTB/RIF. All children with clinically negative TB showed negative result on GeneXpert MTB/RIF. McNemar test showed no significant difference between the results of GeneXpert MTB/RIF and MTB sputum culture with positive clinical manifestations of TB. Kappa test showed significant equivalence between GeneXpert MTB/RIF and MTB sputum culture with positive clinically manifestation of TB.

This study showed that GeneXpert MTB/RIF had 100 % sensitivity (95% CI 100-100), specificity of 95% (95% CI 85-100), PPV 89% (95%CI 68-100), NPV 100% (95% CI 100-100), LR +20 (95% CI 2.82-128), LR -0.

## DISCUSSION

In our study, more than 50% of samples had a history of close contact to adult patients with positive AFB smear, cough > 3 weeks, fever of unknown origin > 2 weeks, positive TST results and X-rays which showed PTB process. More than 50% of the sample had a value of TB score > 6. Positive TST results were not always found in children suspected of pulmonary TB. In our

study, Positive TST results obtained in 66.67% of the samples, of which only 55.5% of the samples were GeneXpert MTB/RIF positive. Sekadde et al (2013) and Nicol et al (2011) obtained positive results only at  $\pm 30\%$  of the samples<sup>13,14</sup>, while Nataprawira et al (2001) get positive TST results only in 9.7% of children who had close contacts with adult TB patients or adults suspected of having TB in Bandung.<sup>15</sup> There are several factors that influence the results of TST, such as malnutrition which effect on phagocytosis, cellular immunity and cytokine production.<sup>16</sup> Malnutrition leads to lymphoid tissue atrophy, thus affecting the development, differentiation and cause a decrease in lymphocytes. Moderate and severe malnutrition lead to decreasing delayed-type hypersensitivity reactions and recall process.<sup>17</sup>

The positive results of microscopic and molecular examination in this study is quite high when compared to previous studies. Giang et al reported positive results of GeneXpert MTB/RIF on 8.6% of samples,<sup>18</sup> Nicol et al reported 12.8%,<sup>14</sup> Sekadde et al reported 14%,<sup>13</sup> Singh et al and Nhu et al reported a respective 16.9% and 16.2%.<sup>19,20</sup> This condition is likely due to several factors, such as the number of MTB in children (paucibacillary) and sputum production capabilities that are lacking in children. Different inclusion criteria with previous studies may also cause these differences. In our study, most of the sample had more than 4 clinical manifestations as well, while these other studies established inclusion criteria of children aged  $\leq 14$  years with at least 2 clinical manifestations of cough  $\geq 2$  weeks and one of the symptoms of weight

loss or fever  $\geq 2$  weeks with unknown origin, or a history of contact with adult TB patients, or a positive result on TST or positive X-ray for TB process.<sup>13,14,20</sup>

Microscopic examination (Olympus CH-20, Olympus Corp., Japan; 1000x magnification) is able to detect 66.7% of specimens with positive GeneXpert MTB/RIF and 75% of specimens positive by MTB sputum culture. Positive smear cases were found mainly in the  $> 5$  years old age group. This occurs because children  $> 5$  years old or adolescents have pathological features of "adult-type" TB that is not paucibacillary with more bacilli accumulated and generally give more positive results on microscopic examination.<sup>21,22</sup> Previous studies reported the same result. Marlowe et al in the US collected 217 sputum specimens and showed that microscopic examination is able to detect 73% of GeneXpert MTB/RIF positive results. Lawn et al (2011) only reported 45% positive smear result of TB cases. This can be explained due to the colony of microscopic detection capabilities is less sensitive than the other two examination modalities. The detection capability of the colony smear microscopy is  $5 \times 10^3$  to  $5 \times 10^4$  bacilli/ml, the detection capability of GeneXpert MTB RIF is 102-107 CFU/ml and the culture detection capability is 10-100 CFU/ml.<sup>12,23</sup> However, microscopic examination is not specific to diagnose TB because there are several other bacteria that are resistant to acid staining which are *Rhodococcus* spp, *Nocardia* spp, *Legionella micdadei*, cysts and isospores of *Cryptosporidium* spp, that will give a false-positive smear result.<sup>24</sup>

In this study, the sensitivity of GeneXpert MTB/RIF is 100% and specificity was 94.7%. There is one GeneXpert MTB/RIF positive result that is not detected by MTB sputum culture. PCR concept used in GeneXpert MTB/RIF sequence all of MTB DNA without the capability to detect the viability of the MTB. GeneXpert MTB/RIF false-positive may result from patients who had been treated as well.

Systematic reviews conducted by WHO in 2013 among 13 studies involving 2,603 participants mention pooled sensitivity of GeneXpert MTB/RIF TB was 66% (95% CI 52-77) and pooled

specificity was 98%.<sup>25</sup> Meta-analysis conducted in 2012 of 18 studies involving 10,224 specimens reported sensitivity of GeneXpert MTB/RIF amounted to 90.4% (95% CI 89.2 to 91.4) and specificity of 98.4% (95% CI 98-98, 7).<sup>26</sup> Recent meta-analysis of the ability of GeneXpert MTB/RIF in the diagnosis of childhood PTB reported pooled sensitivity of 62% (95% CI 51-73) and a pooled specificity of 98% (95% CI 97-99).<sup>27</sup>

Bates et al reported no significant differences between specimens derived from sputum or liquid gastric washings in the GeneXpert MTB/RIF examination and concluded the use of liquid gastric washings can replace sputum specimens if they are not available.<sup>28</sup> In the study conducted by Nhu et al and Singh et al, stored sputum specimens were used instead of fresh sputum specimens.<sup>19,20</sup> In a sputum that was kept frozen and then thawed, the DNA will be damaged and affect the viscosity of sputum, thus giving bias.<sup>29</sup> Performing GeneXpert MTB/RIF examination twice on one specimen reportedly do not increase the rate of case detection. Repeated examination of GeneXpert MTB/RIF will increase the cost, even though there are still other supporting diagnostic examination. BBLK Surabaya's policy is to do single sputum GeneXpert MTB/RIF examination for each patient.<sup>14,20</sup>

The existence of GeneXpert MTB/RIF machines is not widely available in primary and secondary health facilities, therefore Sekkade et al conducted a study in Uganda and analyze the clinical characteristics associated with GeneXpert MTB/RIF positive results. It is intended to help health workers in limited medical care facilities to predict the likelihood of TB in a suspected TB children. Researchers reported some characteristics of the sample that has a tendency to get positive result of GeneXpert MTB/RIF, such as age group of  $> 5$  years, a positive TST result and a positive TB contacts.<sup>13</sup>

All sputum specimens with positive result of GeneXpert MTB/RIF and MTB sputum culture show sensitive result to rifampicin in this study. A study by Carriquiry et al on 130 patients aged  $> 18$  years in Peru in 2012 reported 100% (95% CI 61-100) and 91% (95% CI 88.7 to 100) for sensitivity and specificity respectively. Predictive result were

66.7% (95% CI 35.4 -87.9) and 100% (95% CI 88.7 to 100) for PPV and NPV respectively.<sup>29</sup> Some researchers have assessed the ability of GeneXpert MTB/RIF in detecting MTB with Rifampicin resistance, but the samples were too small and therefore cannot be assessed.<sup>14,20,30</sup>

<sup>32</sup> Tuberculosis is more progressive and fatal in children aged < 5 years old, while those aged ≥ 5 years old was associated with disease progression being "adult-type TB". Other than that this age group is the most common group of contracting TB in countries with high TB prevalence. This type of "adult-type TB" has the potential to cause extensive damage to lung parenchyma due to calcification and formation of cavities, and this age group is potentially infectious to the community.<sup>21</sup> In our study, the statistical test shows significant equivalence between GeneXpert MTB/RIF and MTB sputum cultures in both age groups. This means that GeneXpert MTB/RIF can be used interchangeably with MTB sputum culture to diagnose PTB in both age groups if there is no MTB sputum culture examination facilities. Beside, molecular methods with GeneXpert MTB/RIF also gives advantage of reading the results quickly (± 2 hours) so clinical decisions to initiate TB treatment can be accelerated. Sekadde et al and Nhu et al reported the same result for sensitivity and specificity of GeneXpert MTB/RIF for > 5 years old age group higher than group age < 5 years old.<sup>19,20</sup>

In our study, statistical analysis suggested there is a link between clinical manifestations and GeneXpert MTB/RIF or MTB sputum culture results. There is equivalence between GeneXpert MTB/RIF results and MTB culture result in the clinical TB group. GeneXpert MTB/RIF and MTB sputum cultures had lower sensitivity in diagnosing TB children clinically than molecularly. This is because children naturally had paucibacillary MTB although demonstrated clinical manifestations of TB and have symptoms improvement after treatment.

Symptoms of TB in children are not specific and more than 50% of children with TB are asymptomatic. Children with TB exhibiting clinical symptoms mostly will experience lung disorders, while 25% - 35% have extra-

<sup>27</sup> pulmonary disorders. Systemic disorders such as fever, night sweats, anorexia may also occur. The most common clinical symptoms are cough, body fatigue and weight loss. Specificity of clinical symptoms depends on the tightness of operational definitions used. However there is no cut-off for clinical symptoms that have been validated until now.<sup>22,31,32</sup>

The diagnosis of PTB in children cannot be established by clinical symptoms alone. Laboratory tests need to be done in children with or without clinical symptoms of PTB. However, the negative results of bacteriological examination does not exclude the possibility of TB disease.<sup>22</sup> Patients aged < 5 years, with a positive TST result and history of close contact with adult TB patients but do not show symptoms of PTB, were given INH prophylaxis of 7-15 mg /kg / day, once daily for 6 months, while patients that show symptoms of PTB were given TB drugs according to standard procedures.<sup>33</sup>

## CONCLUSION

GeneXpert MTB/RIF has a good sensitivity and specificity to diagnose pulmonary tuberculosis (PTB) in children which give parallel results with MTB sputum culture methods in aiding the diagnosis of PTB in children aged ≥ 3 months old - 14 years old with suspected PTB.

## <sup>26</sup> CONFLICT OF INTEREST

There is no conflict of interest of this study.

## ACKNOWLEDGEMENT

The author sincerely thank all members of Pediatric Department and staff in Dr. Soetomo Hospital for their permission, facilities and financial support.

## REFERENCES

1. Cuevas LE, Browning R, Bossuyt P, Casenghi M, Cotton MF, Cruz AT, et al. Evaluation of tuberculosis



- diagnostics in children: 2. methodological issues for conducting and reporting research evaluations of tuberculosis diagnostics for intrathoracic tuberculosis in children. consensus from an expert panel. *Journal of Infectious Diseases* 2012;205:209-15.
2. Edwards DJ, Kitelete F, Rie AV. Agreement between clinical scoring systems used for the diagnosis of pediatric tuberculosis in the HIV era. *Int J Tuberc Lung Dis* 2007;11:263-9.
  3. Starke JR. Tuberculosis (*Mycobacterium tuberculosis*). In: Kilegman RM, Stanton BF, Schor NF, Geme JWS, Behrman RE, editors. *Nelson textbook of pediatrics* 19th edition. Philadelphia: Elsevier, 2011;996-1011.
  4. WHO. Automated real-time nucleic acid amplification technology for rapid and simultaneous detection of tuberculosis and rifampicin resistance: Xpert MTB/RIF assay for the diagnosis of pulmonary and extrapulmonary TB in adults and children. policy update. 2013;23-33
  5. Kemenkes RI. Petunjuk Teknis Manajemen TB Anak. 2013;2-79.
  6. Marais BJ, Pai M. New approaches and emerging technologies in the diagnosis of childhood tuberculosis. *Paediatric Respiratory Reviews* 2007;8:124-33.
  7. Kulkarni S, Singh P, Memon A, et al. An in-house multiplex PCR test for the detection of *Mycobacterium tuberculosis*, its validation & comparison with a single target TB-PCR kit. *Indian J Med Res* 2012;135:788-94.
  8. Caws M, Wilson SM, Clough C, Drobniewski F. Role of IS6110-targeted PCR, culture, biochemical, clinical, and immunological criteria for diagnosis of tuberculous meningitis. *J Clin Microbiol* 2000;38:3150-5.
  9. Narayanan S, Parandaman V, Narayanan PR, Venkatesan P, Girish C, Mahadevan S, et al. Evaluation of PCR using TRC4 and IS6110 primers in detection of tuberculous meningitis. *J Clin Microbiol* 2001;39:2006-8.
  10. Chakravorty S, Tyagi JS. Novel multipurpose methodology for detection of mycobacteria in pulmonary and extrapulmonary specimens by smear microscopy, culture, and PCR. *J Clin Microbiol* 2005;43:2697-702.
  11. WHO. Rapid implementation of the Xpert MTB-RIF diagnostic test. Technical and operational 'how-to' practical considerations. Geneva, Switzerland: World Health Organization 2011;5-25
  12. Lawn SD, Nicol MP. Xpert®MTB/RIF assay: development, evaluation and implementation of a new rapid molecular diagnostic for tuberculosis and Rifampicin resistance. *Future Microbiol* 2011;6:1067-82.
  13. Sekadde MP, Wobudeya E, Joloba ML, et al. Evaluation of the Xpert MTB/RIF test for the diagnosis of childhood pulmonary tuberculosis in Uganda: a cross-sectional diagnostic study. *BMC Infect Dis* 2013;13:133-41.
  14. Nicol MP, Workman L, Isaacs W, et al. Accuracy of the Xpert MTB/RIF test for the diagnosis of pulmonary tuberculosis in children admitted to hospital in Cape Town, South Africa: a descriptive study. *The Lancet* 2011;1-6.
  15. Nataprawira HMD, Kartasasmita CB, Rosmayudi O, et al. Diagnosis of pediatric tuberculosis using the Indonesian national consensus for pediatric tuberculosis. *Pediatr Indones* 2001;41:185-90
  16. Chandra RK. Nutrition and the immune system : an introduction. *Am J clin Nutr* 1997;66:460-3
  17. Cunningham-Rundles S, Moon A, McNeeley DF. Malnutrition and host defense. In *Nutrition in Pediatrics*. 4th ed. Hamilton, Ontario, Canada: BC Decker Inc, 2008:261-72
  18. Giang DC, Duong TN, Ha DTM, et al. Prospective evaluation of GeneXpert for the diagnosis of HIV-negative pediatric TB cases. *BMC Infect Dis* 2015;15:70-80
  19. Singh S, Singh A, Prajapati S, et al. Xpert MTB/RIF assay can be used on archived gastric aspirate and induced sputum samples for sensitive diagnosis of pediatric tuberculosis. *BMC Microbiol* 2015;15:191-201
  20. Nhu NTQ, Ha DTM, Anh ND, et al. Evaluation of Xpert MTB/RIF and MODS assay for the diagnosis of pediatric tuberculosis. *BMC Infectious Diseases* 2013;13:31-40.
  21. Marais BJ, Gie RP, Schaaf HS, et al. The natural history of childhood intra-thoracic tuberculosis: a critical review of literature from the pre-chemotherapy era. *Int J Tuberc Lung Dis* 2004;8:392-402
  22. WHO. Guidance for national tuberculosis programmes on the management of tuberculosis in children. 2014;2:21-74
  23. Marlowe EM, Novak-Weekley SM, Cumpio J, et al. Evaluation of the cepheid Xpert MTB/RIF assay for direct detection of *Mycobacterium tuberculosis* complex in respiratory specimens. *J Clin Microbiol* 2011;49:1621-3.
  24. Kemenkes RI. Petunjuk teknis pemeriksaan biakan, identifikasi, dan uji kepekaan *Mycobacterium tuberculosis* pada media padat. 2012:10-44
  25. WHO. Automated real-time nucleic acid amplification technology for rapid and simultaneous detection of tuberculosis and rifampicin resistance: Xpert MTB/RIF assay for the diagnosis of pulmonary and extrapulmonary TB in adults and children. policy update. 2013:23-33
  26. Chang K, Lu W, Wang J, et al. Rapid and effective diagnosis of tuberculosis and Rifampicin resistance with Xpert MTB/RIF assay: a meta-analysis. *J Infect* 2012;64: 580-8.
  27. Detjen AK, Keil T, Roll S, et al. Interferon-g release assays improve the diagnosis of tuberculosis and nontuberculous mycobacterial disease in children in a

- country with a low incidence of tuberculosis. *Clinical Infectious Diseases* 2007;45:322-28.
28. Bates M, O'Grady J, Maeurer M, et al. Assessment of the *Xpert MTB/RIF* assay for diagnosis of tuberculosis with gastric lavage aspirates in children in sub-Saharan Africa: a prospective descriptive study. *The Lancet Infect Dis* 2013;13:36-42.
  29. Carriquiry G, Otero L, Gonzalez-Lagos E, et al. A diagnostic accuracy study of *xpert MTB/RIF* in HIV positive patients with high clinical suspicion of pulmonary tuberculosis in Lima, Peru. *PLoS ONE* 2012;7:1-7.
  30. Theron G, Peter J, Zyl-Smit R, et al. Evaluation of the *Xpert MTB/RIF* assay for the diagnosis of pulmonary tuberculosis in a high HIV prevalence setting. *Am J Respir Crit Care Med* 2011;184:132-40.
  31. Hesselting AC, Schaaf HS, Gie RP, et al. A critical review of diagnostic approaches used in the diagnosis of childhood tuberculosis. *Int J Tuberc Lung Dis* 2002;6:1038-45.
  32. Shingadia D, Burgner D. *Mycobacterial infections*. In: Taussig LM, Landau LI, Souef PNL, Martinez FD, Morgan WJ, Sly PD, editors. *Pediatric respiratory medicine* 2nd edition. Philadelphia: Elsevier 2008:597-614.
  33. Kemenkes RI. *Petunjuk Teknis Manajemen TB Anak*. 2013:2-79.

# GeneXpert MTBRIF and Mycobacterium tuberculosis Sputum Culture in Establishing the Diagnosis of Pulmonary Tuberculosis

## ORIGINALITY REPORT

24%

SIMILARITY INDEX

18%

INTERNET SOURCES

19%

PUBLICATIONS

0%

STUDENT PAPERS

## PRIMARY SOURCES

- 1** Matthew Bates, Justin O'Grady, Markus Maeurer, John Tembo et al. "Assessment of the Xpert MTB/RIF assay for diagnosis of tuberculosis with gastric lavage aspirates in children in sub-Saharan Africa: a prospective descriptive study", *The Lancet Infectious Diseases*, 2013  
Publication 2%
- 2** [www.biomedcentral.com](http://www.biomedcentral.com)  
Internet Source 2%
- 3** Bates, Matthew, Justin O'Grady, Markus Maeurer, John Tembo, Lophina Chilukutu, Chishala Chabala, Richard Kasonde, Peter Mulota, Judith Mzyece, Mumba Chomba, Lukundo Mukonda, Maxwell Mumba, Nathan Kapata, Andrea Rachow, Petra Clowes, Michael Hoelscher, Peter Mwaba, and Alimuddin Zumla. "Assessment of the Xpert MTB/RIF assay for diagnosis of tuberculosis with gastric lavage aspirates in children in sub-Saharan Africa: a

prospective descriptive study", The Lancet Infectious Diseases, 2013.

Publication

---

4	<a href="https://www.123dok.com">123dok.com</a> Internet Source	1%
5	<a href="https://www.ncbi.nlm.nih.gov">www.ncbi.nlm.nih.gov</a> Internet Source	1%
6	Tingyu Tang, Fang Liu, Xiaoling Lu, Qingdong Huang. " Evaluation of GeneXpert MTB/RIF for detecting in a hospital in China ", Journal of International Medical Research, 2017 Publication	1%
7	<a href="https://pt.scribd.com">pt.scribd.com</a> Internet Source	1%
8	<a href="https://e-journal.unair.ac.id">e-journal.unair.ac.id</a> Internet Source	1%
9	Ghulam Rasool, Arif Muhammad Khan, Raza Mohy-Ud-Din, Muhammad Riaz. " Detection of in AFB smear-negative sputum specimens through MTB culture and GeneXpert MTB/RIF assay ", International Journal of Immunopathology and Pharmacology, 2019 Publication	1%
10	<a href="https://hdl.handle.net">hdl.handle.net</a> Internet Source	1%

---

11 Youness Mechal, Elmostafa Benaissa, Nadia El mrimar, Yassine Benlahlou et al. "Evaluation of GeneXpert MTB/RIF system performances in the diagnosis of extrapulmonary tuberculosis", BMC Infectious Diseases, 2019  
Publication 1%

---

12 [apps.who.int](https://apps.who.int)  
Internet Source 1%

---

13 [www.nice.org.uk](https://www.nice.org.uk)  
Internet Source 1%

---

14 [eprints.soton.ac.uk](https://eprints.soton.ac.uk)  
Internet Source 1%

---

15 [repository.unair.ac.id](https://repository.unair.ac.id)  
Internet Source 1%

---

16 Nicol, M.P.. "Accuracy of the Xpert MTB/RIF test for the diagnosis of pulmonary tuberculosis in children admitted to hospital in Cape Town, South Africa: a descriptive study", The Lancet Infectious Diseases, 201111  
Publication 1%

---

17 N. N. You, L. M. Zhu, G. L. Li, L. Martinez, W. Lu, Q. Liu, H. T. Yang. "A tuberculosis school outbreak in China, 2018: reaching an often overlooked adolescent population", Epidemiology and Infection, 2019  
Publication <1%

---

18

[idoc.pub](http://idoc.pub)

Internet Source

<1%

---

19

[doaj.org](http://doaj.org)

Internet Source

<1%

---

20

[journals.lww.com](http://journals.lww.com)

Internet Source

<1%

---

21

Elizabeth Whittaker, Heather J Zar. "Promising directions in the diagnosis of childhood tuberculosis", Expert Review of Respiratory Medicine, 2014

Publication

<1%

---

22

[creativecommons.org](http://creativecommons.org)

Internet Source

<1%

---

23

[www.apjpch.com](http://www.apjpch.com)

Internet Source

<1%

---

24

K. Myo, M. Zaw, T. L. Swe, Y. Y. Kyaw, T. Thwin, T. T. Myo, K. O. Aye, A. A. Myint. "Evaluation of Xpert MTB/RIF assay as a diagnostic test for pulmonary tuberculosis in children in Myanmar ", The International Journal of Tuberculosis and Lung Disease, 2018

Publication

<1%

---

25

Ogechukwu Egini, Edouard Guillaume, Titilope Adeyemo, Chiemeziem Nwanyanwu, Fnu Shweta, Eric Jaffe. "In vitro regulation of reactive oxygen species formation in red blood

<1%

cells of homozygous sickle cell patients using  
Vitamin C", F1000Research, 2017

Publication

26

[banglajol.info](http://banglajol.info)

Internet Source

<1%

27

[vdocuments.mx](http://vdocuments.mx)

Internet Source

<1%

28

[www.hpv2009.org](http://www.hpv2009.org)

Internet Source

<1%

29

Chisti, Mohammad Jobayer, Stephen M. Graham, Trevor Duke, Tahmeed Ahmed, Hasan Ashraf, Abu Syed Golam Faruque, Sophie La Vincente, Sayera Banu, Rubhana Raqib, and Mohammed Abdus Salam. "A Prospective Study of the Prevalence of Tuberculosis and Bacteraemia in Bangladeshi Children with Severe Malnutrition and Pneumonia Including an Evaluation of Xpert MTB/RIF Assay", PLoS ONE, 2014.

Publication

<1%

30

Mercedes C Becerra, Soumya Swaminathan. "Commentary: A targets framework: Dismantling the invisibility trap for children with drug-resistant tuberculosis", Journal of Public Health Policy, 2014

Publication

<1%

31

[docplayer.info](http://docplayer.info)

Internet Source

&lt;1%

32

[malariajournal.biomedcentral.com](http://malariajournal.biomedcentral.com)

Internet Source

&lt;1%

33

[www.msh.org](http://www.msh.org)

Internet Source

&lt;1%

34

Kamon Kawkitinarong, Gompol Suwanpimolkul, Pairaj Kateruttanakul, Weerawat Manosuthi et al. "Real-Life Clinical Practice of Using the Xpert MTB/RIF Assay in Thailand", *Clinical Infectious Diseases*, 2017

Publication

&lt;1%

35

Luxia Kong, Bangruan Xie, Qian Liu, Lijuan Hua, Anish Bhusal, Chen Bao, Jiannan Hu, Shuyun Xu. "Application of acid-fast staining combined with GeneXpert MTB/RIF in the diagnosis of non-tuberculous mycobacteria pulmonary disease", *International Journal of Infectious Diseases*, 2021

Publication

&lt;1%

36

Mohammed Ahmed Yassin, Luelseged Takele, Sahlemariam Gebresenbet, Emebet Girma, Meskele Lera, Ersido Lendebo, Luis E. Cuevas. "HIV and Tuberculosis Coinfection in the Southern Region of Ethiopia: A Prospective Epidemiological Study", *Scandinavian Journal of*

&lt;1%



37

[erj.ersjournals.com](http://erj.ersjournals.com)

Internet Source

<1%

38

[old.biomedcentral.com](http://old.biomedcentral.com)

Internet Source

<1%

39

M Held, M Laubscher, L Workman, H J Zar, R Dunn. "Diagnostic accuracy of GeneXpert MTB/RIF in musculoskeletal tuberculosis: High sensitivity in tissue samples of HIV-infected and HIV-uninfected patients", South African Medical Journal, 2017

Publication

<1%

40

Marilda Casela, Silvânia Maria Andrade Cerqueira, Thais de Oliveira Casela, Mariana Araújo Pereira et al. "Rapid molecular test for tuberculosis: impact of its routine use at a referral hospital", Jornal Brasileiro de Pneumologia, 2018

Publication

<1%

41

Toyin Togun, Beate Kampmann, Madhukar Pai. "Diagnosis of Childhood Tuberculosis", Elsevier BV, 2017

Publication

<1%

42

[academic.oup.com](http://academic.oup.com)

Internet Source

<1%

43	<a href="http://archive.lstmed.ac.uk">archive.lstmed.ac.uk</a> Internet Source	<1%
44	<a href="http://article.sapub.org">article.sapub.org</a> Internet Source	<1%
45	<a href="http://test.dovepress.com">test.dovepress.com</a> Internet Source	<1%
46	<a href="http://www.futuremedicine.com">www.futuremedicine.com</a> Internet Source	<1%
47	<a href="http://www.health.gov.au">www.health.gov.au</a> Internet Source	<1%
48	<a href="http://www.hpv.min-saude.pt">www.hpv.min-saude.pt</a> Internet Source	<1%
49	<a href="http://www.redalyc.org">www.redalyc.org</a> Internet Source	<1%
50	Anna Cristina Calçada Carvalho, Claudete Aparecida Araújo Cardoso, Terezinha Miceli Martire, Giovanni Battista Migliori et al. "Epidemiological aspects, clinical manifestations, and prevention of pediatric tuberculosis from the perspective of the End TB Strategy", <i>Jornal Brasileiro de Pneumologia</i> , 2018 Publication	<1%
51	Jose Domínguez, Irene Latorre, Neus Altet, Lourdes Mateo, Malú De Souza-Galvão, Juan	<1%

Ruiz-Manzano, Vicente Ausina. "IFN- $\gamma$ -release assays to diagnose TB infection in the immunocompromised individual", Expert Review of Respiratory Medicine, 2014

Publication

---

52

"Tuberculosis of the Central Nervous System", Springer Science and Business Media LLC, 2017

Publication

---

<1%

53

L. Rigouts. "Clinical practice", European Journal of Pediatrics, 04/25/2009

Publication

---

<1%

---

Exclude quotes Off

Exclude matches Off

Exclude bibliography Off

# GeneXpert MTBRIF and Mycobacterium tuberculosis Sputum Culture in Establishing the Diagnosis of Pulmonary Tuberculosis

---

## GRADEMARK REPORT

---

FINAL GRADE

**/100**

GENERAL COMMENTS

**Instructor**

---

PAGE 1

---

PAGE 2

---

PAGE 3

---

PAGE 4

---

PAGE 5

---

PAGE 6

---

PAGE 7

---

PAGE 8

---

PAGE 9

---