

Cytotoxic effect of natural cuttlefish bone xenograft: an in vitro and in vivo study

by Komang Agung Irianto Suryaningrat

Submission date: 17-Mar-2021 12:02PM (UTC+0800)

Submission ID: 1535109043

File name: 3_Cytotoxic_effect_of_natural_cuttlefish.pdf (1.09M)

Word count: 4270

Character count: 22597

Cytotoxic effect of natural cuttlefish bone xenograft: an *in vitro* and *in vivo* study

Komang Agung Irianto, Suyenci Limbong



4

pISSN: 0853-1773 • eISSN: 2252-8083
<https://doi.org/10.13181/mji.oa.203141>
Med J Indones. 2020;29:136–42

Received: September 15, 2018

Accepted: May 13, 2020

Authors' affiliations:

Department of Orthopaedic and
Traumatology, Faculty of Medicine,
Universitas Airlangga, Dr. Soetomo
Hospital, Surabaya, Indonesia

Corresponding author:

Komang Agung Irianto
Department of Orthopaedic and
Traumatology, Dr. Soetomo Hospital,
Jalan Prof. Dr. Moestopo No. 6–8,
Surabaya 60286, Indonesia
Tel/Fax: +62-31-5501481
E-mail: komang168@yahoo.com

ABSTRACT

BACKGROUND Commercialized synthetic bone grafts are commonly used to replace the bone defect. Cuttlefish bone is naturally available and widely studied, but the specific cytotoxicity test has not been conducted. This study aimed to evaluate the cytotoxicity of the xenograft compared to commercial grafts.

METHODS We performed an *in vitro* test evaluating the viability of human mesenchymal stem cells (hMSCs) when cultured for 48 hours with the tested materials (cuttlefish bone graft and fabricated PerOssal®). The trypsinized mitochondrial activity of the viable hMSC was assayed based on colorimetry of the formazan color change. The tested material was considered nontoxic if >70% of the hMSCs were viable. The *in vivo* cytotoxic effect was evaluated by implanting the graft material in the femoral muscle of New Zealand (NZ) white rabbits. Nine rabbits were used in each test (cuttlefish bone, PerOssal®, and NaCl 0.9%). The systemic acute pyrogenic effect was evaluated based on 72 hours body weight changes and rectal temperature changes every 30 min in the first 3 hours and 72 hours post-implantation.

RESULTS The mean percentage of hMSC viability when cultured with cuttlefish bone graft and PerOssal® was comparable (93.47% and 105.37%, respectively, $p = 0.240$). The *in vivo* cytotoxicity on NZ rabbit was similar between all tested materials, as shown by the minor changes in body weight (<10% body weight, $p = 0.391$) and rectal temperature (<0.5°C, $p = 0.127$).

CONCLUSIONS Cuttlefish bone xenograft and fabricated PerOssal® have a similar non-cytotoxic effects on hMSCs and non-pyrogenic systemic effects on rabbits.

KEYWORDS bone substitute, cuttlefish bone, cytotoxicity tests, xenograft

Bone grafting is an operative procedure performed to replace a bone defect.¹ This procedure is the second most frequently performed procedure after blood transfusion.² According to a previous study, a bone graft may be an autograft, allograft, and xenograft.³ Currently, an autograft is still the gold standard for repairing bone defects, as it has osteoconductive, osteoinductive, and osteogenesis

properties and provides growth factors with a relatively low risk of infection.¹ However, bone grafting has several disadvantages, including postoperative pain, hematoma formation, infection, neurovascular trauma, and esthetic deformity, and bone grafts are limited to small defects.⁴ Moreover, allograft and xenograft are developed because of the limited availability of autograft. Allograft requires a

1

Copyright © 2020 Authors. This is an open access article distributed under the terms of the Creative Commons Attribution-NonCommercial 4.0 International License (<http://creativecommons.org/licenses/by-nc/4.0/>), which permits unrestricted non-commercial use, distribution, and reproduction in any medium, provided the original author and source are properly cited. For commercial use of this work, please see our terms at <https://mji.ui.ac.id/journal/index.php/mji/copyright>.

rigorous preparation to reduce a possible transmission of disease, thus losing most of its osteoconductive, osteoinductive, and osteogenesis capabilities.^{3,5}

Bovine xenograft has been widely used.^{5,6} It plays a role in the development of subsequent xenografts. One of the origins of xenografts is cuttlefish bone; it has a lower production cost, high worldwide availability, and naturally procured. Cuttlefish bone xenograft acts as a scaffold with a high degree of porosity, which is ideal for bone grafting. It has appropriate interconnectivity, suitable for the biologic activity of bone growth, and revascularization.⁷ Following the application of cuttlefish extract, the callus grows thicker with higher osteoblast proliferation, ultimately resulting in better bone healing process.⁶⁻⁸ However, before the natural xenograft could be introduced for a human clinical trial, the toxicity of cuttlefish bone xenograft *in vitro* and *in vivo* needs to be evaluated to prove the safety of its application.

Synthetic bone grafts are currently used owing to their availability and high reproducibility.^{3,4} The bioavailability of synthetic bone graft is widely investigated, which largely depends on the architectural features of the composite material.^{1,5} PerOssal® is a synthetic bone replacement material that consists of nanocrystalline hydroxyapatite and calcium sulfate. It possesses osteostimulation and osteoinductive properties, which may promote bone healing. In addition, this material may be used to administer antibiotics into a bone defect.

Several studies are evaluating the *in vitro* and *in vivo* potency of cuttlefish bone xenograft as a bone replacement or as a scaffold.⁶⁻⁹ Nonetheless, before a medical device is allowed for human application, several biocompatibility safety assessments according to U.S. Food and Drug Administration, Japanese guidelines, and International Organization for Standardization should be met.¹⁰ The *in vitro* test would be assessed by the viability of human mesenchymal stem cells (hMSCs) after cultured with the tested material. New Zealand (NZ) white rabbit is commonly used as an animal *in vivo* cytotoxicity test model because of its docility and easy-to-maintain health.¹¹ The systemic cytotoxic effect was evaluated based on acute systemic responses such as an alteration in body weight due to diarrhea and loss of appetite and on the pyrogenicity of the implanted material based on rectal temperature changes before and after culture.^{11,12} The rabbit pyrogen test was

used for biocompatibility because rabbits respond to pyrogens regardless of their source.¹¹ Thus, this study was aimed to compare the *in vitro* and *in vivo* toxicity of cuttlefish bone xenograft with a commercially available bone substitute (PerOssal®).

METHODS

This was an experimental study performed at the Laboratory of Tropical Disease Center of Universitas Airlangga, from February through August 2018.

Study design

The *in vitro* cell viability test was performed using hMSC cell line incubated for 48 hours in 96-well plates with the tested material.¹³ Cell viability was assayed by colorimetry which measured the optical density (OD) of the cell suspension and compared with the blank well as the control as described by 3-(4,5-dimethylthiazol-2-yl)-2,5-diphenyltetrazolium bromide (MTT) assay.^{14,15} This process was repeated in triplicate.

The *in vivo* cytotoxic test was a controlled animal laboratory experiment. The minimum number of required samples was calculated using a sample size formula.¹⁶ In this study, 27 NZ white rabbits (*Oryctolagus cuniculus*) were divided into three groups, with nine rabbits for each group (cuttlefish bone, PerOssal®, and control). This study was fully compliant with animal research on the reporting of *in vivo* experiments (ARRIVE) criteria.¹² The study protocol was approved by the Animal Care and Use Committee, Universitas Airlangga, Indonesia (certificate number: 0326/KEPK/V/2018).

Handling of animals

All rabbits were housed in the animal care laboratory and were well-taken care according to the standards of the National Institute of Health.¹² The inclusion criteria were as follows: male rabbits, age 6–9 months, healthy without any disability, ambulatory, and body weight of 2,500–3,000 g. The exclusion criteria were surgical wound infection and death. The rabbits were housed individually in separate cages (100 × 60 × 75 cm) with following environmental conditions: temperature of 21 ± 2°C, humidity of 60 ± 10%, and lighting of 350 lux intensity with 12-hours dark and light cycle. All rabbits were fed regularly and provided with an unlimited supply of water.

Preparation of hMSCs

The human bone marrow was extracted from the femoral shaft of a patient undergoing femoral nailing. The patient has consented to have his/her bone marrow tissue extracted for further research purposes. The bone marrow was put into a 50-ml conical tube, phosphate-buffered saline (PBS) 1:1 was added, suspended to be homogenous, and then carefully layered on 5 ml of Ficoll in a conical tube. It was then centrifuged for 30 min at 1,600 rpm, at 26°C until separated into four layers. The second layer was a buffy coat, which looked like a cloudy ring. The buffy coat was put in a 15-ml conical tube, washed with 10 ml of PBS, and centrifuged for 5 min at 1,600 rpm. The pellet was resuspended in culture media until becoming homogenous, placed in a 10-cm plate, labeled, and then incubated in a CO₂ incubator.¹³

Preparation of cuttlefish bone xenograft

Aragonite (CaCO₃) was obtained from the lamellae of cuttlefish bone powder, using high-energy milling three-dimensional motion or mortar. It was heated in a furnace at 350°C for 3 hours. The CaCO₃ content was quantified using X-ray diffraction (XRD). Then, 100 g of CaCO₃ was added to 1 liter of distilled water to acquire a 1 M solution. The 0.6 M NH₄H₂PO₄ solution was prepared by dissolving 69 g of the compound in 1 liter of distilled water.^{7,16} The hydrothermal reaction to produce hydroxyapatite (HA) was as follows: 10 CaCO₃ + 6 NH₄H₂PO₄ + 2H₂O → Ca₁₀(PO₄)₆(OH)₂ + 3 (NH₄)₂CO₃ + 7 H₂CO₃. The reaction was carried out by mixing two

previous solutions using a magnetic stirrer for 30 min. The mixture was then heated at 200°C for 12 hours. The mixture was then cooled at room temperature and washed with a distilled water using a magnetic stirrer until it was separated from the water and a neutral pH was achieved. The final washing process contained methanol to limit the agglomeration of HA during the drying process.¹⁷ The washed sample was filtered through a filter paper and dried in an electric oven at 50°C until completely dry. The HA sample that was formed during this process was labeled. The crystal structures were characterized using XRD (Figure 1). The morphology and the ratio Ca/PO₄ were quantified using scanning electron microscope-energy dispersive X-ray spectroscopy.^{7,17}

The synthetic bone graft was PerOssal® (AAP biomaterials GmbH, Germany; product number: 42000000-AKS-004699091). It is in powder form and ready to use. Several other similar products are available from different brands. However, this product is already approved by the Indonesian Food and Drug Regulator (registration number of AKL 21302802278).

Surgical procedure and animal sample preparation

The animals were anesthetized with ketamine (40 mg/kg) and xylazine (5 mg/kg) intramuscularly and placed prone on a warm pad. The right thigh was disinfected and draped aseptically. The right quadriceps muscle of the control group was injected with NaCl 0.9% alone (0.5 cc), while the cuttlefish bone group was injected with cuttlefish bone extract + NaCl 0.9% mixture (0.5 cc), and the Perossal® group was injected with Perossal® + NaCl 0.9% mixture. The bone graft material was injected in the muscle to introduce the possible clinical cytotoxic effect (rectal temperature and body weight changes to represent the health status of the animal sample) as shown in Figure 2. Given the aseptic process, within 72 hours of evaluation, the injection site should not show an inflammatory reaction and functional lesions. If a wound infection is present, the sample will be excluded.

A reduction in body weight is one of the indicators of systemic health status in response to cytotoxicity, which usually appears with diarrhea or loss of appetite.¹² Body weight was measured before implantation and at 72 hours after implantation. Rectal temperatures may signify immunological systemic reaction in the *in vivo* models.¹² The rabbit

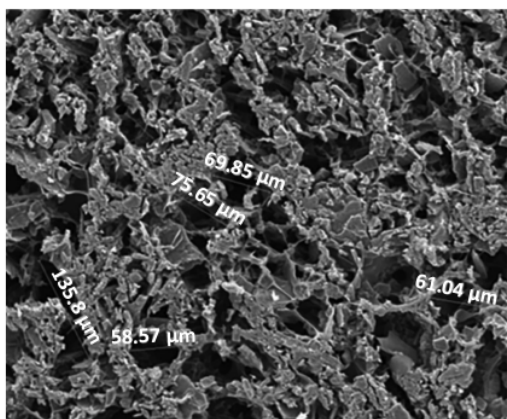


Figure 1. Crystal structure of the cuttlefish hydroxyapatite under the scanning electron microscope

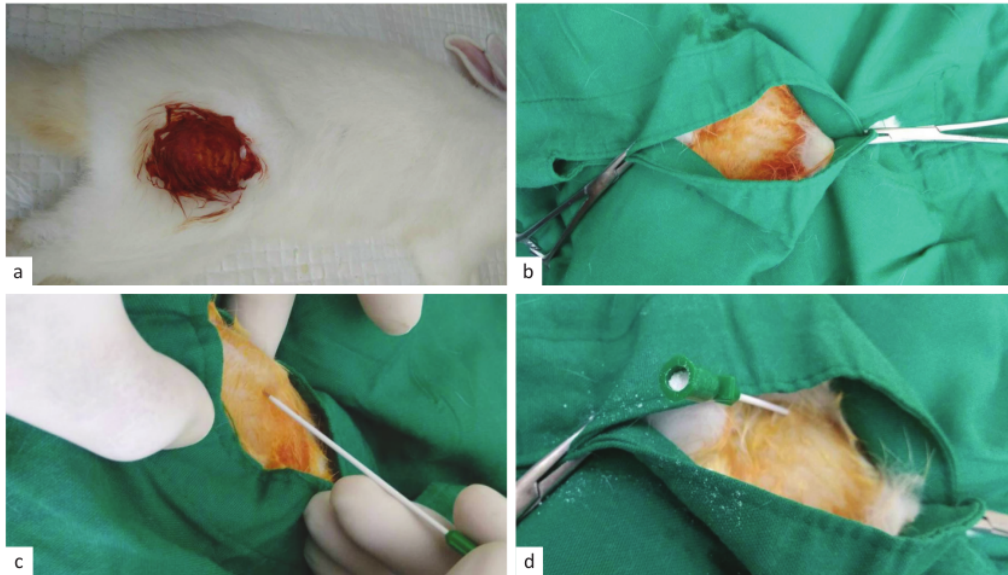


Figure 2. Surgical procedure (sequential clockwise from the upper left). (a) Disinfection; (b) draping; and (c and d) injection of the mixtures (NaCl 0.9% for the control group, cuttlefish bone extract + NaCl mixture for the cuttlefish bone group, and PerOssal® + NaCl mixture for the PerOssal® group)

pyrogenicity test was conducted by recording the rectal temperatures at baseline (30 min before injection) and 30-min intervals up to 3 hours post-implantation and 72 hours thereafter. The tested material is not pyrogenic if no rabbit showed an individual rise in temperature of $\geq 0.5^{\circ}\text{C}$ above its respective control temperature. The requirements for the absence of pyrogenicity are met if not more than three of the nine rabbits have an individual rise in temperature of $\geq 0.5^{\circ}\text{C}$ or if the sum of the maximum rise in temperature of the nine rabbits does not exceed 3.3°C .¹⁰⁻¹² During housing, the health status of the animals were monitored three times daily. No adverse events should be observed in any of the animals.¹²

Cytotoxicity test using MTT assay

The viability of hMSCs *in vitro* was evaluated by MTT assay methods on each group. Cell viability was measured based on formazan color change. The cells at 80% confluent were trypsinized and plated in 96-well culture plates (1×10^4 cells/well). Each well contained 100 μl of cell suspension, and the plates were incubated for 24 hours at 37°C in a 5% CO_2 environment. The media was removed from each well after 24 hours. Subsequently, 100 μl of eluent

from the 0.5, 1, 5, 25, 50, 100, or 200 $\mu\text{g/ml}$ cuttlefish bone powder, PerOssal® powder, positive/negative control, and blank wells was placed into the 96-well culture plates (3 \times 8 wells/test material), as shown in Figure 3.

Cell viability was assessed after the incubation. The experiments were repeated three times. Then, 100% viability was determined by the mean OD of the control group (hMSC with media). The OD of the blank well was considered 0% cell viable. The results for the experimental, positive, and negative control groups were normalized to the blank. The relative cell count ratio was calculated using the following formula (Equation 1):¹³

$$\% \text{ Cell viability} = \frac{\text{optical density sample} - \text{optical density blank}}{\text{optical density control} - \text{optical density blank}} \times 100\% \quad (1)$$

Statistical analysis

Acquired data were tested by normality test using the Shapiro-Wilk test. If the data were distributed normally, the one-way analysis of variance was used. Otherwise, nonparametric Mann-Whitney test was performed. For repeated measurements, we used the general linear model. If data were not distributed normally, Kruskal-Wallis test was performed. All

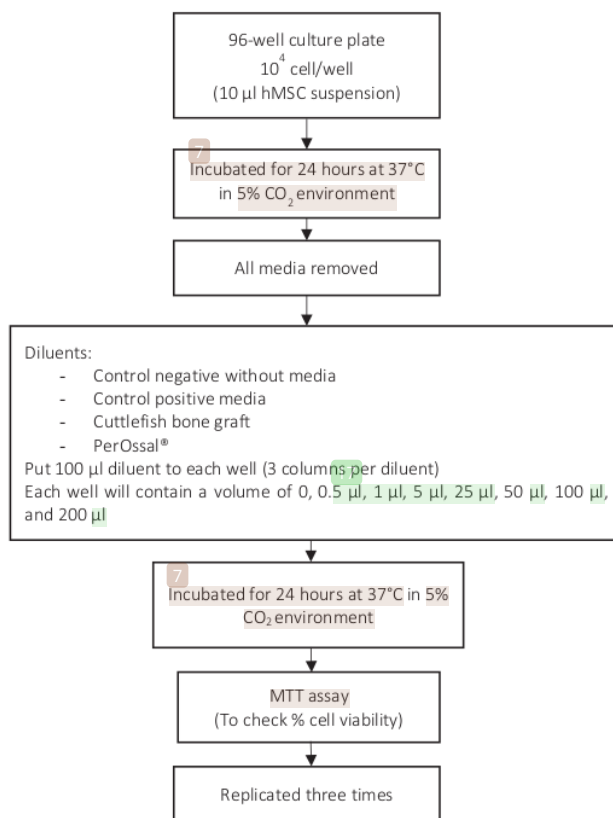


Figure 3. Flow chart of the cytotoxic test using hMSCs and MTT assay. hMSCs=human mesenchymal stem cells; MTT=3-(4,5-dimethylthiazol-2-yl)-2,5-diphenyltetrazolium bromide

statistical analysis was done using the SPSS software version 18.0 (IBM Corp., USA).

RESULTS

As shown in Table 1, the mean percentage of cell viability in the cuttlefish bone and PerOssal® groups were similar. All nine tests of both materials showed >70% viability which meant that cuttlefish bone and PerOssal® grafts were considered nontoxic for medical device assessment, mean (standard deviation) cell viability: cuttlefish bone at 0.9347 (±0.09) and PerOssal® at 1.0573 (±0.32); $p = 0.240$ as shown in Figure 4.¹

Post-hoc analysis using Bonferroni indicated that cuttlefish bone and PerOssal® bone grafts endured similar acute systemic outcome ($p = 0.864$). Body weights shifted less than 5% (100 g of 2,500 g).

The average rectal temperature of the nine rabbits before implantation was within the healthy normal range. Thus, the pyrogenicity test may be

performed. The average rectal temperatures were done after implantation taken at 30-min interval for the first 3 hours and 72 hours. The three tested implant materials showed similar rectal temperature changes before and after the implantation. Since the temperature changes were within the range of 0.5°C, all cuttlefish bone, PerOssal®, and control groups (NaCl 0.9%) passed the pyrogenicity test ($p = 0.127$) (Table 1).

DISCUSSION

In vivo cytotoxic test was conducted to ensure the safety of the implant for human application. The short-term *in vitro* study is not comparable to an *in vivo* study since most measurements of the toxic effects was only conducted within 12–24 hours after exposure to a toxic substance.^{14,18} *In vivo* reactions tend to be more complex and possibly performed beyond 24 hours. Thus, it is imperative to use animal samples to test bonegrafts before human clinical

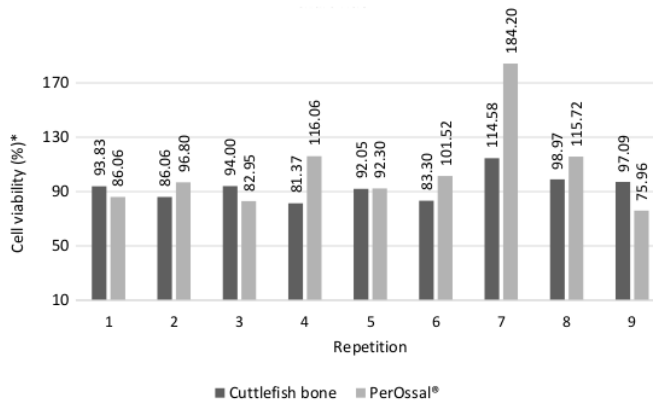


Figure 4. Human mesenchymal stem cell (hMSC) viability test using cuttlefish bone and PerOssal® materials

*The viability of the cells was quantified based on optical density of the cell solution by colorimetry assay and the reference was normal control as 100%, hence it may be exceeded 100%

Table 1. Toxicity effect of cuttlefish bone and PerOssal®

	<i>In vitro</i> hMSC		<i>In vivo</i> (9 NZ rabbits per test)*	
	Cell viability (%), mean (SD)	Body weight (g), mean (SD)	Temperature before intervention (°C), mean (SD)	Temperature after intervention (°C), mean (SD)
Control	100	46.67 (10.00)	39.55 (0.31)	39.46 (0.25)
PerOssal®	105.73 (32)	45.55 (11.30)	39.36 (0.35)	39.21 (0.29)
Cuttlefish	93.47 (9)	40.00 (11.18)	39.48 (0.86)	39.43 (0.17)
<i>p</i>	0.240 [†]	0.391	0.488 [‡]	0.127 [†]

hMSC=human mesenchymal stem cell; NZ=New Zealand; SD=standard deviation

*Every rabbit weighted 2,500–3,000 g; [‡]Kruskal–Wallis test; [†]one-way analysis of variance

trials.⁸ The animal selection criteria require using animals with similar physiology and pathology with humans.^{9,11,12} In this study, all NZ white rabbits survived until the end of the study, which is in accordance with the nature of NZ white rabbit.¹¹

In this study, cuttlefish bonegraft as xenograft has similar cytotoxicity with PerOssal® when incubated with hMSC from human femoral bone marrow. A similar result using osteoblast cell line from rat calvaria was reported by Vajrabhaya et al¹⁴ where the decline in viability percentage of the cell was similar in the cuttlefish bone powder group and control group. In this study, the proliferative effect of the cuttlefish bone also showed good results. This fact may be explained by the properties of HA from cuttlefish bone that is quite similar to those found in calcium-based tissues in humans.¹⁹ Other studies found that the cuttlefish bonegraft has the potency to induce osteoblast differentiation and proliferation.^{7,20,21} This circumstance is postulated because of the presence of Mg²⁺ on the material, which was proven to promote osteogenesis.^{22,23} Thus, cuttlefish bone powder poses low toxicity and ability

to promote osteogenesis, prompting its possible application as bone graft material.^{14,24}

The cytotoxic effect *in vivo* in laboratory animals may present as diarrhea or severe loss of appetite which are usually related to reduced body weight or even death. The rectal temperature signifies an immunological inflammatory reaction in *in vivo* models. The acute systemic toxicity examination in this study for PerOssal®, cuttlefish bone, and control groups showed minor acute systemic toxicity reaction and pyrogenicity, which is consistent with the previous experiment with rat fracture model.^{6,7,19}

Dogan and Okumus⁸ reported that cuttlefish bone xenograft did not affect the physiological measurement of the NZ rabbit after 24 weeks of implantation, allowing the cuttlefish bone to become bone xenograft potential material. They also showed no increase in inflammatory cells histologically; thus, the graft does not cause inflammatory reactions.^{20,22} This is in agreement to the implantation of cuttlefish bone extract into *Rattus norvegicus* with fractured tibial bone that resulted in thicker callus formation and higher proliferation of osteoblast without

significant inflammatory reaction.⁶ This information signified that an accelerated bone healing process may be affected by the application of cuttlefish bone extract.^{24,25}

This study has several limitations. The authors only examine the toxicity effect of cuttlefish bone graft in the acute phase (72 hours). The only aspects explored were the pyrogenicity and the effect of the graft to the animals' body weight, and the authors did not evaluate chronic toxicity. Further bioavailability tests need to be conducted to evaluate the effect of cuttlefish bone xenograft in terms of its genotoxicity, sensitization, irritation, chronic toxicity, implant rejection, and carcinogenicity. In conclusion, the cuttlefish bone xenograft has similar nontoxic effect on hMSCs from bone marrow to the fabricated synthetic bone graft PerOssal®, and both passed the pyrogenicity rabbit test.

Conflict of Interest

The authors affirm no conflict of interest in this study.

Acknowledgment

The authors would like to acknowledge Oen Sindrawati, Yudhistira Pradnyan Klopung, and William PS for editing and revising the final text.

Funding Sources

None.

REFERENCES

- Hung NN. Basic knowledge of bone grafting. In: Zorzi A, editor. Bone grafting. London: IntechOpen; 2012. p. 11–38. Available from: <https://www.intechopen.com/books/bone-grafting/basic-knowledge-of-bone-grafting>.
- Kinaci A, Neuhaus V, Ring DC. Trends in bone graft use in the United States. *Orthopedics*. 2014;37(9):e783–8.
- Finkemeier CG. Bone-grafting and bone-graft substitutes. *J Bone Joint Surg Am*. 2002;84(3):454–64.
- Pape HC, Evans A, Kobbe P. Autologous bone graft: properties and techniques. *J Orthop Trauma*. 2010;24Suppl 1:S36–40.
- Sarin P, Lee SJ, Apostolov ZD, Kriven WM. Porous biphasic calcium phosphate scaffolds from cuttlefish bone. *J Am Ceram Soc*. 2011;94(8):2362–70.
- Klopung LP, Purwati, Edward M. The healing effect of cuttlefish bone on fractured bone in rat model. *Bali Med J*. 2016;5(2):193–6.
- Irianto KA, Pribadi A, Irsyam IA, Klopung YP, Sindrawati O. A comparison of osteoblast cell proliferation and osteocalcin expression in cuttlefish bone and bovine bone xenograft. *Mol Cell Biomed Sci*. 2019;3(2):75–80.
- Dogan E, Okumus Z. Cuttlebone used as a bone xenograft in bone healing. *Vet Med*. 2014;59(5):254–60.
- Cadman J, Zhou S, Chen Y, Li Q. Cuttlebone: characterisation, application and development of biomimetic materials. *J Bionic Eng*. 2012;9(3):367–76.
- Medical Device and Diagnostic Industry (MDDI) Online. Biocompatibility safety assessment of medical devices: FDA, ISO, and Japanese guidelines [Internet]. 2000 [cited 2018 Sep 18]. Available from: <https://www.mddionline.com/biocompatibility-safety-assessment-medical-devices-fda-iso-and-japanese-guidelines>.
- Mapara M, Thomas BS, Bhat KM. Rabbit as an animal model for experimental research. *Dent Res J*. 2012;9(1):111–8.
- Kilkenny C, Browne WJ, Cuthill IC, Emerson M, Altman DG. Improving bioscience research reporting: the ARRIVE guidelines for reporting animal research. *PLoS Biol*. 2010;8(6):e1000412.
- Rantam FA, P Ferdiansyah, A Purwati. Mesenchymal stem cell, hematopoietic, and application model. Surabaya: Universitas Airlangga; 2014. Indonesian.
- Vajrabhaya LO, Korsuwannawong S, Surarit R. Cytotoxic and the proliferative effect of cuttlefish bone on MC3T3-E1 osteoblast cell line. *Eur J Dent*. 2017;11(4):503–7.
- BioVision Incorporated. Ready-to-use cell proliferation colorimetric reagent, WST-1 (catalog #K304-2500; assays; store at -20°C). Milpitas: BioVision Incorporated.
- Supranto J. Simple random sampling. In: Supranto J. Sampling technique in survey and experimental study. Jakarta: Rineka Cipta; 2000. Indonesian.
- Aminatun, Siswanto, Penga YM, Istifarah, Apsari R. The effect of sintering process on the characteristics of hydroxyapatite from cuttlefish bone (*Sepia Sp.*). *Res J Pharm Biol Chem Sci*. 2013;4(4):1431–42.
- Dimarino AM, Caplan AI, Bonfield TL. Mesenchymal stem cells in tissue repair. *Front Immunol*. 2013;4:201.
- Park JW, Bae SR, Suh JY, Lee DH, Kim SH, Kim H, et al. Evaluation of bone healing with eggshell-derived bone graft substitutes in rat calvaria: a pilot study. *J Biomed Mater Res A*. 2008;87(1):203–14.
- Guda T, Walker JA, Pollot BE, Appleford MR, Oh S, Ong JL, et al. In vivo performance of bilayer hydroxyapatite scaffolds for bone tissue regeneration in the rabbit radius. *J Mater Sci Mater Med*. 2011;22(3):647–56.
- Blokhuis TJ, Arts JJ. Bioactive and osteoinductive bone graft substitutes: definitions, facts and myths. *Injury*. 2011;42 Suppl 2:S26–9.
- Cerrolaza M, Duarte V, Garzón-Alvarado D. Analysis of bone remodeling under piezoelectricity effects using boundary elements. *J Bionic Eng*. 2017;14(4):659–71.
- Checa AG, Cartwright JH, Sánchez-Almazo I, Andrade JP, Ruiz-Raya F. The cuttlefish *Sepia officinalis* (Sepiidae, Cephalopoda) constructs cuttlebone from a liquid-crystal precursor. *Sci Rep*. 2015;5:11513.
- Oryan A, Alidadi S, Moshiri A, Maffulli N. Bone regenerative medicine: classic options, novel strategies, and future directions. *J Orthop Surg Res*. 2014;9(1):18.
- Nazirkar G, Singh S, Dole V, Nikam A. Effortless effort in bone regeneration: a review. *J Int Oral Health*. 2014;6(3):120–4.

Cytotoxic effect of natural cuttlefish bone xenograft: an in vitro and in vivo study

ORIGINALITY REPORT

12%

SIMILARITY INDEX

10%

INTERNET SOURCES

10%

PUBLICATIONS

0%

STUDENT PAPERS

PRIMARY SOURCES

- 1** erepo.unud.ac.id 2%
Internet Source
- 2** Barile, . "Acute dermal toxicity testing", Principles of Toxicology Testing Second Edition, 2013. 1%
Publication
- 3** Aminatun Aminatun, D.E. Fadhilah Handayani, Prihartini Widiyanti, Dwi Winarni, Siswanto Siswanto. "In vivo approach on femur bone regeneration of white rat (*Rattus norvegicus*) with the use of hydroxyapatite from cuttlefish bone (*Sepia spp.*) as bone filler", *Veterinary World*, 2019 1%
Publication
- 4** Jamal Zaini, Agus Dwi Susanto, Erlang Samoedro, Vonni Christiana Bionika, Budhi Antariksa. "Health consequences of thick forest fire smoke to healthy residents in Riau, Indonesia: a cross-sectional study", *Medical Journal of Indonesia*, 2020 1%
Publication

5	Eka Jusup Singka, Innes Ericca. "Hajj health management in Indonesia", Medical Journal of Indonesia, 2020	1 %
Publication		
6	La-ongthong Vajrabhaya, Vanida Sangalungkarn, Ratchapin Srisatjaluk, Suwana Korsuwannawong, Chareerut Phruksaniyom. "Hypochlorite solution for root canal irrigation that lacks a chlorinated odor", European Journal of Dentistry, 2019	1 %
Publication		
7	www.revistafarmacia.ro	1 %
Internet Source		
8	Neng Nenden Mulyaningsih, Ariadne Lakshmidewi Juwono, Djarwani Soeharso Soejoko, Dewi Apri Astuti. "Multi-hole spherical CT scan method to characterize large quantities of bones in rats", Medical Journal of Indonesia, 2021	1 %
Publication		
9	www.dovepress.com	<1 %
Internet Source		
10	1library.net	<1 %
Internet Source		
11	Nurfairuz Abdul Malik, Mahaneem Mohamed, Mohd Zulkifli Mustafa, Azalina Zainuddin. "In vitro modulation of extracellular matrix genes	<1 %

by stingless bee honey in cellular aging of human dermal fibroblast cells", Journal of Food Biochemistry, 2019

Publication

12	balimedicaljournal.org Internet Source	<1 %
13	synapse.koreamed.org Internet Source	<1 %
14	www.ncbi.nlm.nih.gov Internet Source	<1 %
15	bio-protocol.org Internet Source	<1 %
16	www.degruyter.com Internet Source	<1 %
17	www.freepatentsonline.com Internet Source	<1 %
18	academic.oup.com Internet Source	<1 %
19	escholarship.org Internet Source	<1 %
20	www.karger.com Internet Source	<1 %
21	www.phcog.com Internet Source	<1 %

Exclude quotes Off

Exclude matches < 10 words

Exclude bibliography On

Cytotoxic effect of natural cuttlefish bone xenograft: an in vitro and in vivo study

GRADEMARK REPORT

FINAL GRADE

/0

GENERAL COMMENTS

Instructor

PAGE 1

PAGE 2

PAGE 3

PAGE 4

PAGE 5

PAGE 6

PAGE 7
