

Valsartan

by Suciati Suciati

Submission date: 18-Dec-2019 12:07PM (UTC+0800)

Submission ID: 1236356411

File name: C-4_extract.pdf (1.32M)

Word count: 18827

Character count: 91363



Valsartan

Febry Ardiana*, Suciati[†], Gunawan Indrayanto^{†,1}

*Research & Development, Bemofarm Pharmaceutical Company, Buduran-Sidoarjo, Indonesia

[†]Faculty of Pharmacy, Airlangga University, Dharmawangsa Dalam, Surabaya, Indonesia

¹Corresponding author: e-mail address: gunawanindrayanto@yahoo.com

Contents

1. General Information	432
1.1 Solubility	432
1.2 Chemical Name	432
1.3 Synonym	432
1.4 CAS Number	433
1.5 pK_a and Partition Coefficient	433
1.6 Proprietary Preparations	433
1.7 Structural Formula	433
1.8 Molecular Formula, Molecular Weight, and Elementary Analysis	433
2. Thermal Analysis	433
3. Polymorph	435
4. Related Compounds and Impurities	436
5. Spectroscopy	440
5.1 Ultraviolet Spectroscopy	440
5.2 Infrared Spectroscopy	446
5.3 Mass Spectroscopy	446
5.4 Nuclear Magnetic Spectroscopy	449
6. Compendial Method	450
6.1 Identification Test	450
6.2 Related Compounds	453
6.3 Assay	454
6.4 Dissolution Test	454
7. Method of Analysis	455
7.1 Spectrophotometer	455
7.2 Capillary Electrophoresis	456
7.3 High-Performance Liquid Chromatography	457
7.4 Thin-Layer Chromatography	457
8. Analysis of Valsartan in Biological Fluids and Herbal Preparation	469
Acknowledgments	486
References	486



1. GENERAL INFORMATION

Valsartan is an antihypertensive drug which selectively inhibits angiotensin receptor type II [1,2]. This tetrazole derivative was first developed by Novartis and marketed under brand name Diovan[®] [3]. This compound is orally active and is rapidly absorbed after oral doses, having a bioavailability of approximately 23% [4].

Valsartan appears as a white or almost white hygroscopic powder [5]. This compound must be kept in an air-tight container and should be protected from light and heat [6].

It is available in film-coated tablets containing valsartan 40, 80, 160, or 320 mg, and capsules with dosage of 80 or 160 mg. Tablet combinations of valsartan with hydrochlorothiazide or amlodipine are also available [3,7].

1.1 Solubility

Valsartan is practically insoluble in water, sparingly soluble in methylene chloride, and freely soluble in anhydrous ethanol and methanol [3,5]. Solubility of valsartan in water is affected by pH; solubility in pH 4.07, 7.02, 9.18, and 10.06 are 56.6, 62.8, 71.6, and 100.0 mg/100 mL, respectively [8].

1.2 Chemical Name

N-(1-Oxopentyl)-*N*-[[2'-(1*H*-tetrazol-5-yl)[1,1'-biphenyl]-4-yl]methyl]-*L*-valine

L-Valine, *N*-(1-oxopentyl)-*N*-[[2'-(1*H*-tetrazol-5-yl)[1,1'-biphenyl]-4-yl]methyl]

p-(*o*-1*H*-Tetrazol-5-ylphenyl)benzyl]-*N*-valeryl-*L*-valine

(2*S*)-3-Methyl-2-[pentanoyl[[2'-(1*H*-tetrazol-5-yl)biphenyl-4-yl]methyl]amino]butanoic acid

3-Methyl-2-[pentanoyl-[[4-[2-(2*H*-tetrazol-5-yl)phenyl]phenyl]methyl]amino]butanoic acid

N-pentanoyl-*N*-[[2'-(1*H*-tetrazol-5-yl)biphenyl-4-yl]methyl]-*L*-valine

(*S*)-*N*-(1-Carboxy-2-methylprop-1-yl)-*N*-pentanoyl-*N*-[2'-(1*H*-tetrazol-5-yl)-biphenyl-4-ylmethyl]amine [2,4-6,9].

1.3 Synonym

CGP-48993, Diovan[®], and Tareg[®] [1,4,9]

1.4 CAS Number

137862-53-4 [4–6]


1.5 pK_a and Partition Coefficient

The pK_a of valsartan is 4.9 [10]; partition coefficient in octanol–water system is 22.2; *n*-octanol–phosphate buffer is 0.033 [8,9]

1.6 Proprietary Preparations

Diovan[®], Co-Diovan[®] (combination with hydrochlorothiazide), and Exforge (combination with amlodipine) [3,4]

1.7 Structural Formula

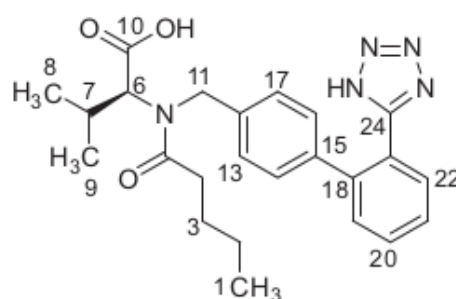

1.8 Molecular Formula, Molecular Weight, and Elementary Analysis

Molecular formula: $C_{24}H_{29}N_5O_3$ [4–6]

Molecular weight: 435.52 [6]

Elementary analysis: HRMS (calculated): 436.23432 $[M + H]^+$ [11].

C: 66.19%, H: 6.71%, N: 16.08%, and O: 11.02% [9].



2. THERMAL ANALYSIS

The melting point of valsartan crystallized from diisopropyl ether is 116–117 °C [9]. Unfortunately, the melting point measurement of the USP valsartan RS (lot L0L195) and valsartan raw material from TEVA, India (Control No. 761120512; using a Mettler Toledo MP 70 system) did not show sharp light intensity curves (see Figure 1), so their melting points cannot be determined accurately.

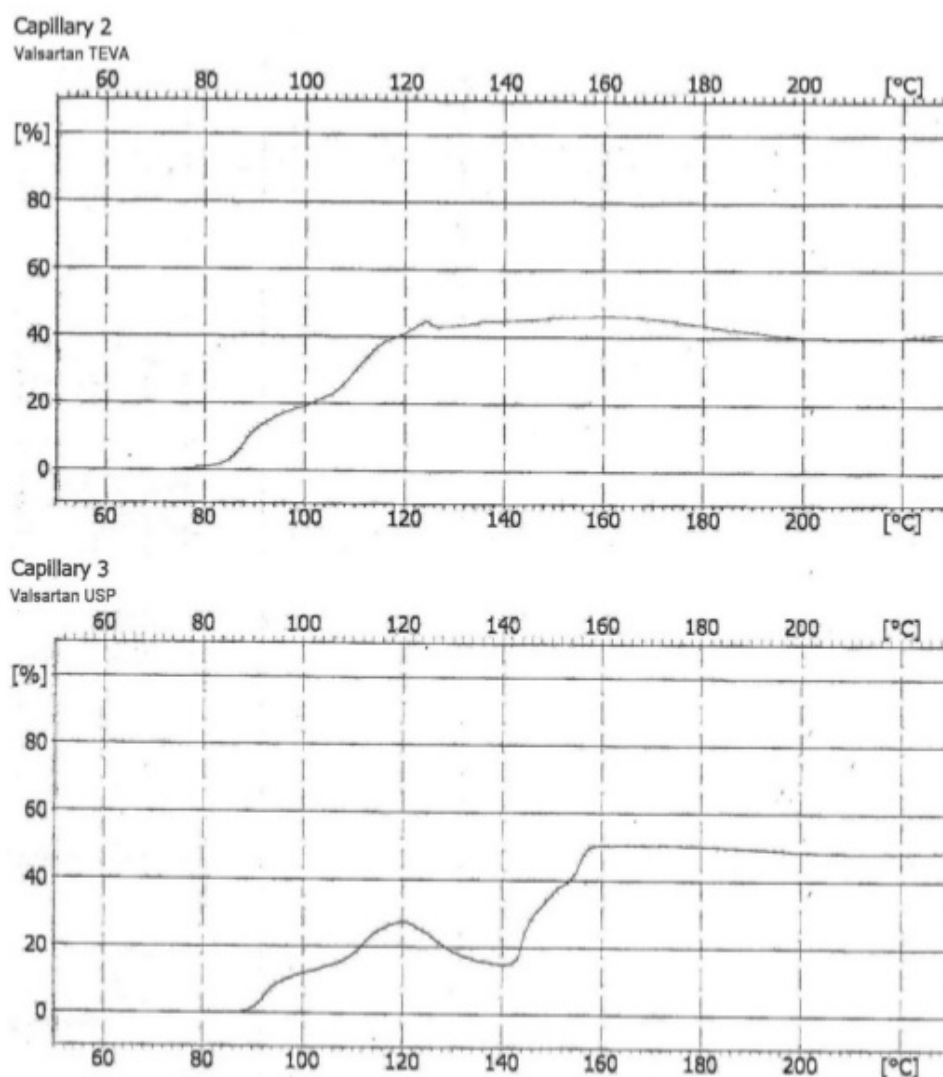


Figure 1 Light intensity curves of USP valsartan RS and valsartan bulk drug (TEVA, India) recorded on a Mettler Toledo MP 70. Y-axis: light intensity (%); X-axis: temperature (°C).

Nalluri *et al.* [12] reported that recrystallization of valsartan from methanol, ethanol, isopropanol, and acetonitrile could yield different endothermic (broad) peaks according to their DSC thermograms (75.2, 65.9, 64.6, and 80.6 °C, respectively). These results indicated that crystallization of valsartan from different solvents could yield different crystal forms. The broad endothermic peak of valsartan was reported at 102.8 °C [13] and 100.8 °C [14], while a sharp endothermic peak was reported at around 118 °C [15], 115 °C [16], and 115.77 °C [17]. Skotnicki *et al.* [18] reported that valsartan gave two endothermic peaks at approximately 80 °C (a small

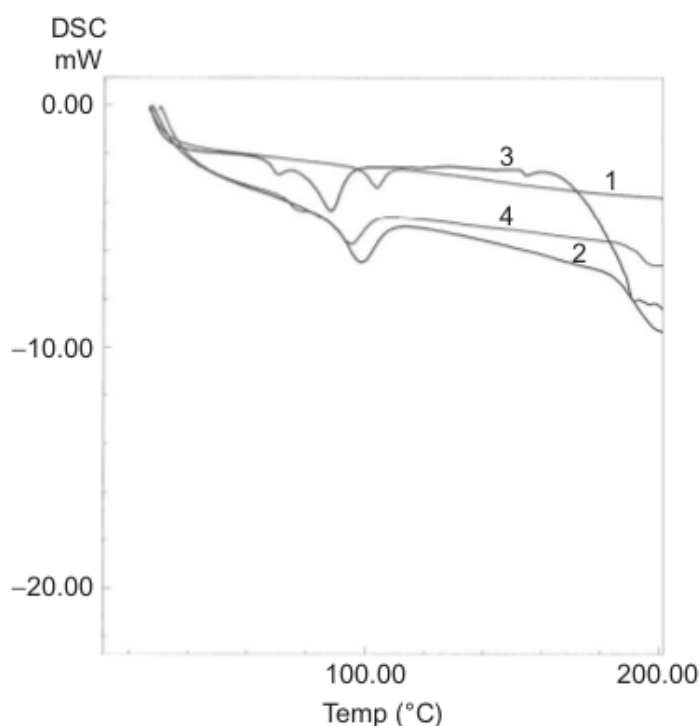


Figure 2 DSC curves of USP valsartan RS and valsartan bulk drug (TEVA, India) recorded on a DSC-60A Shimadzu Thermal Analyzer. (1) USP valsartan RS rate 20 °C/min, (2) USP valsartan RS rate 40 °C/min, (3) valsartan TEVA rate 20 °C/min, and (4) valsartan TEVA rate 40 °C/min.

broad peak) and at 100.6 °C (a sharp peak). These peaks are related to water evaporation and enthalpy relaxation. The glass transition of valsartan appeared at around 76 °C. [Figure 2](#) shows the DSC thermograms of USP valsartan RS and valsartan TEVA measured by using a DSC-60A Shimadzu Thermal Analyzer. The DSC of valsartan TEVA exhibits two endothermic peaks (2 and 3), which are similar to that of DSC reported by Skotnicki *et al.* [18]. Based on our data and several published results, it can be suggested that crystal forms of valsartan available in the market are not identical.



3. POLYMORPH

Valsartan may occur in 12 different polymorphs and an amorphous form. X-Ray data of various crystal forms are presented in [Table 1](#).

Burgbacher *et al.* [21] reported the invention of a highly crystalline form of valsartan. The crystal was characterized by its XRPD pattern that

Table 1 X-Ray Powder Diffraction of Valsartan Polymorphs

Polymorph	Scattering Angle (Degree 2θ)	Reference
Form I	(5.4; 13.0; 16.3; 19.5; 20.7; 23.4) $\pm 0.2^\circ$	[19]
Form II	(5.8; 12.7; 14.0; 17.6; 20.8; 22.5) $\pm 0.2^\circ$	[19]
Form III	(5.1; 10.1; 15.3; 18.6) $\pm 0.2^\circ$	[19]
Form IV	(6.2; 10.7; 14.5; 15.7; 19.0; 23.5; 24.8) $\pm 0.2^\circ$	[19]
Form VI	(5.5; 13.3; 14.3; 17.7; 21.1; 22.3) $\pm 0.2^\circ$	[19]
Form VII	(5.2; 15.2; 15.9; 18.6; 22.8; 23.6) $\pm 0.2^\circ$	[19]
Form VIII	(5.7; 13.6; 18.0)	[19]
Form IX	(6.3; 14.0; 17.9) $\pm 0.2^\circ$	[19]
Form X	(5.6) $\pm 0.2^\circ$ broad peaks: (15.0; 20.6) $\pm 0.2^\circ$	[19]
Form XI	(5.2; 10.5; 12.9; 13.9; 18.8) $\pm 0.2^\circ$	[19]
Form XIII	(5.1; 11.6; 15.8; 18.6; 26.2) $\pm 0.2^\circ$	[19]
Amorphous	Very broad peak	[20]
Valsartan highly crystalline form	(9.308; 10.74; 11.643; 13.854; 15.136; 16.056; 16.686; 17.643; 18.561; 19.186; 20.024; 20.567; 21.335; 21.595; 21.858; 22.879; 24.597; 25.051; 26.292; 31.032) $\pm 0.2^\circ$	[21]

had a scattering angle peak ($2^\circ\theta$) about $(31.0 \pm 10.2)^\circ$; its melting point was $(140.8 \pm 3)^\circ\text{C}$. The invention stated that this highly crystalline form of valsartan can be easily dried compared to other forms of valsartan. This crystal has less residual solvent, higher stability, and better purity compared to the other forms of valsartan. The highly crystalline form of valsartan can also be characterized with SEM that showed very low water content. The quasi flower-like conglomerates and the formation of spheroid conglomerates were assumed for high flowability of the highly crystalline valsartan.



4. RELATED COMPOUNDS AND IMPURITIES

Valsartan has one chiral center in the valine moiety, and usually the pure (*S*)-enantiomer is used for pharmaceutical purposes. This *S* enantiomer

is derived from the L-valine that is used in the synthesis [22]. Method for testing enantiomer impurity is described in BP 2013 [5].

There are three related compounds or impurities described in the USP 36 [6] and BP 2013 [5]. These compounds are impurity 15 A, B, and C in BP 2013 [5] which are identical to USP valsartan-related compound A RS, USP valsartan-related compound B RS, and USP valsartan-related compound C RS.

Nie *et al.* [23] isolated three impurities from crude valsartan, and one of the impurities was identified as impurity A. This compound may be formed by hydrolysis 30 of a methyl ester intermediate in the last step of valsartan synthesis; the methyl group can react with the nitrogen of the tetrazole acid isostere to form impurity A. Approximately 80% of valsartan could be recovered by dissolving impurity A in methanol (pH 11) over 24 h. Impurity A described by Nie *et al.* [23] is different from impurity A in BP 2013 [5].

Sampath *et al.* [24] detected five impurities, i.e., impurities I, II, III, IV, and V from crude drug valsartan. The characterization and structure elucidation of the impurities were achieved by spectroscopic methods (IR, NMR, and MS) as well as by synthesis. The mechanism of impurities formation has been described in detail.

Mehta *et al.* [25] reported three impurities (DP-1, DP-2, and DP-3) from the results of degradation studies using acidic (in 1 N HCl), alkaline (in 2 N NaOH), neutral (in water), oxidative (in 30% H₂O₂), and photolytic studies (exposure to UV radiation). The structure of DP-1 reported by Mehta and coworkers [25] was identical to that of impurity I described by Sampath *et al.* [24].

Bianchini *et al.* [26] reported two degradation products of valsartan (DP-1 and DP-2) from exposure of valsartan to UV-VIS radiation (320 nm). DP-1 was formed by light-induced decarboxylation of valsartan, while DP-2 could be formed by further decomposition 30 of tetrazole moiety from the loss of nitrogen and cyclization. The chemical structures of DP-1 and DP-2 are not identical to those of DP-1 and DP-2 reported by Mehta *et al.* [25].

Table 2 summarizes known impurities and related substances of valsartan, and their chemical structures are presented in Figure 3. The compendial method of separation and analysis of the impurities are summarized in Table 3. Raw pharmaceutical substances are required to have related substances or impurities less than 0.05% [24].

Table 2 Impurities and Related Compounds of Valsartan

Compound Number	Chemical Name	Other Name	Chemical Formula	Degradation Process or Formation	Ref.
1	22-N-Valeryl-N-{[2'-(1H-tetrazole-5-yl)biphenyl-4-yl]-methyl} valine	Compound A RS, impurity A	C ₂₄ H ₂₉ N ₅ O ₃		[5,6]
2	22-N-Butyryl-N-{[2'-(1H-tetrazole-5-yl)biphenyl-4-yl]-methyl} valine	Compound B RS, impurity C, CGP 55390	C ₂₃ H ₂₇ N ₅ O ₃		[5,6,27]
3	7-N-Valeryl-N-{[2'-(1H-tetrazole-5-yl)biphenyl-4-yl]-methyl} valine benzyl ester	Compound C RS, impurity B	C ₃₁ H ₃₅ N ₅ O ₃		[5,6]
4	(S)-N-Valeryl-N-{[2'-(1-methyl-tetrazole-5-yl)biphenyl-4-yl]-methyl} valine	Impurity A ^a	C ₂₅ H ₃₁ N ₅ O ₃	1 Hydrolysis of the methyl ester in the last step of crude synthesis	[23]
5	8 (S)-N-(1-Carboxy-2-methylprop-1-yl)-N-[2'-(1H-tetrazol-5-yl)-biphenyl-4ylmethyl]amine	Impurity I, DP-1 ^b	C ₁₉ H ₂₁ N ₅ O ₂	Intermediate product of valsartan Simple amide hydrolysis under acidic and neutral conditions	[24,25]
6	8 (S)-N-(1-Carboxy-2-methylprop-1-yl)-N-(5-phenylthio)pentanoyl-N-[2'-(1H-tetrazol-5-yl)-biphenyl-4ylmethyl]amine	Impurity II	C ₃₀ H ₃₃ N ₅ O ₃ S	Intermediate product of valsartan	[24]

7	⁸ (S)-N-(1-Carboxy-2-methylprop-1-yl)-N-(5-phenyl)pentanoyl-N-[2'-(1H-tetrazol-5-yl)-biphenyl-4ylmethyl]amine	⁸⁶ Impurity III	C ₃₀ H ₃₃ N ₅ O ₃	Condensation of impurity I with 5-phenylvaleroyl chloride in the presence of pyridin	[24]
8	(S)-N-(1-Carboxy-2-methylprop-1-yl)-N-4-pentenoyl-N-[2'-(1H-tetrazol-5-yl)-biphenyl-4ylmethyl]amine	Impurity IV	C ₂₄ H ₂₇ N ₅ O ₃	Condensation of impurity I with 4-pentanoic acid	[24]
9	(S)-N-(1-Carboxy-2-methylprop-1-yl)-N-(5-hydroxy)pentanoyl-N-[2'-(1H-tetrazol-5-yl)-biphenyl-4ylmethyl]amine	Impurity V	C ₂₄ H ₂₉ N ₅ O ₄	Condensation of impurity I with 5-chlorovaleric acid	[24]
10	⁶⁰ n/a	DP-2 ^b	C ₁₉ H ₂₁ N ₅ O ₃	N-Dealkylation in photoneutral condition	[25]
11	n/a	DP-3 ^b	C ₂₄ H ₂₈ N ₅ O	Cyclization of tetrazole ring in photo acidic condition	[25]
12	N-[2'-(1H-Tetrazol-5-yl)-biphenyl-4ylmethyl]-N-isobutylpentanamide	DP-1	C ₂₃ H ₃₀ N ₅ O [M + H] ⁺	Photodegradation of valsartan	[26]
13	⁴¹ N-(Diazirinol[1,3-f]phenanthridin-4ylmethyl)-N-isobutylpentanamide	DP-2	C ₂₃ H ₂₇ N ₃ O	Photodegradation of valsartan	[26]

^aDegradation product from Nie *et al.* [23].^bDegradation products from Mehta *et al.* [25].

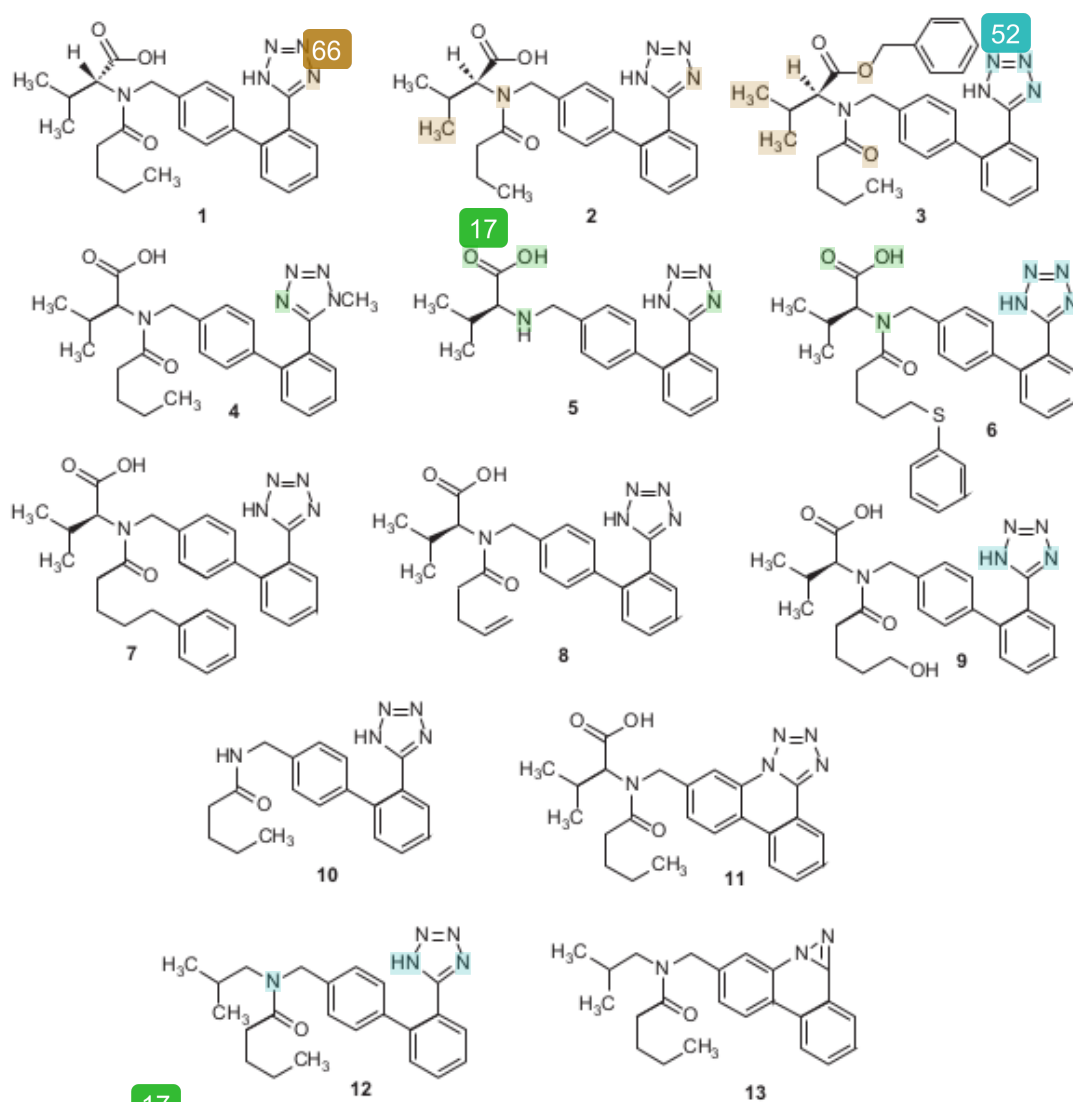


Figure 3 Chemical structures of impurities and related compound of valsartan. Compound numbers were described in Table 2.

5. SPECTROSCOPY

5.1 Ultraviolet Spectroscopy

The UV spectra of valsartan in methanol (10 ppm, scanned 200–400 nm) are presented in Figure 4. These spectra were recorded using a Spectrophotometer Shimadzu UV-1800 series; valsartan showed emission wavelength at 399 nm (excitation 259 nm), and the intrinsic fluorescence was decreased in basic pH [32].

Table 3 Compendial LC Methods and Their Validation Report

Analyte(s)	Sample	Test	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	Ref.
Valsartan	Raw material	Compendial USP assay method	Water:ACN:glacial acetic acid (500:500:1, v/v/v) ²⁷	Column: L1 (125 mm × 3 mm; 5 μm) Temperature: n/a Detector: 273 nm	Water:ACN:glacial acetic acid (500:500:1, v/v/v) ⁴ Flow rate: 0.4 mL/min Injection volume: 5 μL	[6]
Valsartan	Raw material	Assay method in USP/MC validation report	Water:ACN:glacial acetic acid (500:500:1, v/v/v) ⁹	Column: L10 Zorbax SB-CN (250 mm × 4.6 mm; 5 μm) ⁶ Temperature: 25 °C Detector: 273 nm	Water:ACN:glacial acetic acid (500:500:1, v/v/v) ⁴ Flow rate: 0.7 mL/min Injection volume: 10 μL	[28]
Valsartan-related compound A	Raw material	Compendial USP method for related compound A	<i>n</i> -Hexane:2-propanol:trifluoroacetic acid (85:15:0.1, v/v/v)	Column: L40 (250 mm × 4.6 mm; 5 μm) ⁴⁰ Temperature: n/a Detector: 230 nm	<i>n</i> -Hexane:2-propanol:trifluoroacetic acid (85:15:0.1, v/v/v) ⁴⁷ Flow rate: 0.8 mL/min Injection volume: 10 μL	[6]
Valsartan-related compounds B, C, and other related compounds	Raw material	Compendial USP method for related compounds B, C, and other related compounds	Water:ACN:glacial acetic acid (500:500:1, v/v/v)	Column: L1 (125 mm × 3 mm; 5 μm) Temperature: n/a Detector: 225 nm	Water:ACN:glacial acetic acid (500:500:1, v/v/v) ⁴ Flow rate: 0.4 mL/min Injection volume: 5 μL	[6]

Continued

Table 3 Compendial LC Methods and Their Validation Report—cont'd

Analyte(s)	Sample	Test	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	Ref.
Valsartan	Raw material	Organic impurities test in USP/MC validation report	Water:ACN:acetic glacial acid (500:500:1, v/v/v)	35 μ m: L10 Zorbax SB-CN (250 mm \times 4.6 mm; 5 μ m) Column oven temp: 25 $^{\circ}$ C Auto sampler temp: 4 $^{\circ}$ C Detector: 225 nm MS condition: Detector: MS Source: ES Scan (+ and -) Capillary (kv): 3.00 Cone (v): 20.0 Extractor (v): 2.0 44 lens (v): 0.1 Source temperature: 80 $^{\circ}$ C Desolvation temperature: 400 $^{\circ}$ C	Water:ACN:26 acetic glacial acid (500:500:1, v/v/v) Flow rate: 0.7 mL/min Injection volume: 10 μ L	[28]
Valsartan	Tablet, organic impurities	Compendial USP assay method	7 N:water (1:1, v/v)	Column: L1 (250 mm \times 4.6 mm; 11 μ m) Column temperature: 30 $^{\circ}$ C Detector: 230 nm	Water:ACN:glacial acetic acid (50:50:1, v/v/v) Flow rate: 1.0 mL/min Injection volume: 20 μ L	[6]
Valsartan	Tablet	Assay method in USP/MC validation report	9 Water:ACN:glacial acetic acid (500:500:1, v/v/v)	59 Column: L10 Zorbax SB-CN (250 mm \times 4.6 mm; 5 μ m) Column oven temp: 25 $^{\circ}$ C Auto sampler temp: 4 $^{\circ}$ C Detector: 225 nm	Water:ACN:4 glacial acetic acid (500:500:1, v/v/v) Flow rate: 0.7 mL/min Injection volume: 10 μ L	[29]

Valsartan	Tablet	Organic impurities test in USP/MC validation report	Water:ACN: glacial acetic acid (500:500:1, v/v/v)	<p>35 μm: L10 Zorbax SB-CN (250 mm × 4.6 mm; 5 μm) Column oven temp: 25 °C Auto sampler temp: 4 °C Detector: 225 nm MS condition: Detector: MS Source: ES Scan (+ and -) MS condition: Capillary (kv): 3.00 Cone (v): 20.0 Extractor (v): 2.0 44 lens (v): 0.1 Source temperature: 80 °C Desolvation temperature: 400 °C</p>	<p>Water:ACN:26 glacial acetic acid [29] (500:500:1, v/v/v) Flow rate: 0.7 mL/min Injection volume: 10 μL</p>
Valsartan, hydrochlorothiazide	Tablet, organic impurities	Compendial USP assay method	ACN:water (1:1, v/v)	<p>Column: L1 (125 mm × 3.0 mm; 5 μm) Column temperature: n/a Detector: 265 nm</p>	<p>Solution A: ACN:water: [6] trifluoroacetic acid (10:90:0.1, v/v/v) Solution B: ACN:water: 4 trifluoroacetic acid (90:10:0.1, v/v/v) (gradient) Flow rate: 0.4 mL/min Injection volume: 10 μL</p>

Continued

Table 3 Compendial LC Methods and Their Validation Report—cont'd

Analyte(s)	Sample	Test	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	Ref.
Valsartan, hydrochlorothiazide	Tablet	Assay method in USP/MC validation report	Solution A: 0.2% acetic acid in water Solution B: ACN Solvent: solution A: solution B (1:1, v/v)	Column: L10 Zorbax SB-CN (250 mm × 4.6 mm; 5 μm) Column oven temp: 25 °C Auto sampler temp: 10 °C Detector: 272 nm	Solution A: 0.2% acetic acid in water Solution B: ACN (gradient) Flow rate: 1 mL/min Injection volume: 10 μL	[30]
Valsartan, hydrochlorothiazide	Tablet	Organic impurities test in USP/MC validation report	Solution A: 0.2% acetic acid in water Solution B: ACN Solvent: solution A: solution B (1:1, v/v)	Column: L10 Zorbax SB-CN (250 mm × 4.6 mm; 5 μm) Column oven temp: 25 °C Auto sampler temp: 10 °C Detector: 272 nm MS condition: Detector: MS Source: ES Scan (+ and -) MS condition: Capillary (kv): 3.00 Cone (v): 25.0 Extractor (v): 2.0 RF lens (v): 0.1 Source temperature: 90 °C Desolvation temperature: 500 °C	Solution A: 0.2% acetic acid in water Solution B: ACN (gradient) Flow rate: 1 mL/min Injection volume: 10 μL	[30]

Valsartan, amlodipine	Tablet	Assay method in USP/MC validation report	Water:MeOH (3:7, v/v)	Column: 11 L1 Purosphere Star RP18 (150 mm × 4.6 mm; 5 μm) Column oven temp: 40 °C Auto sam 10 temp: 10 °C Detector: 237 nm	Solution A: 25 mM ammonium acetate pH 7.5 Solution B: ACN (gradient) Flow rate: 1 mL/min Injection volume: 20 μL	[31]
Valsartan	Raw material	Organic impurities test in USP/MC validation report	Water:MeOH (3:7, v/v)	Column: L 13 Purosphere Star RP18 (150 mm × 4.6 mm; 5 μm) Column oven temp: 40 °C Auto sampler temp: 10 °C Detector: 237 nm MS condition: Detector: MS Source: ES Scan (+ and -) MS condition: Capillary (kv): 3.00 Cone (v): 25.0 Source temperature: 90 °C Desolvation temperature: 500 °C	Solution A: 25 mM ammonium acetate pH 7.5 12 Solution B: ACN (gradient) Flow rate: 1 mL/min Injection volume: 20 μL	[31]
Valsartan	Raw material	Compendial BP method of enantiomer purity	Trifluoroacetic acid:2- propanol: hexane (0.1:15:85, v/v/v)	Column: silica gel OD for chiral separations (250 mm × 4.6 mm) Temperature: n/a Detector: 230 nm	Trifluoroacetic acid:2- 18 propanol:hexane (0.1:15:85, v/v/v) Flow rate: 0.8 mL/min Injection volume: 5 μL	[5]
Valsartan	Raw material	Compendial BP method for related substances	Glacial acetic acid:ACN: water (1:500:500, v/v/v)	Column: end-capped 85 decylsilyl silica gel (125 mm × 3.0 mm; 5 μm) Temperature: n/a Detector: 230 nm	17 Glacial acetic acid:ACN:water (1:500:500, v/v/v) Flow rate: 0.4 mL/min Injection volume: 10 μL	[5]

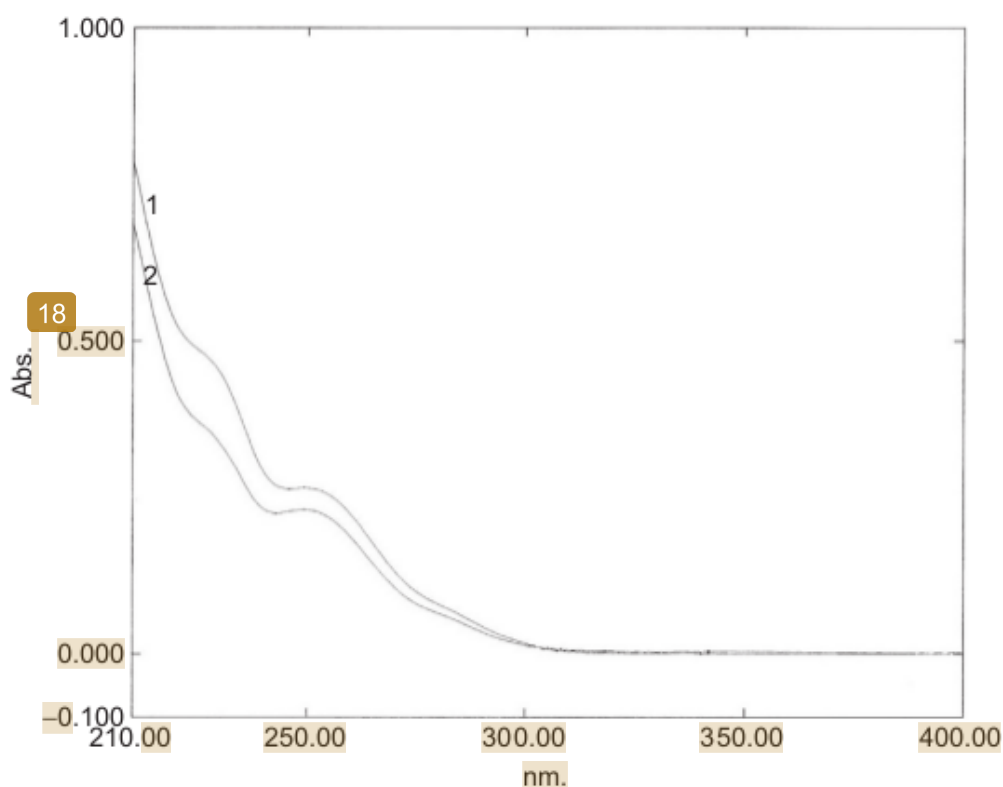


Figure 4 UV spectra of valsartan in methanol (10 ppm) recorded on a Shimadzu UV-1800 series spectrophotometer. (1) USP valsartan RS and (2) valsartan (TEVA, India).

5.2 Infrared Spectroscopy

The IR spectra of USP valsartan RS in KBr (2 mg/200 mg) pellets are shown in Figure 5. IR spectra were obtained by using a Shimadzu IR Prestige-21 Spectrophotometer. The principle peaks and their assignments are listed in Table 4.

Several authors used IR spectrum to characterize nanoparticles of valsartan with Eudragit[®] L 100 [33], solid dispersion of valsartan and cyclodextrin [34,35], and hydrogel formulations [13]. Data from these references suggested that the IR spectrum of valsartan bulk drug is not fully identical, which confirmed the results of DSC (see Section 2).

5.3 Mass Spectroscopy

The ESI-MS of valsartan (TEVA) is shown in Figure 6A. The ESI-M-S (positive ion mode) is measured by using a UPLC Waters H Class and MS: XEVO TQD (direct injection to MS; cone voltage 30 V), while MS/MS (collision; 20 V) is presented in Figure 6B. The important

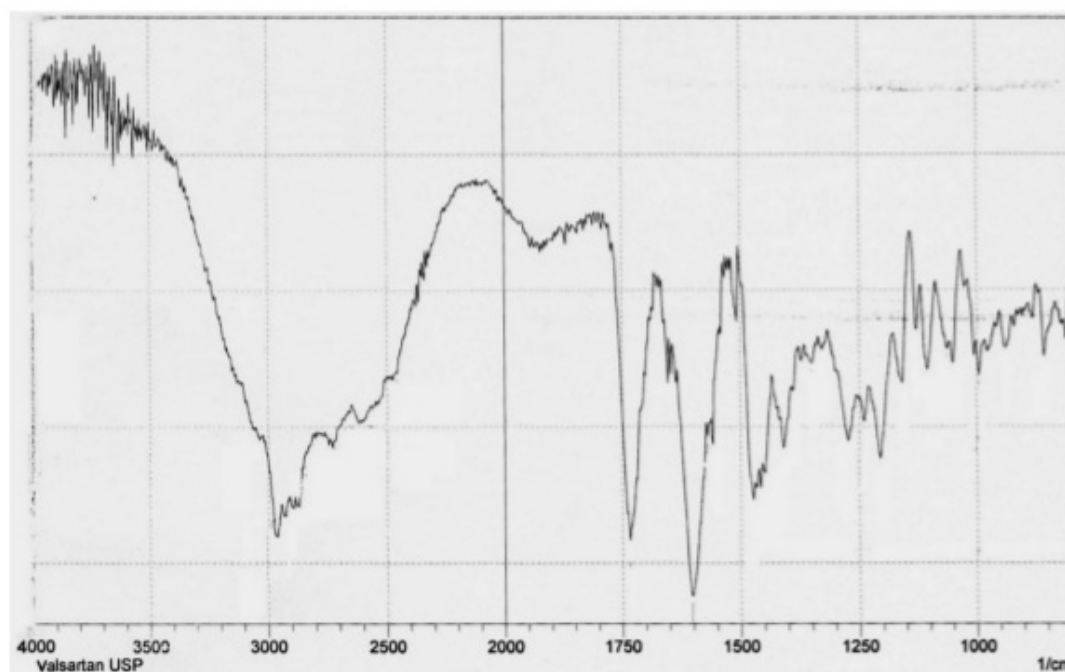


Figure 5 IR spectrum of USP valsartan RS recorded on a Shimadzu IR Prestige-21 spectrophotometer using KBr pellet sampling (2 mg/200 mg).

Table 4 IR Characteristic Band of USP Valsartan RS Lot L0L195^a

Frequency (cm ⁻¹)	Assignment
3035.96	O—H stretch
2873.94	C—H stretch
2927.94	C—H stretch
1734.01	C=O stretch, acid
1602.85	C=O stretch, amide
1570.06	N—N bending
1473.62	C—OH plane bend
1448.54	C—OH plane band
1390.68	C—O stretch
995.27	C—N stretch

^aAssignments were compared to previous published data [13,24].

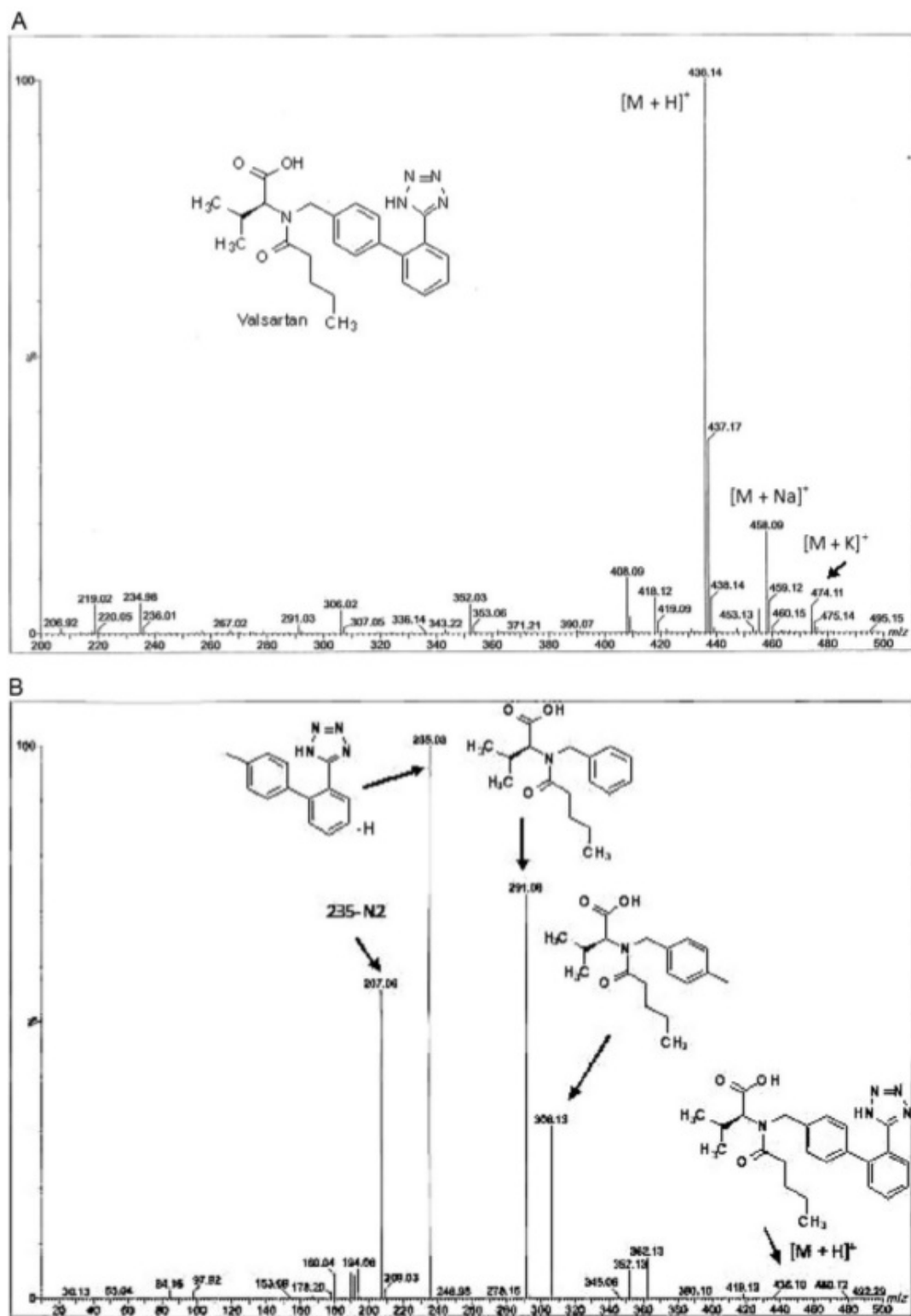


Figure 6 (A) ESI-MS spectrum of valsartan (TEVA) recorded on a UPLC Waters H Class and MS: XEVO TQD (direct injection to MS; cone voltage 30 V). (B) MS/MS spectrum of valsartan (precursor ion m/z 436, collision energy 20 V).

fragments are also described. The MS data are similar to that of described in the literature [36–39]. The ESI negative ion mode is also reported; the important fragments (m/z) are 434.25 ($M-H^-$), 390.1, 350.1, and 179.2 [40,41]. Both positive and negative ion modes are used for determination of valsartan by LC–MS/MS (see Section 8).

The EI-MS of valsartan dimethyl ester was reported by Maurer *et al.* [42]. The important fragments (m/z) were 463 ($[M]^+$, valsartan dimethyl ester), 378 [$M-57$], 320 [$M-115$], 268 [$M-115-57$], and 250 [$M-115-78$].

5.4 Nuclear Magnetic Spectroscopy

Li *et al.* [43] reported two conformations for valsartan based on the result of 2D NOESY experiment and quantum chemistry calculations. The two conformers, namely conformers A (major) and B (minor), are presented in Figure 7 (ChemBio3D Ultra 13, energy minimization with MM2 software to an RMS gradient of 0.100). These *cis/trans* isomer existed as the result of the conformational interchange via rotation around the peptide C(O)—N bond. Potamitis and coworkers [44] as well as Li *et al.* [43] suggested that the two conformers can be distinguished by NOESY experiment. A NOESY experiment was conducted on valsartan sample. In the major conformer A, NOESY cross peaks were observed from H-2, H-3, and H-4 to the methylene protons at C-11 (Figures 7 and 8). Other NOESY correlation observed was between H-4 and the aromatic proton H-13/H-17.

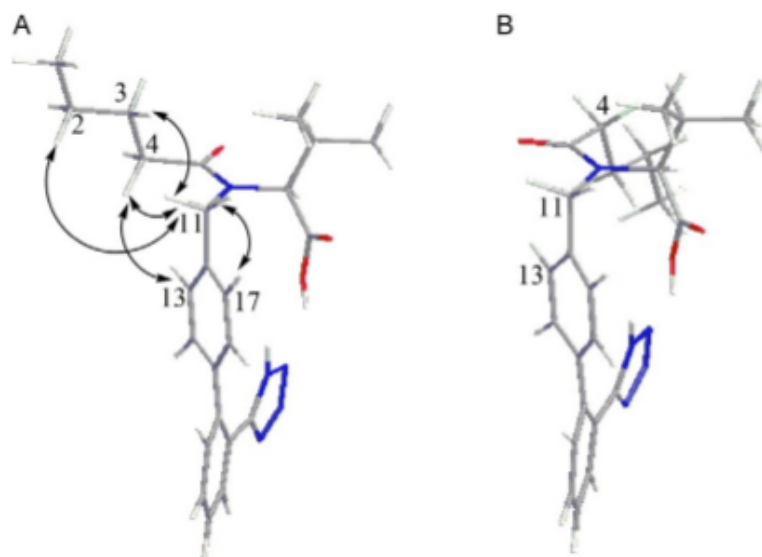


Figure 7 Conformers of valsartan and their NOESY correlation. Structures were drawn by ChemBio3D Ultra 13, energy minimization with MM2 software to an RMS gradient of 0.100: (A) major and (B) minor.

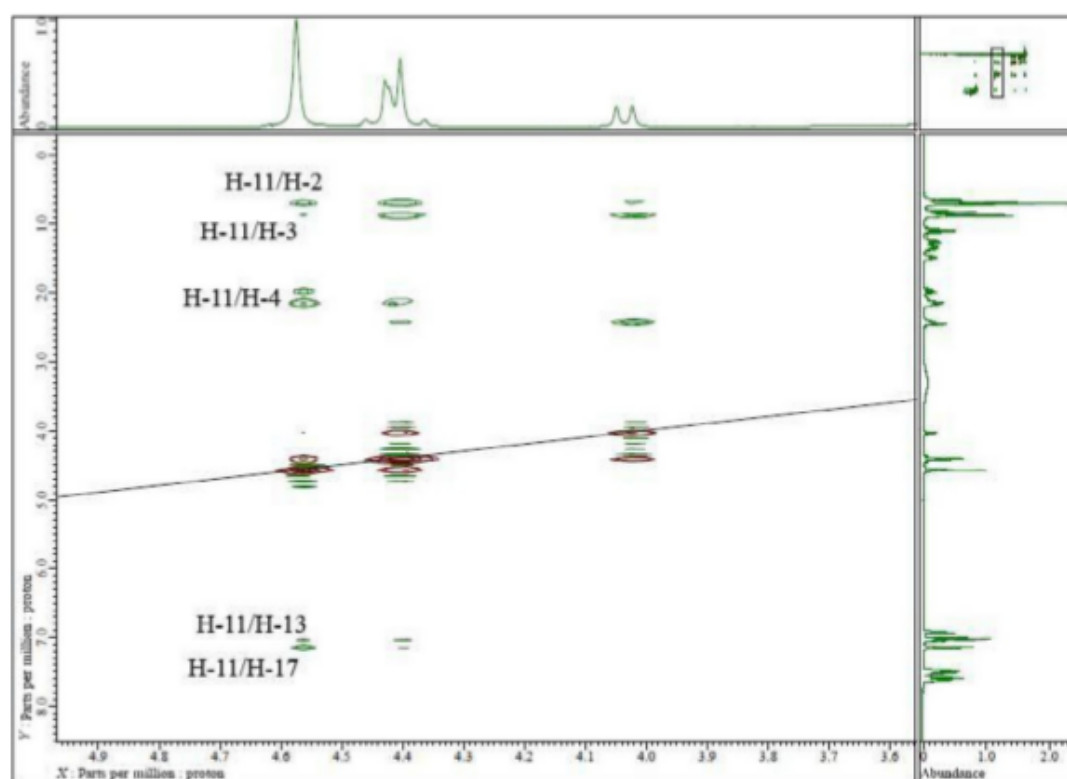


Figure 8 Selected 2D NOESY spectrum for valsartan (TEVA, India) recorded on a JEOL ECS-400, NMR spectrometer in DMSO- d_6 at 294 K with mixing time 500 ms.

These data indicated the closeness of the aliphatic *n*-butyl side chain to the aromatic ring. In the minor conformer, these NOESY correlations were absent. The composition of the two conformers was calculated at different temperature ranging from 293 to 313 K by using ^1H NMR spectroscopy. The data suggested that population of the minor conformer B increased as the temperature increased. Potamitis *et al.* [44] suggested that between 50 and 65 °C the ratio of isomer A:B is $\sim 60:40$. The ^1H and ^{13}C NMR spectra also indicated two sets of data which are presented in Table 5.



6. COMPENDIAL METHOD

6.1 Identification Test

USP 36 [6] described two methods for identification of valsartan bulk drug, i.e., IR absorption and LC system (by comparison of the retention time). This LC method has also been applied for identification of valsartan in tablets (valsartan tablets, valsartan, and hydrochlorothiazide tablets;

Table 5 NMR Data of Two Conformers of Valsartan

#	Major Conformer						Minor Conformer			
	δ_C			δ_H			δ_C		δ_H	
	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆ ^b	DMSO- <i>d</i> ₆ ^c	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆ ^b	DMSO- <i>d</i> ₆ ^c	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆ ^b	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆ ^b
1	14.5	13.6	14.1	0.84 (<i>t</i> , 7.4)	0.75	0.74 (<i>t</i> , 7.2)	14.5	13.6	0.95 (<i>t</i> , 7.4)	0.88
2	23.7	21.6	22.1	1.24	1.15	1.10–1.18	23.8	21.7	1.38	1.31
3	28.8	26.7	28.0	1.51	1.37; 1.39	1.52 (<i>sextet</i> , 6.9)	28.9	26.9	1.66	1.54
4	34.8	32.4	32.9	2.19; 2.33	2.04; 2.20	2.41–2.51 (<i>m</i>)	34.7	32.4	2.50; 2.64	2.46
5	177.5	173.4	172.4	–	–	–	177.3	173.4	–	–
6	65.2	62.9	63.4	4.58 (<i>d</i> , 10.1)	4.45	4.03 (<i>d</i> , 6.5)	68.2	65.7	4.14 (<i>d</i> , 10.7)	4.08
7	29.6	27.5	27.2	2.31	2.21	2.17 (<i>septet</i> , 6.5)	29.5	27.5	2.24	2.13
8	19.7	18.7	19.2	0.84 (<i>d</i> , 6.5)	0.75	0.73 (<i>d</i> , 6.5)	19.6	18.4	0.79 (<i>d</i> , 6.7)	0.70
9	20.9	20.0	20.5	1.00 (<i>d</i> , 6.5)	0.93	0.91 (<i>d</i> , 6.5)	20.3	19.3	1.01 (<i>d</i> , 6.5)	0.93
10	173.9	171.8	174.2	–	–	–	173.3	171.5	–	–
11	50.9	48.7	49.1	4.62; 4.80	4.62	4.41(<i>d</i> , 13.6) 4.49 (<i>d</i> , 13.6)	47.6	45.5	4.60	4.47
12	139.1	137.7	138.1	–	–	–	139.8	137.1	–	–
13	128.1	126.2	126.7	7.24	7.20	7.05 (<i>br d</i> , 8.0)	129.0	126.9	7.18	7.09
14	130.6	128.7	128.2	7.10	7.06	7.18 (<i>d</i> , 8.0)	130.1	128.2	7.02	6.97

Continued

Table 5 NMR Data of Two Conformers of Valsartan—cont'd

#	Major Conformer						Minor Conformer			
	δ_C		δ_H				δ_C		δ_H	
	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆	DMSO- <i>d</i> ₆	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆	DMSO- <i>d</i> ₆	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆	CD ₃ OD ^{a,b}	DMSO- <i>d</i> ₆
15	140.0	137.7	138.2	–	–	–	139.1	138.2	–	–
16	130.6	128.7	128.7	7.10 (<i>d</i> , 8.3)	7.06	6.95 (<i>d</i> , 8.0)	130.1	128.2	7.02 (<i>d</i> , 8.3)	6.97
17	128.1	126.2	128.7	7.24 (<i>d</i> , 8.3)	7.20	7.08 (<i>d</i> , 8.0)	129.0	126.9	7.18 (<i>d</i> , 8.3)	7.09
18	143.4	141.1	141.6	–	–	–	143.6	141.3	–	–
19	132.1	130.5	129.2	7.55	7.53	7.48–7.59 (<i>m</i>)	132.1	130.5	7.53	7.53
20	131.9	130.5	131.0	7.67	7.63	7.59–7.71 (<i>m</i>)	132.8	130.5	7.66	7.63
21	129.3	127.6	127.4	7.56	7.57	7.48–7.59 (<i>m</i>)	129.2	127.6	7.54	7.57
22	132.8	130.9	131.6	7.65	7.68	7.59–7.71 (<i>m</i>)	132.0	130.9	7.64	7.68
23	124.7	123.5	123.8	–	–	–	124.5	123.5	–	–
24	156.9	155.0	155.0	–	–	–	157.0	155.0	–	–
NH	–	–	–	–	–	1.97 (<i>s</i>)				

19

^aData from Li *et al.* [43].^bData from Potamitis *et al.* [44].^cData from Bianchini *et al.* [26].

see Table 3). In addition, a TLC method has also been used for identification of valsartan and hydrochlorothiazide tablets (201). The sample ground tablets (the weight is equivalent to a single tablet) was dissolved in a 2 mL of acetone by sonication for 5 min, followed by centrifugation. The mobile phase was a mixture of EtOAc:dehydrated alcohol:3.6 M ammonium hydroxide. The R_f of the samples should be identical to that of the standard solutions.

British Pharmacopeia 2013 [5] described three methods (A, B, and C) for identification of valsartan, i.e., (A) IR absorption, (B) enantiomeric purity, and (C) specific optical rotation. The identification must be conducted by using either methods A and B or A and C. Method B should be performed by using LC system (see Table 3). The sample used in this method contains impurity A, and the area of the impurity should not be more than the area of the principle peak in the standard solution. For method C, 0.2000 g sample was dissolved in 20 mL methanol, and the value should be between -64.0 and -69.0 .

6.2 Related Compounds

USP 36 [6] described two limit tests for valsartan-related compounds using LC methods. Test 1 is the limit test for related compound A (RS A), while test 2 is for related compound B (RS B), related compound C (RS C), and other related compounds. The percentage of the related compounds can be calculated by using the formula:

$$100(C_s/C_u)(r_u/r_s)$$

in which C_s is the concentration (mg/mL) of the valsartan-related compound A RS in the standard solutions, while C_u is the concentration of valsartan (mg/mL) in the test solutions; r_u and r_s are the peak responses of RS A from the test and standard solutions, respectively. The acceptability limit is 1%. Details of LC methods are described in Table 3.

British Pharmacopeia 2013 [5] described an LC method for determination of the impurity in bulk drugs (see Table 3). The area of impurity C should not be more than twice the area of the principle peak in the chromatogram of reference solution A (0.2%). The area of each unspecified impurity should be less than that of reference solution A (0.1%), and the area of total impurities must not be more than the reference solution A (0.3%).

6.3 Assay

USP 36 [6] described LC methods for the determination of valsartan bulk drugs, valsartan tablets, and valsartan and hydrochlorothiazide tablets. All three LC methods used the same material for column (L1) but with different dimensions and particle sizes, while mobile phases are identical. A gradient elution was applied for analysis of valsartan and hydrochlorothiazide. Details of LC conditions are presented in Table 3.

British Pharmacopeia 2013 [5] used potentiometry for assay of valsartan in bulk drugs. The sample (0.170 g) in 70 mL of 2-propanol R was then titrated with 0.1 M tetrabutylammonium hydroxide in 2-propanol. Analyses were performed under nitrogen, and the end point was determined potentiometrically. In this method, 1 mL of 0.1 M tetrabutylammonium hydroxide in 2-propanol is equivalent to 21.78 mg $C_{24}H_{29}N_5O_3$.

6.4 Dissolution Test

6.4.1 Valsartan Tablets

USP 36 [6] described the determination of valsartan in dissolution medium (phosphate buffer pH 6.8) by using spectrophotometer at UV 250 nm against blank (medium). The percentage of dissolved valsartan can be calculated by using the equation:

$$\text{Result} = (A_u/A_s) \times (C_s/L) \times V \times 100$$

A_u = absorbance of sample solution; A_s = absorbance of standard solution; C_s = concentration of USP valsartan RS in standard solution (mg/mL); L = label claim (mg/tablet); V = volume medium (1000 mL); tolerance is not less than 80% (Q value).

6.4.2 USP Valsartan and Hydrochlorothiazide Tablets

6.4.2.1 Spectrophotometer

The determination of valsartan and hydrochlorothiazide in dissolution medium (phosphate buffer pH 6.8) was performed by using spectrophotometer at UV 250 nm (valsartan) and at 272 nm (hydrochlorothiazide) against blank (medium). The percentage of the dissolved valsartan and hydrochlorothiazide can be calculated using equations:

$$\text{Valsartan : Result} = [(AT_2 \times D) - (AT_1 \times E)/(C \times D) - (B \times E)] \times 12,500$$

$$\text{Hydrochlorothiazide : Result} = [(AT_1 \times C) - (AT_2 \times B)/(D \times C) - (E \times B)] \times 80,000$$

AT_1 = absorbance of sample solution at 272 nm; AT_2 = absorbance of sample solution at 250 nm; $B = A1\% V_{272}$ absorptivity (1%, 0.2 cm, 272 nm) valsartan in medium; $C = A1\% V_{272}$ absorptivity (1%, 0.2 cm, 250 nm) valsartan in medium; $D = A1\% V_{272}$ absorptivity (1%, 0.2 cm, 272 nm) hydrochlorothiazide in medium; $E = A1\% V_{272}$ absorptivity (1%, 0.2 cm, 250 nm) hydrochlorothiazide in medium.

6.4.2.2 Liquid Chromatography

Valsartan and hydrochlorothiazide can be analyzed in dissolution medium using an LC method. The sample diluents are ACN:water (1:1). Details of the LC condition are described in Table 3. The percentage of dissolved valsartan and hydrochlorothiazide can be calculated by equation:

$$\text{Result} = (r_u/r_s) \times (C_s/L) \times D \times V \times 100$$

r_u = peak response valsartan or hydrochlorothiazide in the sample solution; r_s = peak response valsartan or hydrochlorothiazide in the standard solution; C_s = concentration of valsartan or hydrochlorothiazide in standard solution; L = label claim of valsartan or hydrochlorothiazide in tablet (mg/tablet); D = dilution factor (if applicable); V = volume of medium (1000 mL). Tolerance is 80% (Q value).



7. METHOD OF ANALYSIS

7.1 Spectrophotometer

Gupta *et al.* [45] described analysis of valsartan in bulk drugs and tablets. The standard and samples were dissolved in methanol, and the determination of analytes in samples was performed by comparing with standards and value of A (1%, 1 cm). Comparison with standards was performed by using two methods, i.e., zero-order absorption (250 nm) and second-order spectra (at 241 nm). Calibration ranges were linear in the range of 10–50 $\mu\text{g/mL}$ ($r=0.999$). Recovery studies were in the range of 99.08–100.31%. A statistical evaluation showed no significant difference among three methods of estimation. A similar spectrophotometric method was proposed by Kumar *et al.* [46]; in this case, the analysis was performed at λ 249 nm (in methanol), and recovery showed good values (97.77–101.4%).

Din *et al.* [47] reported two simultaneous spectrophotometric determination of valsartan and hydrochlorothiazide in tablets. The first method applied ratio derivative method; in this case, first amplitude derivatives were at λ 231.5 and 260.5 nm (valsartan) and 270.6 nm (hydrochlorothiazide);

mean recover for valsartan was 100.4% (RSD = 1.76%). Second method applied inverse least square absorbance matrix of zero-order spectra in the range of 225–280 nm ($\Delta\lambda = 5$ nm, at 12 wavelengths). Spectra were measured in various binary mixtures; in this case, recovery of valsartan was 101.2% (RSD = 1.58%).

Chaudhary *et al.* [48] determined valsartan and hydrochlorothiazide simultaneously in tablets by using absorption ratio method. Molar absorptivity was calculated at both wavelengths 270.5 nm (λ_{\max} of hydrochlorothiazide) and 231.5 nm (isoabsorptive point). The recovery showed satisfactory results (valsartan: 99.05–102.23%; hydrochlorothiazide: 97.42–100.22%). RSD of in- and interday precision was below 2%.

Gupta *et al.* [49] reported analysis of valsartan and amlodipine in tablets by using simultaneous equation method and absorption correction method. A simultaneous equation method was developed by using λ_{\max} of valsartan (250 nm, MeOH) and amlodipine (238 nm, MeOH), while absorption correction method applied at iso-bestic wavelength of the drugs (236 nm). Recoveries were in the range of 99.48–100.13% for both spectrophotometric methods.

Ramachandran *et al.* [50] developed method for simultaneous determination of valsartan and ezetimibe in tablets by using complex formation with bromophenol blue (BPB) and bromocresol green (BCG). Valsartan–dye complex shows λ_{\max} at 425 and 428 nm, respectively, while ezetimibe cannot react with both BPB and BCG. The analysis was then performed by using two wavelengths for each method of complexation (λ 425 or 428 for valsartan) and at λ_{\max} for ezetimibe (250 nm). Good recovery was achieved using this method (99.3–100.3%).

Birajdar *et al.* [51] reported simultaneous determination of valsartan and nebivolol in tablets using Q method of analysis and simultaneous equations (at λ 280.2 and 246.6 nm). Isoabsorptive point of valsartan and nebivolol (at λ 275 nm) and λ 246.6 nm were selected for Q method. Recoveries for both analytes were in the range of 98.80–101.2% by using both spectrophotometric methods.

7.2 Capillary Electrophoresis

Hillaert and van den Bossche [52] described the separation of six ARA-IIs, i.e., valsartan, candesartan, eprosartan mesylate, irbesartan, losartan potassium, and telmisartan, by using capillary zone electrophoretic method

(CZE). Best separation was achieved by using 60 mM phosphate buffer pH 2.5. In 2003, these authors reported the separation of the same six ARA-II by using micellar electrokinetic capillary chromatography (MEKC) in a 55 mM sodium phosphate buffer solution (pH 6.5) containing 15 mM SDS. Both CZE and MEKC methods were performed in a fused silica capillary column and detected by using UV 214 nm.

Hillaert and van den Bossche [53] also reported simultaneous analysis of combination of the six ARA-IIs and hydrochlorothiazide by using CZE and MEKC methods. The experiments were conducted on a fused silica capillary (total length was 85 cm, 33 cm to the detector; 50 μm i.d.) and by UV 214 nm detection. The running buffers were similar to that of in the previous publications [52,54]. Sulfanilamide (CZE) and eprosartan mesylate (MEKC) were used as internal standard for a quantitative evaluation. The recovery of valsartan, losartan potassium, and irbesartan at three-level concentrations for both methods was in the range of 96.1–102.3%.

7.3 High-Performance Liquid Chromatography

Table 6 summarizes HPLC analysis of valsartan and other drugs in pharmaceutical preparations and raw materials. Most of the methods used C18 as the stationary phase.

7.4 Thin-Layer Chromatography

Silica gel has been reported as the stationary phase in most of the analysis of valsartan by (HP)TLC. The mobile phase usually contains a small amount of acid or base (ammonia), and UV was used for detection. Summary of these works is presented in Table 7.

Inglot and coworkers [102,103] reported the retention behavior of several antihypertension drugs on diol, silica, RP-8, RP-18, and CN plates in various mobile phases. The best separation on diol plate was achieved by using hexane:dioxane:formic acid (3:7:0.1, v/v/v), while a mixture of hexane:isopropanol:formic acid (4:6:0.1, v/v/v) was the best mobile phase for CN plate. An effective separation on RP-8 and RP-18 plates can be achieved by using simple mobile phase containing DMSO:phosphate buffer pH 5 (8:2, v/v). The detection limit of valsartan determined by using densitometry at 254 nm and video densitometry were 0.1 and 0.2 $\mu\text{g}/\text{spot}$, respectively.

Table 6 Summary of HPLC Analysis of Valsartan

Analyte (s)	Sample	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	LOD, LOQ, Recovery (Rec) Valsartan	Ref.
Valsartan, hydrochlorothiazide	Tablet	ACN:acetate buffer (pH 4.0, 0.1 M; 40:60, v/v)	Column: RP Hypersil ODS (200 mm × 4.6 mm, 5 μm) Temperature: (20–25) °C Detector: 220 nm	ACN:acetate buffer (pH 4.0, 0.1 M; 40:60, v/v) Flow rate: 1.0 mL/min Injection volume: n/a	LOD: 1.0 μg/mL LOQ: n/a Rec: 99.75%	[55]
Valsartan, hydrochlorothiazide	Tablet, synthetic mixtures	MeOH	Column: Supelcosil C18 (150 mm × 4.6 mm; 5 μm) Temperature: n/a Detector: 225 nm	0.02 M Phosphate buffer pH 3.2:ACN (55:45, v/v) Flow rate: 0.9 mL/min Injection volume: 20 μL	LOD: 0.017 μg/mL LOQ: 0.058 μg/mL Rec: Tablet: 102.2% Synthetic mixtures: 100.8%	[56]
Valsartan,	Capsule	Std: ACN, diluted with mobile phase Sample: ACN	Column: Shim-pack C18 (250 mm × 4.6 mm; 10 μm) Temperature: room Detector: 265 nm	ACN:PO ₄ buffer pH 2.7 (5:55, v/v) Flow rate: 1.3 mL/min Injection volume: 20 μL	LOD: 0.2 μg/mL LOQ: 1 μg/mL Rec: 99.91%	[57]
Valsartan	Crude valsartan	MeOH	Column: ODS Hypersil C18 (200 mm × 4.6 mm; 5 μm); Temperature: ambient Detector: 254 nm MS: ESI (–) ion mode Scan range: 80–600 <i>m/z</i> 434.2 <i>m/z</i> valsartan	Acetic acid 0.2%:water (55:45, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: n/a LOQ: n/a Rec: n/a	[23]
Valsartan, hydrochlorothiazide	Tablet, impurities	ACN:water (4:6; v/v)	Column: Hypersil 120-5 ODS (150 mm × 4.6 mm; 5 μm) Column temperature: 25 °C Detector: 256 nm	Mobile phase A: ACN:water = 10:90 (v/v, pH 2.5) Mobile phase B: ACN:water = 90:10 (pH 2.5) Gradient time program Flow rate: 1.0 mL/min Injection volume: 50 μL	LOD: 0.02 μg/mL LOQ: 0.004 μg/mL Rec: (103.7–105.3)%	[27]

Valsartan, amlodipine	Capsule	ACN:phosphate buffer (0.02 M, pH 3.0; 56:44, v/v)	Column: Kromasil RP18 (250 mm × 4.6 mm) Column temperature: 25 °C Detector: 234 nm	0.02 M, pH 3.0; 56:44, v/v Flow rate: 1.0 mL/min Injection volume: 20 µL	LOD: 0.018 µg/mL LOQ: 0.054 µg/mL Rec: (99.09–101.79)%	[58]
Valsartan, nebivolol hydrochloride	Capsule	MeOH:1-hexanesulfonic acid monohydrate sodium salt (80:20, v/v)	Column: HIQ sil C18 ODS (250 mm × 4.6 mm, 5 µm) Temperature: ambient Detector: 289 nm	MeOH:1-hexanesulfonic acid monohydrate sodium salt (80:20, v/v) Flow rate: 1.0 mL/min Injection volume: 20 µL	LOD: 17.58 µg/mL LOQ: 53.29 µg/mL Rec: (99.08 ± 1.05)%	[59]
Valsartan, hydrochlorothiazide	Tablet	MeOH:ACN:water:isopropyl alcohol (22:18:68:2, pH 8.0, v/v/v/v)	Column: Diamonsil C18 (250 mm × 4.6 mm; 5 µm) Temperature: n/a Detector: 270 nm	MeOH:ACN:water:isopropyl alcohol (22:18:68:2, pH 8.0, v/v/v/v) Flow rate: 1.0 mL/min Injection volume: 20 µL	LOD: n/a LOQ: n/a Rec: n/a	[60]
Valsartan	Raw material	n/a	Column: Kromasil C18 (250 mm × 4.6 mm, 5 µm) Temperature: n/a Detector: 265 nm	ACN:water:0.5% orthophosphoric acid (40:60, v/v) Flow rate: 1.0 mL/min Injection volume: 20 µL	LOD: 0.25 µg/mL LOQ: 1.0 µg/mL Rec: (98.76–101.30)%	[61]
Valsartan, amlodipine	Tablet	Water:ACN (50:50, v/v)	Column: Xterra RP18 (150 mm × 4.6 mm, 5 µm) Temperature: 25 °C Detector: Valsartan: UV 265 nm	A: 1000 mL water + 0.2 mL trifluoroacetic acid B: Water:ACN:trifluoroacetic acid (400:600:1, v/v/v) Flow rate: 1.5 mL/min Injection volume: 10 µL	LOD: 0.95 µg/mL LOQ: 1.6 µg/mL Rec: (99.76–99.91)%	[62]

Continued

Table 6 Summary of HPLC Analysis of Valsartan—cont'd

Analyte (s)	Sample	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	LOD, LOQ, Recovery (Rec) Valsartan	Ref.
Valsartan	Tablet	NH ₄ H ₂ PO ₄ 0.01 M, pH 3.5:MeOH (50:50, v/v)	Column: Phenomenex mini C18 (250 mm × 4.6 mm; 5 μm) Temperature: 25 °C Detector: 210 nm	NH ₄ H ₂ P ₂ 0.01 M, pH 3.5: MeOH = 50:50 Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.0388 μg/mL LOQ: 0.1176 μg/mL Rec: (99.63–100.587)%	[63]
Valsartan	Degradation product	MeOH	Column: HIQ sil C18 (250 mm × 4.6 mm; 5 μm) Temperature: ambient Detector: 250 nm	MeOH:water (pH 7.2; 70:30, v/v) Flow rate: 1.2 mL/min Injection volume: 20 μL	LOD: n/a LOQ: n/a Acid hydrolysis (valsartan) (98.50–100.20)% Oxidation (valsartan): (98.70–101.50)%	[64]
Valsartan, ramipril	Synthetic mixture, degradation studies	ACN:water (55:45, v/v) pH 3.6	Column: Hypersil C18 (250 mm × 4.6 mm, 5 μm) Temperature: (25 ± 2) °C Detector: 215 nm	ACN:water = 55:45 pH 3.6 Flow rate: 1 mL/min Injection volume: 20 μL	LOD: 0.0280 μg/mL LOQ: 0.0849 μg/mL Rec: (98.42–99.34)%	[65]
Valsartan, amlodipine besylate	Tablet	ACN:phosphate buffer (50:50, v/v)	Column: Zorbax ODS (250 mm × 4.6 mm, 5 μm) Temperature: room Detector: 210 nm	ACN:phosphate buffer (50:50, v/v) Flow rate: 1.0 mL/min Injection volume: n/a	LOD: 0.36 μg/mL LOQ: 1.28 μg/mL Rec: 100.48%	[66]
Valsartan	Capsule	0.1 M phosphate buffer:ACN (20:80, v/v)	Column: Venusil XBP C18 (250 mm × 4.6 mm, 5 μm) Temperature: ambient Detector: 273 nm	0.1 M phosphate buffer:ACN (20:80, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.5 μg/mL LOQ: 1.5 μg/mL Rec: 100.57–100.80%	[67]
Valsartan	Tablet	Water:ACN (50:50, v/v)	Column: Thermo-Hypersil ODS (150 mm × 4.6 mm; 5 μm) Temperature: 25 °C Detector: 273 nm	Water:ACN:glacial acetic acid (200:500:01, v/v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 2.72 μg/mL LOQ: 8.25 μg/mL Rec: (99.0–100.4)%	[68]

Valsartan, nebivolol HCl	Capsule	ACN:water (1:1, v/v)	Column: Phenomenex RP C18 (250 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 275 nm	Ammonium acetate 5 mM, pH 3.5:ACN (30:70, v/v) Flow rate: 0.8 mL/min Injection volume: 50 μL	LOD: 0.032 μg/mL LOQ: 0.095 μg/mL Rec: 99.58%	[51]
Valsartan, hydrochlorothiazide	Tablet	MeOH	Column: Phenomenex Luna C18 (150 mm × 4.6 mm, 5 μm) Temperature: (20 ± 1) °C Detector: 250 nm	ACN:MeOH:PO ₄ buffer 50 mM, pH (3 ± 0.1) 20:50:30, v/v/v Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.2077 μg/mL LOQ: 0.6294 μg/mL Rec: (100.73–102.22)%	[69]
Valsartan, atorvastatin	Tablet	Diluent 1: ACN: MeOH (50:50, v/v) Diluent 2: Water: MeOH (50:50, v/v)	Column: Hypersil BDS C18 (250 mm × 4.6 mm, 5 μm) Temperature: 40 °C Detector: 225 nm	0.1% Acetic acid:ACN (50:50, v/v) Flow rate: 2.0 mL/min Injection volume: 10 μL	LOD: n/a LOQ: n/a Rec: 99.2%	[70]
Valsartan	Pure, tablet	Phosphate buffer pH 3:ACN (50:50, v/v)	Column: Xterra C18 (100 mm × 4.6 mm, 5 μm) Temperature: ambient Detector: 210 nm	Phosphate buffer pH 3:ACN (20:50, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.012 μg/mL LOQ: 0.040 μg/mL Rec: 100.2%	[71]
Valsartan	Tablet	MeOH	Column: C18 (250 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 265 nm	Ammonium dihydrogen phosphate:MeOH (33.5:66.5, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 6 ng/mL LOQ: 18 ng/mL Rec: (99.4–100.6)%	[72]
Valsartan, amlodipine	Tablet	MeOH	Column: Microbondapak C18 (250 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 210 nm	MeOH:potassium hydrogen phosphate 0.1 M, pH 3.0 (65:35, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.02 μg/mL LOQ: 0.06 μg/mL Rec: (98.95–99.76)%	[73]

Continued

Table 6 Summary of HPLC Analysis of Valsartan—cont'd

Analyte (s)	Sample	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	LOD, LOQ, Recovery (Rec) Valsartan	Ref.
Valsartan	Capsule	Ammonium dihydrogen orthophosphate pH 7.5:MeOH (50:50, v/v)	Column: Inertsil C18 (250 mm × 4.6 mm, 5 μm) Temperature: room temp Detector: 210 nm	Ammonium dihydrogen orthophosphate pH 3.5: MeOH (50:50, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.056 μg/mL LOQ: 0.156 μg/mL Rec: (100.01–100.04)%	[74]
Valsartan, hydrochlorothiazide	Tablet	MeOH and mobile phase	Column: C18 Inertsil (250 mm × 4.6 mm; 10 μm) Temperature: room Detector: 259 nm	0.02 M potassium dihydrogen orthophosphate:MeOH: triethylamine (25:75:0.2, v/v/v, pH 6.0) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 16 μg/mL LOQ: 48 μg/mL Rec: (99.39–100.05)%	[75]
Valsartan	Tablet	Phosphate buffer: ACN = 55:45	Column: Kromasil C18 (250 mm × 4.6 mm, 5 μm) Temperature: ambient Detector: 233 nm	Phosphate buffer:ACN 55:45, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.034 μg/mL LOQ: 0.104 μg/mL Rec: 99.97%	[76]
Valsartan, aliskiren hemifumarate	Tablet	ACN:0.05 M potassium dihydrogen PO ₄ buffer pH 3.5 (45:55, v/v)	Column: perchrom ODS BP (200 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 224 nm	ACN:0.05 M Potassium dihydrogen PO ₄ buffer pH 3.5 (45:55, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.93 μg/mL LOQ: 2.84 μg/mL Rec: (99.32–99.92)%	[77]
Valsartan, ramipril	Capsule	MeOH	Column: RP 18 (250 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 225 nm	Phosphate buffer 1%:ACN 60, v/v, pH 3.2) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 1 μg/mL LOQ: 3 μg/mL Rec: (97.67–100.23)%	[78]

Valsartan, irbesartan, losartan potassium, telmisartan	Tablet	MeOH	Column: ACE RP-C18 (250 mm × 4.6 mm, 5 μm) Temperature: 40 °C Detector: 220 nm	Potassium dihydrogen phosphate 0.025 M, pH 6.0 (5:35, v/v) Flow rate: 1.5 mL/min Injection volume: 50 μL	LOD: 0.13 μg/mL [79] LOQ: 0.4 μg/mL Rec: (99.87–101.35)%
Valsartan, amlodipine, hydrochlorothiazide	Tablet	ACN:phosphate buffer (0.05 M) pH 2.8 (40:60, v/v)	Column: Phenomenex Kinetex (150 mm × 4.6 mm, 5 μm) Temperature: (22–25) °C Detector: 227 nm	ACN:phosphate buffer (0.05 M) pH 2.8 (40:60, v/v) Flow rate: 0.8 mL/min Injection volume: 20 μL	LOD: 1.42 μg/mL [80] LOQ: 4.31 μg/mL Rec: (98.9–101.0)%
Valsartan	Tablet	MeOH	Column: Phenomenex Luna C18 (250 mm × 4.6 mm; 5 μm) Temperature: n/a Detector: 210 nm	MeOH:PO ₄ buffer pH 3.0 (65:35, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.02 μg/mL [81] LOQ: 0.06 μg/mL Rec: (98.79–99.81)%
Valsartan	Tablet	MeOH	Column: Kromasil C18 (250 mm × 4.6 mm, 5 μm) Temperature: n/a Detector: 250 nm	Gradient condition of ACN and PO ₄ buffer pH 3.5 Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: n/a [82] LOQ: n/a Rec: (99.27–99.92)%
Valsartan, hydrochlorothiazide	Tablet	ACN:water (1:1, v/v)	Column: Xterra (250 mm × 4.6 mm, 5 μm) Temperature: room Detector: 265 nm	Ammonium acetate 0.20 M, pH 5.6:ACN (gradient elution) Flow rate: 1.5 mL/min Injection volume: 20 μL	LOD: 0.0375 μg/mL [83] LOQ: 0.064 μg/mL Rec: (100.1–102.3)%
Valsartan	Tablet	Water:ACN (1:1, v/v)	Column: Phenomenex C18 (250 mm × 4.6 mm; 5 μm) Temperature: ambient Detector: 240 nm	ACN:5 mM ammonium acetate pH 4.5 (75:25, v/v) Flow rate: 1.0 mL/min Injection volume: 50 μL	LOD: 10 ng/mL [84] LOQ: 28 ng/mL Rec: (99.85 ± 1.03)%

Continued

Table 6 Summary of HPLC Analysis of Valsartan—cont'd

Analyte (s)	Sample	Solvent	Column and Condition	Mobile Phase, Flow Rate, Injection Volume	LOD, LOQ, Recovery (Rec) Valsartan	Ref.
Valsartan, amlodipine besylate	Capsule	MeOH	10 μ m: Waters Symmetry C18 (150 mm \times 4.6 mm; 5 μ m) Temperature: 25 $^{\circ}$ C Detector: 238 nm	49 OH:KH ₂ PO ₄ buffer 0.01 M, pH 2.5 (60:40, v/v) Flow rate: 1.0 mL/min	LOD: 44 ng/mL LOQ: 300 ng/mL Rec: (98.97–99.97)%	[85]
Valsartan, hydrochlorothiazide	Tablet	0.02 M Phosphate buffer:ACN: MeOH = 50:40:10	10 Column: Supelcosil C18 (150 mm \times 4.6 mm; 5 μ m) Temperature: ~25 $^{\circ}$ C Detector: 225 nm	48 0.02 M Phosphate buffer pH 2.9:ACN:MeOH (50:40:10, 25 /v) Flow rate: 1.4 mL/min Injection volume: 50 μ L	LOD: n/a LOQ: n/a Rec: 99.94%	[86]
Valsartan	Tablet	35 KH ₂ PO ₄ buffer: ACN (pH 3.0 \pm 0.1, 55:45, v/v)	33 Column: Xterra C18 (150 mm \times 4.6 mm; 5 μ m) Temperature: 23 \pm 1 $^{\circ}$ C Detector: 286 nm	1 KH ₂ PO ₄ buffer: 0.1 N (pH 3.0 \pm 0.1, 55:45, v/v) Flow rate: 0.7 mL/min Injection volume: 20 μ L	20 LOD: 0.0180 μ g/mL LOQ: 0.1016 μ g/mL Rec: (99.40–99.68)%	[87]
Valsartan, ramipril	Capsule	n/a	1 Column: Hypersil C18 (150 mm \times 4.6 mm, 5 μ m) Temperature: n/a Detector: 220 nm	18 20 mM Phosphate buffer: ACN (35:65, v/v) Flow rate: 0.8 mL/min Injection volume: 20 μ L	LOD: n/a LOQ: n/a Rec: (98.4–99.6)%	[88]
Valsartan, amlodipine besilate, olmesartan medoxomil, hydrochlorothiazide	Tablet	MeOH diluted with ACN MeO 28 water (7:13:80, v/v/v)	Column: CN (200 mm \times 4.6 mm, 5 μ m) Temperature: 30 $^{\circ}$ C Detector: 235 nm	6 ACN:MeO 10 mM phosphoric acid (pH 2.5, 2 13:80, v/v/v) Flow rate: 1.0 mL/min Injection volume: 20 μ L	LOD: 0.1 μ g/mL LOQ: 0.3 μ g/mL Rec: 95.8%	[89]

Valsartan, amlodipine hydrochlorothiazide	Tablet	MeOH	Column: ⁵⁶ Zorbax SB-C8 (250 mm × 4.6 mm; 5 μm) Temperature: 25 °C Detector: 225 nm	0.025 M phosphoric acid: ² ACN (gradient) Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.24 μg/mL [90] LOQ: 0.80 μg/mL Rec: 100.65 ± 1.23%
Valsartan, amlodipine besylate, hydrochlorothiazide	Tablet	MeOH	Column: ⁷ Kromasil KR-5 C18 (250 mm × 4.6 mm; 5 μm) Temperature: n/a Detector: 232 nm	50 mM KH ₂ PO ₄ pH 3.7: ACN (56:44, v/v) Flow rate: 1.0 mL/min Injection volume: 20 μL	³⁴ LOD: 1.1 μg/mL [91] LOQ: 3.3 μg/mL Rec: 99.69 ± 0.63% to 100.20 ± 0.07%
Valsartan, propranolol hydrochloride	Raw material, gel formulation	ACN: MeOH:0.01 M Na ₂ PO ₄ pH 3.5 = 50:35:15	Column: ¹ Hypersil ODS C18 (250 mm × 4.6 mm, 5 μm) Temperature: 25 ± 0.2 °C Detector: 250 nm	⁸⁴ ACN:MeOH:0.01 M Na ₂ HPO ₄ pH 3.5 = 50:35:15 Flow rate: 1.0 mL/min Injection volume: 20 μL	LOD: 0.45 ⁴² mL [92] LOQ: 1.39 μg/mL Rec: (99.79–102.93)%

Table 7 Summary of TLC Analysis of Valsartan

Analyte(s)	Sample	Stationary Phase	Mobile Phase, Solvent	Detection	LOD, LOQ, Recovery, Precision (Valsartan)	Ref.
Valsartan, amlodipine	Bulk drug, tablet	Silica gel F254	Mobile phase: toluene: MeOH:acetic acid (7:3:0.1, v/v/v) Solvent: MeOH	UV 244 nm	LOD: 50 ng/spot LOQ: 100 ng/spot Recovery: (98.72 ± 0.5)% RSD precision: Intraday: (0.14–0.39)% Interday: (0.19–0.52)%	[93]
Valsartan, ramipril	Capsules	Silica gel F254	Mobile phase: ethyl acetate: chloroform:glacial acetic acid (8:2:0.2, v/v/v) Solvent: MeOH	UV 220 nm	LOD: 42.4 ng/spot LOQ: 1400.8 ng/spot Recovery: (98.92–101.22)% RSD precision: Intraday: 0.037% Interday: 0.046%	[94]
Valsartan, telmisartan	Tablets	Silica gel F254	Mobile phase: 1,4-dioxane: heptane:formic acid 99% (5:5:0.1, v/v/v) Solvent: MeOH	Video scanning	LOD: n/a LOQ: n/a Recovery: (99.62–101.59)% Intermediate precision (RSD): 5.14%	[95]
Valsartan, ramipril	Capsule	Silica gel 60F ₂₅₄	Mobile phase: chloroform: ethyl acetate:methanol: glacial acetic acid (5:5:1:0.2, v/v/v/v) Solvent: n/a	UV 210 nm	LOD: 100 ng/spot LOQ: 330 ng/spot Recovery: 99.099.7% RSD precision: Intraday: 0.136–0.277% Interday: 0.172–0.259%	[96]

Valsartan, hydrochlorothiazide	Tablet	13 Silica gel 60F ₂₅₄	Mobile phase: CHCl ₃ : MeOH:NH ₄ OH (8:2:1, v/v/v) Solvent: MeOH	UV 225 nm	Recovery: 99.94% RSD precision: Intraday: 1.22% Interday: 1.59%	[86]
Valsartan, telmisartan, potassium losartan	Tablet	Silica gel 60F ₂₅₄	Mobile phase: CHCl ₃ : 6 MeOH:acetone:toluene: acetic acid (7.5:1.5:5:5:0.01:0.003, v/v/v/v/v) Solvent: mobile phase	UV 254 nm	Recovery: 98.74–101.1% Precision (RSD): 0.85%	[97]
Valsartan	Bulk drug, formulation	Silica gel 60F ₂₅₄	Mobile phase: 30 toluene: ethyl acetate:MeOH: formic acid (60:20:20:1, v/v/v/v) Solvent: MeOH	UV 250 nm	LOD: 25 ng/band LOQ: 150 ng/band Recovery: 99.30–101.09% RSD precision: Intraday: 0.76–1.52% Interday: 0.40–1.35%	[98]
Valsartan, hydrochlorothiazide	Tablets	Silica gel 60F ₂₅₄	CHCl ₃ : 7 MeOH:acetone: toluene:acetic acid (6:2:1:0.1, v/v/v/v) Solvent: MeOH	UV 260 nm	LOD: 100 ng/spot LOQ: 300 ng/spot Precision: 0.27–1.19%	[99]

Continued

Table 7 Summary of TLC Analysis of Valsartan—cont'd

Analyte(s)	Sample	Stationary Phase	Mobile Phase, Solvent	Detection	LOD, LOQ, Recovery, Precision (Valsartan)	Ref.
Valsartan, losartan, irbesartan, candesartan, eprosartan mesylate, telmisartan	Tablets	HPTLC F ₂₅₄	CH ₂ Cl ₂ :EtOA:EtOH: Gla ⁷⁶ acetic acid:H ₂ O (45:40:5:1:0.5, v/v/v/v/v) Solvent: n/a	UV 260 nm	LOD: 0.02 µg/spot LOQ: 0.5 µg/spot Recovery: 98–102% Precision: n/a	[100]
		RP- Diphenyl-F	ACN:MeOH:0.1 M ammonium acetate: 4% ammonia (30:20:50:0.5, v/v/v/v/v) Solvent: n/a			
Valsartan, nebivolol HCl	Dosage form	Silica gel	EtOA:MeOH:25% ammonia (12:2:1, v/v/v) Solvent: n/a	UV 280 nm	LOD: n/a LOQ: n/a Recovery: n/a Precision: n/a	[101]
Valsartan, hydrochlorothiazide, amlodipine besylate	Tablets	Silica gel	EtOA:MeOH:10% ammonia (17:4:2, v/v/v) Solvent: n/a	UV 320 nm	LOD: 3200 ng/zone LOQ: 6400 ng/zone Recovery: (98–101)% Precision (RSD): Intraday and interday: <0.8%	[101]



8. ANALYSIS OF VALSARTAN IN BIOLOGICAL FLUIDS AND HERBAL PREPARATION

HPLC is the most common method used in the analysis of valsartan and related compounds, biological fluids, and samples. In this technique, an RP-18 column is normally fitted, and detection was performed by various methods including UV, fluorescent, and MS (SIM and MRM (multiple reaction monitoring)).

Maurer and coworkers [42] reported the analysis of the methylation product of valsartan and other drugs in human urine by using GC-EIMS. The GC conditions used were splitless injection, column 12 mm × 0.2 mm i.d. 10^5 column, injection temperature of 280 °C, column temperature of 100–310 °C at 30 °C/min, carrier gas He 1 mL/min, EI mode (70 eV, source 220 °C, interface 260 °C).

Valsartan is metabolized by CYP-450 isozyme into valeryl-4-hydroxy-valsartan in human. This metabolite has 200-fold lower affinities to the AT1 receptor than valsartan. The renal clearance of valsartan has reported only 30% of the total plasma clearance. The recovery of valsartan are 83% in feces and 13% in urine from oral solution administration of valsartan [2,3,104].

In the quantitative analysis of valsartan in a complex matrix, it is recommended to use an LC–MS technique with MRM due to its low LOQ and specificity. In this method, two daughter ions should be selected: one ion for qualification and the other for quantification purposes. The intensity ratio of qualifier and quantification ions must be compared between standard and samples. The qualifier and quantification ions as described by Ref. [105] in MS–MS Agilent (precursor ion m/z 436.2; ESI positive mode; fragmentation voltage: 98 V; collision energy: 14 eV; dwell time: 20 ms) for valsartan are m/z 235 and 291.1, respectively.

The HPLC methods that applied for analysis of valsartan and related compounds in biological fluids and samples are summarized in Table 8.

Kesting and coworkers [11] investigated the presence of adulterant in the Chinese herbal medicines which were used for the treatment of hypertension. The identification of adulterant was conducted by using LC–HRMS and LC–MS–SPE–NMR. The results showed the presence of valsartan, amlodipine, indapamide, and other drugs used for antihypertension as the commercial soft capsules of herbal drug preparations [11].

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
Valsartan, hydrochlorothiazide	<p>Column 29 Zorbax SB-Aq C-18 4 50 mm × 4.6 mm; 5 μm Mobile phase: ACN:10 mM ammonium acetate (60:40, v/v; pH 4.5) Flow rate: 1.2 mL/min Column temperature: 35 °C Detection: MRM (negative ion mode): valsartan (<i>m/z</i>, quantification: 434.2 → 350.2; identification: 434.2 → 179.1) HCT (<i>m/z</i> quantification: 295.9 → 268.9; identification: 295.9 → 204.9) Probenecid: (<i>m/z</i> 283.9 → 239.9)</p>	Human plasma	Probenecid	<p>Standard solutions: Stock solutions of valsartan (4 mg/mL), HCT (1 mg/mL), and probenecid (5 mg/mL) were prepared by dissolving in ACN:water (80:20, v/v); working solution IS (150 ng/mL) Samples: 100 μL IS + 100 μL plasma sample + 300 μL ACN vortex-mixed 1 min, then centrifuged 10 min (10,000 × <i>g</i>), upper layer collected; 1 11 aliquot was injected into LC-MS/MS</p>	<p>Accuracy: ±7.85% Precision: Intraday: (1.92–4.81)%; Interday: (2.21–7.12)% LOD: n/a LOQ: 4 ng/mL</p>	[106]
Valsartan	<p>32 column (20 mm × 2.1 mm) packed with OASIS M 32 30 μm; SPE was connected 10-port switching valve Samples (40 μL) was flushed from sample loop at 2 mL/min using water:acetic acid (85:15), after loading in SPE, and then SPE was washed with water and THF for 0.5 min. Elution was performed using THF:water:formic acid (95:5:5, v/v/v)</p>	Human plasma and urine	Candesartan	<p>Standard solutions: Stock solutions, valsartan (4 mg/mL) and candesartan (1 mg/mL), were dissolved in MeOH with 0.5% NH₄OH. Calibration and QC samples were prepared by spiking matrices with standard solutions Samples: Samples were centrifuged 5 min (2400 × <i>g</i>); 50 μL sample then deposited in 96-well polypropylene plate; after</p>	<p>Accuracy: Plasma: (101–108)% Urine: (41.7–46.4)% Precision: Intraday: Plasma: (2.0–8.5)% Urine: (0.6–9.3)% Interday: Plasma: (2.1–16.9)% Urine: (2.7–5.0)% LOD: n/a LOQ: n/a</p>	[107]

	Detection: MRM (positive ion mode): valsartan: (m/z 436.2 \rightarrow 291.2.1) Candesartan: (m/z : 441.2 \rightarrow 263.2)			addition of 50 μ L acetic acid (15%, v/v) and 10 μ L IS solutions, the plates were sealed and vortexed at 1000 rpm (10 s), and kept 4 °C until injection		
Valsartan	Column: Phenomenex C18 (10 mm \times 4.6 mm; 5 μ m) equipped with guard column NovaPak C8 Column temperature: ambient Mobile phase: phosphate buffer 5 2.8:ACN (70:30, v/v) Flow rate: 1.3 mL/min Detection: fluorescent Excitation: 265 nm Emission: 378 nm	Human plasma	Losartan	Standard solutions: Stock solutions of valsartan (50 μ g/mL) were prepared by dissolving in 0.1 M KOH and pH was adjusted to 8 with 1 M HCl. IS was prepared with same method (2 μ g/mL) Samples: – 1 mL Human plasma + 100 μ L IS, acidified with 125 μ L phosphoric acid to pH 2.5, add 10 mL methyl- tert-butyl ether (MTBE) then vortex-mixed and centrifuged (1800 \times g for 5 min) – Organic layer was transferred to tube containing 200 μ L 0.05 M NaOH and then the mixture was vortexed and centrifuged (1800 \times g for 5 min) – Organic layer was discarded, aqueous layer neutralized with 75 μ L 0.2 M phosphoric acid, 125 μ L was injected into HPLC	Accuracy and precision: [108] Intraday: (91.8–100.6)% RSD: (0.7–10.8%) Interday: (88.3–100.6)% RSD: (0.8–8.9%) LOD: n/a LOQ: 10 ng/mL	

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
Valsartan, losartan, irbesartan, candesartan cilexetil, and candesartan M1	<p>27 mm: μBondapak C18 (300 mm \times 3.9 mm; 10 μm)</p> <p>+ gr. 19 column Novapak C18 (20 mm \times 3.9 mm; 4 μm)</p> <p>Mobile phase: gradient elution 5 mM acetate buffer pH 4 and ACN</p> <p>Column temperature: room</p> <p>Detection: Fluorescence detector</p> <p>Excitation: 250 nm</p> <p>Emission: 375 nm</p>	Human plasma	Bumetanide	<p>Standard solutions: Stock solutions were prepared in MeOH and ACN</p> <p>Samples:</p> <ul style="list-style-type: none"> – 0.25 mL plasma was spiked with appropriate amount of stock solutions and IS, then acidified with 0.25 mL 1 M H₃PO₄ and the mixture was then shaken and centrifuged 5 min (10,000 \times g at 4 °C) – Clean up was performed 75 g Bond Elut C8 that was conditioned with 2 mL MeOH followed by 1 mL 0.1 M phosphate buffer pH 2 – Samples was slowly passed into the cartridge (0.5 mL) – Column was wa 74 with 0.5 mL MeOH–1 M phosphate buffer pH 2 (50:50, v/v) and dried. – Analyte then eluted with 0.5 mL MeOH – 0.1 mL ethylene glycol 10% (v/v) was added 	<p>Accuracy and precision: Intraday: Relative errors: –2.6% to +3.7% RSD: 3.4–5.4%</p> <p>Interday: Relative error: –4.6% to +6.7% RSD: 5.0–8.0% LOD: n/a LOQ: 51.5 ng/mL</p>	[32]

				<ul style="list-style-type: none"> - The eluate was 10exed and dried under N₂ at 40 °C - The residue was dissolved in starting mobile phase (0.25 mL) - 20 µL solution was injected to the HPLC system 	
Valsartan, amlodipine	Column: HICHROM 13 eosil 100-5 C18 (250 mm × 4.6 mm) Temperature: room Mobile phase: phosphate buffer (0.01 M, pH 3.6): 9 N: MeOH (50:40:10, v/v/v) Flow rate: 1.0 mL/min Detection: UV 240 nm	Rat liver perfusate	n/a	Standard solutions: Standard solutions were prepared in MeOH (100 µg/mL) 73 ples: Spiked samples were prepared by adding 50 µL standard solutions into blank liver perfusate (100 µL) and then making up into 1000 µL by mobile phase, after filtering (0.45 µm) then injected into HPLC	Accuracy and precision [109] Intraday: Bias: -1.03% to +12.67%; RSD: 2.20-4.19% Interday: Bias: +0.89% to +12.33%; 20 D: 2.93-4.74% LOD: 0.02 µg/mL LOQ: 0.05 µg/mL
Valsartan	Column: Nucleosil C18 (50 mm × 4 mm) + guard column Phenomenex (4 mm × 3 mm) Mobile phase 9 ACN:15 mM KH ₂ (PO ₄) ₃ pH 2.0 (45:55, v/v) Flow rate: 1 mL/min Column temperature: 40 °C Detection: fluorescence Excitation: 234 nm Emission: 374 nm	Human plasma	n/a	Standard solutions: Stock solutions were prepared by dissolving standard (20 mg) in 25 mL MeOH; calibration, and QC samples were prepared by spiking standard solutions into drug-free plasma Samples: 0.2 mL plasma + 1 mL MeOH then vortex-mixed 30 s (2000 rpm) and then centrifuged 3 min (2500 × g); 5 µL supernatant injected into HPLC system	Intraday: Bias: +1.8% to +7.2%; RSD: 1.1-3.0% Interday: Bias: -0.8% to +1.4%; RSD: 1.3-4.9% LOD: 25 pg/mL LOQ: 97.9 ng/mL [110]

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
Valsartan	<p>5 Column: Chromolith RP-18e (100 mm × 4.6 mm) Mobile phase: 36 M Na₂HPO₄ buffer:ACN (60:40, v/v, adjusted to pH 3.5) Flow rate: 2 mL/min Column temperature: n/a Detection: fluorescence Excitation: 230 nm Emission: 295 nm</p>	Human plasma	Atorvastatin	<p>72 Standard solutions: Stock solutions (2 mg/mL) were prepared in MeOH Samples: 450 µL plasma + 50 µL IS (12 µg/mL) + 500 µL ACN then mixed and centrifuged for 10 min (8000 rpm); 20 µL was injected into HPLC</p>	<p>Recovery (average): 96.3 ± 1.3% Precision (%CV) Intraday: 2.6–4.3% Interday: 3.0–5.2% LOD: n/a LOQ: 20 ng/mL</p>	[111]
Valsartan, losartan, telmisartan	<p>17 Column: Chromolith 18e monolithic (25 mm × 4.6 mm) Mobile phase: 5 mM phosphate buffer pH 3.8:ACN:MeOH (65:20:15, v/v/v) Flow rate: 3.0 mL/min Temperature: room Detection: fluorescent Excitation: 259 nm Emission: 399 nm</p>	Human urine	n/a	<p>Standard solutions: Individual stock solution (1000 µg/mL) were prepared in MeOH Samples: – 1 mL sample were transferred into 16 mm × 125-mm disposable glass tubes and 3 mL of extraction solvent 1 (2% MeOH in 5 mM phosphate buffer pH 3.8) and then vortexed and centrifuged 10 min (2500 × g) – 2 mL of supernatant was transferred into glass screw-cap vials, then processed by using SPE (L 19 rocart 25-4 LiChrospher RP-18 ADS (25 mm × 4 mm; 25 µm)), and extracted using solvent 1</p>	<p>Recovery: (92.9–102.5)% Precision (CV) Intraday: <3.5% 71 rday: <3.5% LOD: 0.001 µg/mL LLOQ: 0.003 µg/mL</p>	[112]

Valsartan, hydrochlorothiazide	Col 10: Betabasic C8 (50 mm × 4.6 mm; 5 μm) Mobile phase: ACN: MeOH:0.001% aqueous ammonia (75:15:10, v/v/v) Flow rate: 0.5 mL/min Column temperature: 50 °C MRM: Valsartan: (negative ion mode: <i>m/z</i> 434.25 → 179.22) HCT: <i>m/z</i> 295.85 → 204.86 Clonazepam: <i>m/z</i> : 314.0 → 278.27	Human plasma	Clonazepam	Standard solutions: Stock solutions (5 mg/mL) and IS (1 mg/mL) were prepared in MeOH:water (50:50, v/v). This stock solution can be used for spiking in plasma Samples: – 500 μL plasma + 25 μL IS + 500 μL orthophosphoric acid 2% in 1.7-mL microtube then vortexed – The mixture then applied to HLB cartridge (30 mg, 3 mL) previously conditioned with 1 mL MeOH and 2 mL water – The loaded cartridge was washed with 2 mL water, then 1 mL 5% MeOH in water, then eluted with 1.0 mL MeOH:ACN (90:10, v/v), and loaded into autosampler vial	Mean recovery: 90.28% CV: 0.12% Precision: Intraday: <10.0% Interday: <10.0% LOD: n/a LLOQ: 50.0 ng/mL	[41]
Valsartan	Col 6: Zorbax SB-C18 (50 mm × 4.6 mm; 3.5 μm) Mobile phase: ammonium acetate 8 mM, pH 3:MeOH (50:70, v/v) Flow rate: 1 mL/min Temperature: room Detection: MS/MS (no detailed data reported)	Human blood	Valsartan- <i>d</i> ₉	Standards: not reported Samples: – Plasma samples were vortexed (3000 rpm) and centrifuged (1900 × <i>g</i> , 5 min, 4 °C); – 100 μL aliquot + 600 μL solution of IS in MeOH were vortexed (3000 rpm)	Accuracy: 98.66% Precision: 3.51% LLOQ: 20.18 ng/mL	[113]

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
				and then centrifuged (1900 × g, 10 min, 4 °C). – 400 µL supernatant was evaporated into dryness with 65 – Residue was reconstituted in 200 µL in mobile phase and measured with LC–MS/MS		
Valsartan, nebivolol	Column: Gemini-C18 55 (50 mm × 2.0 mm; 3 µm) Guard column C18 17 mm × 3 mm Mobile phase: ACN 4 .05 mM formic acid buffer (50:50, v/v, pH 3.5) Flow rate: 0.25 mL/min Column temperature: 20 °C	Human plasma	Tamsulosin	Standard solutions: Analytes and IS were prepared individually in ACN (100.0 µg/mL). The stocks were diluted with water to give serial of standard solutions; QC samples were prepared by spiking standards into plasma	Accuracy: 101.1–103.7% Precision: Intrarun (RSD): 1.6–2.5% Interrun (RSD): 2.9–4.8% LOD: n/a LLOQ: 1.0 ng/mL 81	[114]
	Detection: Valsartan (negative ion mode) <i>m/z</i> : 434.2 → 179.0 nebivolol (positive ion mode) <i>m/z</i> : 406.1 → 150.9 IS (positive ion mode) <i>m/z</i> 409.4 → 228.1			Samples: 200 µL plasma + 400 µL ACN contained IS (240 ng/mL), then vortexed (1 min) and centrifuged (14,000 rpm; 10 min); aliquots were filtered through 0.2 µm and injected into LC–MS–MS system		

Valsartan, valeryl-4-hydroxy-valsartan	<p>70</p> <p>Column: Waters Atlantis dC18 (100 mm × 3.9 mm; 3.9 μm)</p> <p>Guard 27 umn: μBondapak C18 (10 μm)</p> <p>Mobile phase: gradient elution using 0.025% TFA in ACN and 0.025% TFA in buffer phosphate buffer (5 mM, pH 2.5)</p> <p>Flow rate: 1.3 mL/min</p> <p>Temperature: 40 ± 0.2 °C</p> <p>Detection: fluorescence</p> <p>Excitation: 234 nm</p> <p>Excitation: 378 nm</p>	Human plasma	Candesartan M1	<p>Standard solutions:</p> <p>Stock solutions valsartan, metabolite, and IS were prepared in MeOH with concentrations 107, 186, and 99.26 μg/mL, respectively.</p> <p>Working solutions prepared by diluting the stock solution</p> <p>Samples:</p> <ul style="list-style-type: none"> - 1 mL blank samples were spiked with valsartan, metabolite, and 1 to achieve concentrations 1.1 μg/mL (analytes) and 1.2 μg/mL (IS) and then 1 mL of H₃PO₄ (0.5 M) was added, vortexed, and centrifuged (5 min, 10,000 rpm) - The mixtures 5 were applied to a C8 SPE (conditioned with 2 mL MeOH and 1 mL phosphate buffer pH 2, 60 mM) - SPE was 64 shed by 1 mL MeOH:phosphate buffer (40:60, v/v) - After drying for 8 min, SPE was eluted using solution of 	<p>Recovery: [104]</p> <p>Valsartan: (96.6–101.2)% (RSD: 1.2–3.1%)</p> <p>Metabolite: 94.6–108.8% (RSD: 0.7–1.6%)</p> <p>LOD: n/a</p> <p>LLOQ: 5 ng/mL (valsartan and metabolites)</p>
--	---	--------------	----------------	--	--

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
				10% ethylene glycol in MeOH (0.1 mL) – The compounds were eluted using 0.5 mL diethylether. – Evaporated with N ₂ stream at 60 °C – Residue + 100 µL mobile phase was reconstituted, vortexed, mixed, and filtered with polypropylene membrane Ø = 13 mm, 0.45 µm – 20 µL was injected into HPLC system		
Valsartan, valeryl-4-hydroxy-valsartan	Same with Iriarte 2006	Human plasma	Candesartan M1	Same with Iriarte 2006 (except the nominal concentrations were different). This work reported the application of Iriarte 2006 for bioavailability studies	LOQ for both of analytes was 5 ng/mL Accuracy and precision of analytes: <15%	[115]
Valsartan, amlodipine, hydrochlorothiazide	38 mm: Gemini C18 (250 mm × 4.6 mm; 5 µm), protected with guard column 12 (4 mm × 2 mm) Column temperature: 35 °C Mobile phase: ACN:10 mM ammonium formate pH 3.5, v/v (gradient from 20:80 to 70:30)	Human plasma	Telmisartan	69 dard solutions: Stock solutions (100 µg/mL) of the analytes were prepared in MeOH Samples: 80 ma samples + 1 mL mixture of MeOH:ACN (50:50, v/v) + 15 µL IS solution, vortexed (3 min) and centrifuged 10,000 rpm (10 min),	Recovery: 76.5 ± 2.1%; Precision: Intraday: 101.1–101.8%; Interday: 101.7–107.0% LOD: 7 ng/mL LOQ: 20 ng/mL	[116]

	12 Flow rate: 1 mL/min Detector: PDA (254 nm)			supernatant filtered, and evaporated with N ₂ Dried residue + 0.5 mL mobile phase reconstituted, filtered, and then injected into HPLC		
Valsartan	13 Column: Hypersyl BDS C18 (250 mm × 4.6 mm; 5 μm) Mobile phase: buffer triethylamine pH 3.0:ACN (45:55, v/v) Flow rate: 0.7 mL/min Column temperature: 25 °C Detection: UV 215 nm	Human plasma	Losartan potassium	Standard solutions: not reported Samples: 0.2 mL plasma was transferred into 2-mL centrifuge tube + 25 μL IS (1 mg/mL), vortexed 30 s (2000 rpm) + 1 mL MeOH, vortexed for 10 min (2000 rpm), then centrifuged 10 min (2500 rpm), and 20 μL supernatant injected into HPLC	Accuracy: n/a Precision: n/a LOD: n/a LOQ: n/a	[117]
Valsartan, gliclazide, benazepril	10 Column: Shimadzu VP-ODS C18 (250 mm × 2.0 mm, 5 μm) Temperature: 40 °C Mobile phase: gradient elution of MeOH and 0.05% formic acid Flow rate: 0.2 mL/min Column temperature: n/a Detection: SRM: Valsartan (<i>m/z</i> 436 → 207) Gliclazide (<i>m/z</i> 324 → 110) Benazepril (<i>m/z</i> 425 → 351) IS (<i>m/z</i> 531 → 489)	Human plasma	Ketoconazole	Standard solutions: Stock solutions were prepared in MeOH (valsartan 500 μg/mL; gliclazide and benazepril 50 μg/mL); calibration and QC samples were prepared by spiking into human blank samples Samples: 100 μL + 10 μL IS (125 μg/mL) + 500 μL MeOH, then vortexed and centrifuged for 10 min (40,000 × <i>g</i>), supernatant transferred into another clean test tube, centrifuged again for 10 min (40,000 rpm), 10 μL clear supernatant was injected into LC-MS	Recovery: (93.52–99.94)% Precision: Intraday: 6.2% Bias: 0.1% Interday: 10.4% Bias: 9.0% LOD: n/a LOQ: 20 ng/mL	[38]

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
Valsartan, verapamil, losartan, telmisartan, irbesartan, flecainide, beta blocker drugs	<p>19 mm: Atlantis dC18 (150 mm × 2.1 mm; 3 μm)</p> <p>Mobile phase: gradient elution. Solvent A: 10 mM ammonium formate pH 3.1; B: ACN linear gradient from 90% A and 10% B to 10% A and 90% B in 63 min, then held 3 min</p> <p>Column temperature: 25 ± 0.8 °C</p> <p>Flow rate: 0.3 mL/min</p> <p>Detection: SIM (<i>m/z</i>): Valsartan (207.3); telmisartan (276.3); losartan (207.3); irbesartan (207.3)</p>	Post mortem whole blood	Diazepam-d5	<p>Standard solutions: Two stock solutions were prepared in MeOH (2500 μM); aqueous calibration solutions and aqueous control solutions were prepared from stock solutions; IS was dissolved in ACN (155 μM) and water (18.5 μM). Spiked calibration and QC samples were prepared by adding aqueous solutions into drug-free sodium fluoride whole blood</p> <p>Samples: 0.5 mL whole blood + 50 μL (18.5 μM IS), then this mixture was precipitated with 1 mL ice cold mixture of ACN:MeOH (85:15, v/v), then frozen for 30 min (−20 °C), and centrifuged (2260 × <i>g</i>, 4 °C, 10 min). Supernatant was decanted and mixed with 0.2 mL 4.5 M HCl before diluting with water, and the samples were applied to SPE (Oasis MCX)</p>	<p>Accuracy: n/a</p> <p>Precision: intraday Bias: −7.3%, and −22%; RSD: 21% and 4% for theoretical concentration 5 and 1 μM, respectively</p> <p>Interday: Bias: −20%, and −21%; RSD: 13% and 28% for theoretical concentration 5 and 1 μM</p> <p>LOD: n/a</p> <p>LOQ: n/a</p>	[118]

Valsartan	<p>Column: XTerra MS (50 mm × 2.1 mm; 3.5 μm), guard column: Opti guard-mini C18 (15 mm × 1 mm) Mobile phase: 0.1% trifluoroacetic acid:MeOH:⁶⁸N (45:30:25, v/v/v) Flow rate: 0.2 mL/min Column temperature: 50 °C Detection: SRM (positive ion mode): Valsartan: <i>m/z</i> 436.2 and 291.2 IS: <i>m/z</i> 429.2 and 401.2</p>	Human plasma	CGP48791	<p>Standard solutions: Stock solutions were prepared in MeOH; Calibration standards and QC samples were prepared by spiking into heparinized human plasma Samples: Mixture of 500 μL plasma + 50 μL IS solution (500 ng/mL in MeOH:H₂O (50:50, v/v) + 250 μL 2% trifluoroacetic acid) was placed in 96-well collection plates, then the entire volume was transferred into 96-well Empore SPE that pretreated with 100 μL MeOH and 500 μL 1% trifluoroacetic acid, rinsed with various solvent, then the analytes were eluted using 500 μL mixture of MeOH:H₂O (90:10, v/v) Eluate dried under N₂, residue reconstituted with 100 μL MeOH/0.1% trifluoroacetic acid (50:50, v/v), filtered then centrifuged for 1900 × <i>g</i> (10 min) at 4 °C. Filtrate applied to LC-MS/MS</p>	<p>Recovery: (63.0–73.0)% [37] Precision: Intraday: Precision (RSD): 2.19–5.40% Accuracy: 94.8–107% Interday: Precision (RSD): 1.87–5.67% Accuracy: 93.5–105% LOQ: 2 ng/mL LOQ: n/a</p>
-----------	--	--------------	----------	---	--

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
Valsartan, amlodipine, indapamide	<p>12 Column: Zorbax Eclipse Plus C18 (50 mm × 3 mm; 3 μm) Mobile phase: 62 Solvent A: ACN; B: H₂O with 0.1% formic acid (gradient elution): A:B (10:90) to (95:5) in 5 min, then held 1 min Detection: Qualitative: HRMS and NMR Quantitative: UV (valsartan: 250 nm; amlodipine: 360 nm; and indapamide: 240 nm)</p>	Capsules of Chinese herbal medicine	n/a	<p>Standard solutions: Analytes were dissolved in ACN:H₂O (1:1) Samples: Capsules were cut and emptied into 100 mL volumetric flask + 30 mL H₂O, stirred 30 min + 30 mL ACN stirred 30 min, and then diluted to volume with ACN:H₂O (1:1)</p>	<p>Recovery: 108.51 ± 3.40% Precision: n/a LOD: n/a LOQ: n/a</p>	[11]
Valsartan, amlodipine, hydrochlorothiazide	<p>Col 29: Aquasil C-18 (50 mm × 2.1 mm; 5 μm) 12 μm temperature: 20 °C Flow rate: 0.2 mL/min Mobile phase: ACN:H₂O containing 0.1% formic acid (50:50, v/v, for separating valsartan and amlodipine) ACN:H₂O containing 0.1% glacial acetic acid (60:40, v/v, for separating hydrochlorothiazide) Detection: MRM (positive ion mode): valsartan (<i>m/z</i> 436 → 418); losartan (<i>m/z</i> 423 → 406); amlodipine (<i>m/z</i> 409 → 238) Negative ion mode: hydrochlorothiazide (<i>m/z</i> 296 → 268); furosemide (<i>m/z</i> 328 → 285)</p>	Rat plasma	Losartan Furosemide (only for hydrochlorothiazide)	<p>Standard solutions: Stock solutions were prepared by dissolved analytes in ACN (10 mg/10 mL). Calibration and QC samples were prepared by spiking diluted stock solutions into drug-free rat plasma Samples: Analytes were extracted from plasma by adding ACN as precipitating agent, by vortexing (1 min), and centrifugating 10,000 rpm (10 min) Separated supernatant layer 67 20 μL was injected into LC-MS/MS</p>	<p>Recovery: Intraday: 87.9–106.3% Interday: 85.7–110% Precision: Intra- and interday: <14.3% LOQ: 1 ng/mL (all analytes)</p>	[119]

Valsartan, amlodipine, olmesartan medoxomil, hydrochlorothiazide	<p>Column: CN column (200 mm × 4.6 mm; 5 μm) Mobile phase: ACN: MeOH: 10 mM phosphoric acid pH 2.5 (7:13:80, v/v/v) Flow rate: 1.0 mL/min Detection: UV 235 nm</p>	Human plasma	n/a	<p>Standard solutions: Stock solutions (1 mg/mL) were prepared in MeOH; dilution was performed using a mixture of ACN:MeOH:H₂O (7:13:80, v/v/v) Samples:</p> <ul style="list-style-type: none"> - 1 mL plasma was spiked with various concentrations of working solutions and then basified with 0.5 mL NaOH - Analytes were extracted by vortex mixer (2 min) using mixture of 5 mL <i>n</i>-hexane: EtOA:isoamyl alcohol (88:10:2, v/v/v) - Samples were then centrifuged for 1 min (1500 rpm), and the organic layer was taken, evaporated, and then dissolved with 0.5 mL mobile phase 	<p>Recovery: 75.8–100.3% [89] Precision: <5.8% LOD: 0.1 ng/mL LOQ: 0.3 μg/mL</p>
Valsartan	<p>The system comprised of 3-μL level pumps (A, B, C), an in-line degasser, a sample cooler, a syringe pump, and a switch valve Separation was performed in a micro RP-C18 nano-flow column (150 μm i.d. and 375 μm o.d.; 3 μm).</p>	Human plasma (microliter sample)	Losartan	<p>Standard solutions: Stock solutions of valsartan were prepared in MeOH (1 mg/mL); IS solutions were prepared (10 μg/mL) in solutions of 100 g/L HFBA Samples: Human plasma samples (10 μL) were added in 500-μL vials and</p>	<p>Bias (relative error): [39] Intraday: -4.69% to +4.50% Interday: -2.75% to 3.00% Precision: Intraday: 4.47–6.70%; Interday: 4.71–6.43% LOD: n/a LOQ: 1 μg/mL</p>

Continued

Table 8 Summary of Analysis of Valsartan, Other Compounds, and Their Metabolite in the Biological Fluids and Samples—cont'd

Analyte(s)	Chromatography Conditions	Sample	Internal Standard (IS)	Preparation of Standards, Sample Extraction, and Cleanup	Accuracy, Precision (RSD), LOD, and LOQ Valsartan	Ref.
	Mobile phase: samples loading (pump 25 ACN:1% formic acid (20:80, v/v) Flow rate: 3 $\mu\text{L}/\text{min}$ For sample analysis (pump B: 5 pump C) 1% formic acid:ACN (20:80, v/v) Flow rate: 1 $\mu\text{L}/\text{min}$ Detection: SIR: valsartan m/z 436 ($[\text{M} + \text{H}]^+$) 27			then 5 μL solutions were added. Samples were vortexed (30 s) and centrifuged 45 000 rpm (2 min) 1 μL of supernatant was injected into the nanoscale LC-MS/MS		
Valsartan, losartan, irbesartan, candesartan cilexetil, candesartan cilexetil M1	Column: $\mu\text{Bondapak C18}$ (300 mm \times 3.9 mm; 10 μm) Guard column: Novapak C18 Mobile phase: ACN and sodium acetate 5 mM, pH 4 (gradient elution) Temperature: room Detector: UV 254 nm		Nimodipine	Standard solutions: Stock solutions (1000 $\mu\text{g}/\text{mL}$) prepared in MeOH and ACN Working solution were made daily by dilution with the same solvent Samples: <ol style="list-style-type: none"> 1. Frozen human urine ($-20\text{ }^\circ\text{C}$) thawed 2. 1 mL vortexed and mixed urine + PO_4 buffer (0.5 mL, 0.1 M, pH 2) \rightarrow shaken and centrifuged (3500 rpm) for 5 min 3. BondElut C8 cartridges + MeOH and 1 mL PO_4 buffer (0.1 M, pH 2) 	Recovery (RSD): 1.85–2.27% Precision (RSD): Intraday: 0.82–0.97% Interday: 2.63–4.30% LOD: n/a LOQ: 0.43 $\mu\text{g}/\text{mL}$	[120]

-
4. Cartridges do not dry before application of 1 mL acidified sample under low vacuum (<5 mmHg)
 5. Column conditioned with MeOH:PO₄ buffer pH 2 (40:60, v/v) → dried at vacuum (>200 mmHg) for 20 min
 6. Analytes eluted with 1 mL MeOH
 7. Eluate dry at 40 °C with N₂
 8. Remaining residue dissolved with 0.5 mL mobile phase
 9. 20 µL injected into HPLC system
-

ACKNOWLEDGMENTS

The authors would like to thank Mr. Widi Tri Raharjo, PT. Kromtekindo Utama, Jakarta, and Mr. Tegar Achsendo, Yuniarta Institute of Tropical Diseases, Airlangga University, Surabaya for measuring ESI-MS and NMR of valsartan, respectively.

REFERENCES

- [1] C.H. Czendlik, A. Sioufi, G. Preiswerk, H. Howald, Pharmacokinetic and pharmacodynamic interaction of single dose of valsartan and atenolol, *Eur. J. Clin. Pharmacol.* 52 (1997) 451–459.
- [2] N. Siddiqui, A. Husain, L. Chaudhry, M. Shamser Alam, M. Mitra, P.S. Bhasin, Pharmacological and pharmaceutical profile of valsartan: a review, *J. Appl. Pharm. Sci.* 01 (04) (2011) 12–19.
- [3] Valsartan in DailyMed, <http://dailymed.nlm.nih.gov/dailymed/search.cfm?startswith=VALSARTAN&cx=20&y=13>, 2014 (accessed 06.05.14).
- [4] S.C. Sweetman, *Martindale: The Complete Drug Reference*, 36th ed., The Pharmaceutical Press, London, 2009, pp. 1420–1421.
- [5] *British Pharmacopoeia 2013*, The Stationary Office, London, Data ©Crown Index + ©System Simulation Ltd Copyright 2012. <http://pharmatechbd.blogspot.com/2013/01/british-pharmacopoeia-2013.html> (accessed 07.05.14).
- [6] *The United States Pharmacopeial Convention, United States Pharmacopoeia—36 National Formulary 31*, The United States Pharmacopeial Convention, Rockville, MD, 2012, pp. 5537–5542.
- [7] Xpil, Patient Information Leaflets Online, Valsartan 40mg, 80mg, 160mg, 320mg Film-Coated Tablets. <http://xpil.medicines.org.uk/ViewPil.aspx?DocID=26854>, 2014 (accessed 03.02.14).
- [8] C.J. Mbah, Physicochemical properties of valsartan and the effect of ethyl alcohol, propylene glycol and pH on its solubility, *Pharmazie* 60 (2005) 849–850.
- [9] Merck & Co., Inc., *The Merck Index*, Merck & Co., Inc., Whitehouse Station, NJ, USA, 2008.
- [10] C. Potamitis, P. Chatzigeorgiou, E. Siapi, K. Viras, T. Mavromooustakos, A. Hodzic, G. Pabst, F. Cacho-Nerin, P. Laggner, M. Rappolt, Interactions of the AT1 antagonist with dipalmitoyl-phosphatidylcholine bilayers, *Biochim. Biophys. Acta* 1808 (2011) 1753–1763.
- [11] J.R. Kesting, J. Huang, D. Sorensen, Identification of adulterants in a Chinese herbal medicine by LC-HRMS and LC-MS-SPE/NMR and comparative in vivo study with standards in a hypertensive rat model, *J. Pharm. Biomed. Anal.* 51 (2010) 705–711.
- [12] B.N. Nalluri, R.M. Krishna, T.P. Rao, P.A. Crooks, Effect of recrystallization on the pharmaceutical properties of valsartan for improved therapeutic efficacy, *J. Appl. Pharm. Sci.* 2 (10) (2012) 126–132.
- [13] K.R. Rajeswari, K. Abbulu, M. Sudhakar, R. Karki, B. Rajkumar, Development and characterization of valsartan loaded hydrogel beads, *Pharm. Lett.* 4 (4) (2012) 1044–1053.
- [14] M.S. Sokar, A.S. Hanafy, A.H. El-Kamel, S.S. El-Gamal, Pulsatile core-in-cup valsartan tablet formulations: in vitro evaluation, *Asian J. Pharm. Sci.* 8 (2013) 234–243.
- [15] P. Raviteja, S. Muralidhar, R. Ramesh, T.V. Narayana, P.V. Kumar, G.V. Kumar, Formulation and evaluation of valsartan fast disintegrating tablets using solid dispersion technique, *Int. J. Innov. Pharm. Res.* 4 (1) (2013) 274–280.
- [16] A. Sharma, C.P. Jain, Preparation and characterization of solid dispersions of valsartan with poloxamer 188, *Pharm. Lett.* 2 (2) (2010) 54–63.

- [17] P.D. Maheskumar, G.S. Shantha Kumar, D. Goli, M.P. Talsania, H. Shah, Design, development and evaluation of nanoparticles of poorly soluble valsartan, *World J. Pharm. Pharm. Sci.* 2 (6) (2013) 5737–5749.
- [18] M. Skotnicki, A. Gawel, P. Cebe, M. Pyda, Thermal behaviour and phase identification of valsartan by standard and temperature-modulated differential scanning calorimetry, *Drug Dev. Ind. Pharm.* 39 (10) (2013) 1508–1514.
- [19] I. Rukhman, E. Flyaks, T. Koltai, J. Aronhime, Polymorphs of valsartan, United States Patent US 7,105,557, September 12, 2006.
- [20] I. Rukhman, E. Flyaks, T. Koltai, J. Aronhime, Amorphous form of valsartan, European Patent Specification EP 1950204 A1, July 30, 2008.
- [21] J. Burgbacher, B.T. Hahn, F.A. Rampf, R. Schneberger, Highly crystalline valsartan, United States Patent US 2013/0137737 A1, May 30, 2013.
- [22] Scientific Discussion EMEA. http://www.ema.europa.eu/docs/en_GB/document_library/EPAR_-_Scientific_Discussion/human/000716/WC500033006.pdf, 2007 (accessed 23.01.14).
- [23] J. Nie, B. Xiang, Y. Feng, D. Wang, Mass spectrometry and nuclear magnetic resonance spectroscopy, isolation and identification of process impurities in crude valsartan by HPLC, mass spectrometry, and nuclear magnetic spectroscopy, *J. Liq. Chromatogr. Relat. Technol.* 29 (2006) 553–568.
- [24] A. Sampath, A. Raghupathi Reddy, B. Yakambaran, A. Thirupathi, M. Prabhakar, P. Pratap Reddy, V. Prabhakar Reddy, Identification and characterization of potential impurities of valsartan AT1 receptor antagonist, *J. Pharm. Biomed. Anal.* 50 (2009) 405–412.
- [25] S. Mehta, R.P. Shah, S. Singh, Strategy for identification and characterization of small quantities of drug degradation products using LC and LC-MS: application to valsartan, a model drug, *Drug Test. Anal.* 2 (2010) 82–90.
- [26] R.M. Bianchini, P.M. Castellano, T.S. Kaufman, Characterization of two new potential impurities of valsartan obtained under photodegradation stress condition, *J. Pharm. Biomed. Anal.* 56 (2011) 16–22.
- [27] D. Ivanović, A. Malenović, B. Jančić, M. Medenica, M. Mašković, Monitoring of impurity level of valsartan and hydrochlorothiazide employing an RP-HPLC gradient mode, *J. Liq. Chromatogr. Relat. Technol.* 30 (2007) 2879–2890.
- [28] Valsartan Summary Validation Report. <https://mc.usp.org/sites/default/files/documents/ValidationReport/2012-06-11%20Valsartan%20Summary%20Validation%20Report.pdf>, 2012 (accessed 01.02.14).
- [29] Valsartan tablet summary validation report. <https://mc.usp.org/sites/default/files/documents/ValidationReport/2012-06-12%20Valsartan%20Tablets%20Summary%20Validation%20Report.pdf>, 2012 (accessed 01.02.14).
- [30] Hydrochlorothiazide and valsartan tablets summary validation report. <https://mc.usp.org/sites/default/files/documents/Hydrochlorothiazide%20and%20Valsartan%20Tablets%20Summary%20Validation%20Report%20-%202013-01-020.pdf>, 2012.
- [31] Amlodipine and valsartan tablets summary validation report. <https://mc.usp.org/sites/default/files/documents/Amlodipine%20and%20Valsartan%20Tablets%20Summary%20Validation%20Report%20-%202013-04-20.pdf>, 2013 (accessed 01.02.14).
- [32] L. González, J.A. Lopez, R.M. Alonso, R.M. Jiménez, Fast screening method for the determination of angiotensin II receptor antagonist in human plasma by high-performance liquid chromatography with fluorometric detection, *J. Chromatogr. A* 949 (2002) 49–60.
- [33] M. Bharathi, S.C.M. Prasad, R.L. Eswari, S. Wajim Raja, R.T. Allena, S. Brito Raj, K. Bhaskar Reddy, Preparation and in vitro & in vivo characterization of valsartan loaded eudragit nanoparticles, *Pharm. Sin.* 3 (5) (2012) 516–525.

- [34] A.K. Mahapatra, P.N. Murthy, S. Biswal, A.P.K. Mahapatra, S.P. Pradhan, Dissolution enhancement and physicochemical characterization of valsartan in solid dispersions with β -CD, HP β -CD, and PVP K-30, *Dissolution Technol.* 18 (2011) 39–45.
- [35] C.E.M. Jensen, R.A.S. Santos, A.M.L. Denadai, C.F.F. Santos, A.N.G. Braga, R.D. Sinisterra, Pharmaceutical composition of valsartan: β -cyclodextrin: physicochemical characterization and anti-hypertensive evaluation, *Molecules* 15 (2010) 4067–4084.
- [36] B. Chen, Y. Liang, Y. Wang, Development and validation of liquid chromatography-mass spectrometry method for the determination of telmisartan in human plasma, *Anal. Chim. Acta* 540 (2005) 367–373.
- [37] N. Koseki, H. Kawashita, H. Hara, M. Niina, M. Tanaka, R. Kawai, Y. Nagae, N. Masuda, Development and validation of a method for quantitative determination of valsartan in human plasma by liquid chromatography-tandem mass spectrometry, *J. Pharm. Biomed. Anal.* 43 (2007) 1769–1774.
- [38] X. Hu, Y. Zheng, J. Sun, L. Shang, G. Wang, H. Zhang, Simultaneous quantification of benazepril, glidazide and valsartan in human plasma by LC-MS-MS and application for rapidly measuring protein binding interaction between rhein and these three drugs, *Chromatographia* 69 (2009) 843–852.
- [39] C.Y. Lu, Y. Chang, W. Tseng, C. Feng, C. Lu, Analysis of angiotensin II receptor antagonist and protein markers at microliter level plasma by LC-MS/MS, *J. Pharm. Biomed. Anal.* 49 (2009) 123–128.
- [40] D. Zhang, X. Du, M. Liu, H. Li, Y. Jiang, L. Zhao, J. Gu, Determination of ecabet in human plasma by high performance liquid chromatography-tandem mass spectrometry, *J. Chromatogr. B* 863 (2008) 223–228.
- [41] H.J. Shah, N.B. Kataria, G. Subbaiah, C.N. Patel, Simultaneous LC-MS-MS analysis of valsartan and hydrochlorothiazide in human plasma, *Chromatographia* 69 (2009) 1055–1060.
- [42] H.H. Maurer, T. Kraemer, J.W. Arit, Screening for the detection of angiotensin-converting enzyme inhibitors, their metabolites, and AT II receptor antagonists, *Ther. Drug Monit.* 20 (6) (1998) 706–713.
- [43] F. Li, H. Zhang, L. Jiang, W. Zhang, J. Nie, Y. Feng, M. Yang, M. Liu, Dynamic NMR study and theoretical calculations on the conformational exchange of valsartan and related compounds, *Magn. Reson. Chem.* 45 (2007) 929–936.
- [44] C. Potamitis, M. Zervou, V. Katsiaras, P. Zoumpoulakis, S. Durdagi, M.G. Papadopoulos, J.M. Hayes, S.G. Grdadolnik, I. Kyrikou, D. Argyropoulos, G. Vatougia, T. Mavromoustakos, Antihypertensive drug valsartan in solution and the AT receptor: conformational analysis, dynamic NMR spectroscopy, in silico docking, and molecular dynamics simulations, *J. Chem. Inf. Model.* 49 (3) (2009) 726–739.
- [45] K.R. Gupta, A.R. Wadodkar, S.G. Wadodkar, UV-spectrophotometric methods for estimation of valsartan in bulk and tablet dosage form, *Int. J. ChemTech Res.* 2 (2) (2010) 985–989.
- [46] P.K. Kumar, M.U. Shankar, P. Subasini, P. Ghanshyam, P. Gourishyam, S. Kanhu, Stress degradation studies of valsartan and development of a validated method by UV spectrophotometric in bulk and pharmaceutical dosage form, *J. Pharm. Biomed. Sci.* 8 (2) (2011) 1–8.
- [47] E. Dinç, B. Uslu, S.A. Özkan, Spectral resolution of a binary mixture containing valsartan and hydrochlorothiazide in tablets by ratio spectra derivative and inverse least square techniques, *Anal. Lett.* 37 (4) (2004) 679–693.
- [48] A.B. Chaudhary, R.K. Patel, S.A. Chaudhary, K.V. Gadhvi, Estimation of valsartan and hydrochlorothiazide in pharmaceutical dosage forms by absorption ratio method, *Int. J. Appl. Biol. Pharm. Technol.* 1 (2) (2010) 455–464.

- [49] K.R. Gupta, A.D. Mahapatra, A.R. Wadodkar, S.G. Wadodkar, Simultaneous UV spectrophotometric determination of valsartan and amlodipine in tablet, *Int. J. ChemTech Res.* 2 (1) (2010) 551–556.
- [50] S. Ramachandran, B.K. Mandal, S.G. Navalgund, Simultaneous spectrophotometric determination of valsartan and ezetimibe in pharmaceuticals, *Trop. J. Pharm. Res.* 10 (6) (2011) 809–815.
- [51] A.S. Birajdar, S.N. Meyyanathan, B. Suresh, Simultaneous determination of nebivolol HCl and valsartan in solid dosage form by spectrophotometric and RP-HPLC method, *Int. J. Pharm. Sci. Res.* 2 (2) (2011) 424–431.
- [52] S. Hillaert, W. Van den Bossche, Optimization and validation of a capillary zone electrophoretic method for the analysis of several angiotensin-II-receptor antagonists, *J. Chromatogr. A* 979 (2002) 323–333.
- [53] S. Hillaert, W. Van den Bossche, Simultaneous determination of hydrochlorothiazide and several angiotensin-II-receptor antagonists by capillary electrophoresis, *J. Pharm. Biomed. Anal.* 31 (2003) 329–339.
- [54] S. Hillaert, T.R.M. De Beer, J.O. De Beer, W. Van den Bossche, Optimization and validation of a micellar electrokinetic chromatographic method for the analysis of several angiotensin-II-receptor antagonists, *J. Chromatogr. A* 984 (2003) 135–146.
- [55] G. Carlucci, V. Di Carlo, P. Mazzeo, Simultaneous determination of valsartan and hydrochlorothiazide in tablets by high performance liquid chromatography, *Anal. Lett.* 33 (12) (2000) 2491–2500.
- [56] E. Şatana, Ş. Altınay, N.G. Göğür, S.A. Özkan, Z. Şentürk, Simultaneous determination of valsartan and hydrochlorothiazide in tablets by first-derivative ultraviolet spectrophotometry and LC, *J. Pharm. Biomed. Anal.* 25 (2001) 1009–1013.
- [57] S. Tatar, S. Sağlık, Comparison of UV and second derivative-spectrophotometric and LC methods for the determination of valsartan in pharmaceutical formulation, *J. Pharm. Biomed. Anal.* 30 (2002) 371–375.
- [58] S.S. Chitlange, K. Bagri, D.M. Sakarkar, Stability indicating RP-HPLC method for simultaneous estimation of valsartan and amlodipine in capsule formulation, *Asian J. Res. Chem.* 1 (1) (2008) 15–18.
- [59] S.U. Kokil, M.S. Bhatia, Simultaneous estimation of nebivolol hydrochloride and valsartan using RP HPLC, *Indian J. Pharm. Sci.* 71 (2) (2009) 111–114.
- [60] D.F. Tian, X.L. Tian, T. Tian, Z.Y. Wang, F.K. Mo, Simultaneous determination of valsartan and hydrochlorothiazide in tablets by RP HPLC, *Indian J. Pharm. Sci.* 70 (3) (2008) 372–374.
- [61] V. Agrahari, V. Kabra, S. Gupta, R.K. Nema, M. Nagar, C. Karthikeyan, P. Trivedi, Determination of inherent stability of valsartan by stress degradation and its validation by HPLC, *Int. J. Pharm. Clin. Res.* 1 (2) (2009) 77–81.
- [62] S.B. Patel, B.G. Chaudhari, M.K. Buch, A.B. Patel, Stability indicating RP-HPLC method for simultaneous determination of valsartan and amlodipine from their combination product, *Int. J. ChemTech Res.* 1 (4) (2009) 1257–1267.
- [63] S.K. Patro, S.K. Kanungo, V.J. Patro, N.S.K. Choudhury, Stability indicating RP-HPLC method for determination of valsartan in pure and pharmaceutical formulation, *E-J. Chem.* 7 (1) (2010) 246–252.
- [64] B.M. Sudesh, K.S. Uttamrao, Determination and validation of valsartan and its degradation products by isocratic HPLC, *J. Chem. Metrol.* 3 (1) (2009) 1–12.
- [65] K.S. Lakshmi, L. Vasubramanian, A stability indicating HPLC method for the simultaneous determination of valsartan and ramipril in binary combination, *J. Chil. Chem. Soc.* 55 (2) (2010) 223–226.
- [66] N.H. Al-Shaalan, M.A. Alnowaiser, Simultaneous determination of amlodipine besilate and valsartan in pharmaceutical formulation using high performance liquid chromatography, *J. Chem. Pharm. Res.* 2 (6) (2010) 129–134.

- [67] G. Thanusha, C.J.G. Babu, K.P.C. Basavaraj, V.R. Panditi, C. Sharadha, Quantitative estimation of valsartan in bulk and pharmaceutical dosage forms, *Int. J. ChemTech Res.* 2 (2) (2010) 1194–1198.
- [68] D.U. Vinzuda, G.U. Sailor, N.R. Sheth, RP-HPLC method for determination of valsartan in tablet dosage forms, *Int. J. ChemTech Res.* 2 (3) (2010) 1461–1467.
- [69] D. Jothieswari, D. Priya, S. Brito Raj, E. Mohanambal, S. Wasim Raja, Design and RP-HPLC method for the simultaneous determination of valsartan and hydrochlorothiazide in bulk and in pharmaceutical formulation, *Int. J. Novel Trend. Pharm. Sci.* 1 (1) (2011) 18–22.
- [70] C.V.N. Prasad, C.S. Kumari, J. Sriramulu, A stability indicating RP-HPLC method for simultaneous estimation of valsartan and atorvastatin from their combination drug product, *Int. J. Pharm. Res. Anal.* 1 (1) (2011) 26–31.
- [71] V.B. Raju, A.L. Rao, Reversed phase HPLC analysis of valsartan in pharmaceutical dosage forms, *Int. J. Chem. Environ. Pharm. Res.* 2 (1) (2011) 56–60.
- [72] D.G.T. Parambi, M. Matthew, V. Ganesan, A validated stability indicating HPLC method for the determination of valsartan in tablet dosage forms, *J. Appl. Pharm. Sci.* 1 (4) (2011) 97–99.
- [73] A. Patnaik, M. Shetty, S. Sahoo, D.K. Nayak, S.K. Veliyath, A new RP-HPLC method for the determination of valsartan in bulk and its pharmaceutical formulations with its stability indicative studies, *Pharm. Sci. Monit.* 2 (3) (2011) 43–53.
- [74] K.S. Nataraj, S.S. Kumar, M.B. Duza, K. Reddy, Quantitative estimation of valsartan in pure and capsule dosage forms by RP-HPLC, *Int. J. Res. Pharm. Chem.* 1 (3) (2011) 587–590.
- [75] S. Rahila, P. Manisha, B. Deepali, M. Prabha, Y. Pramod, RP-HPLC method for simultaneous estimation of valsartan and hydrochlorothiazide in solid dosage form, *Int. Res. J. Pharm.* 2 (3) (2011) 162–164.
- [76] P.V. Santosh Kumar, M. Sahu, K.D. Prasad, M.C. Shekhar, Development and validation of analytical method for the estimation of valsartan in pure and tablet dosage form by RP-HPLC method, *Int. J. Res. Pharm. Chem.* 1 (4) (2011) 945–949.
- [77] P.V. Chokshi, K.J. Trivedi, N.S. Patel, Development and validation of RP-HPLC method for analysis of aliskiren hemifumarate and valsartan in their combination tablet dosage form, *Int. J. ChemTech Res.* 4 (4) (2012) 1623–1627.
- [78] N. Dubey, N. Dubey, M. Mandhanya, U. Bhadoriya, P. Sharma, Simultaneous estimation of ramipril and valsartan by RP-HPLC method in combined dosage form, *Asian J. Pharm. Med. Sci.* 2 (1) (2012) 23–28.
- [79] A.A. Ekshanawane, L.M. Abdelaziz, H.M. Hafez, Stability indicating HPLC method for simultaneous determination of several angiotensin-II-receptor antagonists in their dosage forms. *Pharm. Anal. Acta* 3 (8) (2012) 175. <http://dx.doi.org/10.4172/2153-2435.1000175>.
- [80] S.M. El-Gizawy, O.H. Abdelmageed, M.A. Omar, S.M. Deryea, A.M. Abdel-Megied, Development and validation of HPLC method for simultaneous determination of amlodipine, valsartan, hydrochlorothiazide in dosage form and spiked human plasma, *Am. J. Anal. Chem.* 3 (2012) 422–430.
- [81] M.A. Haque, S.H. Amrohi, P. Kumar, G. Nivedita, P. Kumar, D. Mohanty, P.V. Diwan, Stability indicating RP-HPLC method for the estimation of valsartan in pharmaceutical dosage form, *IOSR J. Pharm.* 2 (4) (2012) 12–18.
- [82] M.D. Kendre, S.K. Banerjee, Precise and accurate RP-HPLC method development for quantification of valsartan in tablet dosage form, *Int. J. Pharm. Sci. Drug Res.* 4 (2) (2012) 137–139.
- [83] M. Kharoaf, N. Malkieh, M. Abualhasan, R. Shubitah, N. Jaradat, A.N. Zaid, Tablet formulation and development of a validates stability indicating HPLC method for quantification of valsartan and hydrochlorothiazide combination, *Int. J. Pharm. Pharm. Sci.* 4 (3) (2012) 284–290.

- [84] S. Muralidharan, J.R. Kumar, S.A. Dhanaraj, Simple and validated RP-HPLC method for the estimation of valsartan in pharmaceutical tablet dosage form, *J. Chem. Pharm. Res.* 4 (9) (2012) 4235–4239.
- [85] R.R. Nahire, S.S. Joshi, V. Meghnani, N. Shastri, K.V.S. Nath, J. Satish, Stability indicating RP-HPLC method for simultaneous determination of amlodipine besylate and valsartan combination in bulk and commercial dosage forms, *Asian J. Pharm. Life Sci.* 2 (2) (2012) 281–290.
- [86] M.R. Rezk, N.M. El Remali, A.E.A. Aleem, Simultaneous determination of valsartan and hydrochlorothiazide in their pharmaceutical formulations, *Pharm. Chem.* 4 (1) (2012) 529–537.
- [87] K.N.K. Reddy, G.D. Rao, P.H. Pratyusha, Isocratic RP HPLC method validation of valsartan in pharmaceutical formulation with stress test stability evaluation of drug substance, *J. Chem. Pharm. Sci.* 5 (1) (2012) 8–12.
- [88] B. Brahmaiah, K. Sujana, A.P. Rani, Development and validation of RP-HPLC method for simultaneous determination of ramipril and valsartan in bulk and pharmaceutical dosage forms, *Asian J. Pharm. Clin. Res.* 6 (1) (2013) 23–25.
- [89] S.E. Kepekci Tekkeli, Development of an HPLC-UV method for the analysis of drugs used for combined hypertension therapy in pharmaceutical preparations and human plasma, *J. Anal. Method. Chem.* 2013 (2013) 1–9.
- [90] R.A. Shaalan, T.S. Belal, F.A. El Yazbi, S.M. Elonsy, Validated stability-indicating HPLC-DAD method of analysis for the antihypertensive triple mixture of amlodipine besylate, valsartan and hydrochlorothiazide. *Arab. J. Chem.* (2013). <http://dx.doi.org/10.1016/j.arabjc.2013.04.012>.
- [91] M. Sharma, K. Charmy, O. Sherikar, P. Mehta, Concurrent estimation of amlodipine besylate, hydrochlorothiazide and valsartan by RP-HPLC, HPTLC and UV-spectrophotometry, *J. Chromatogr. Sci.* 52 (2014) 27–35.
- [92] S.S. Imam, A. Ahad, M. Aqil, Y. Sultana, A. Ali, A validated RP-HPLC method for simultaneous determination of propranolol and valsartan in bulk drug and gel formulation, *J. Pharm. Bioallied Sci.* 5 (1) (2013) 61–65.
- [93] S.R. Dhaneswar, N.G. Patre, M.V. Mahadik, Validated TLC method for simultaneous quantitation of amlodipine besylate and valsartan in bulk drug and formulation, *Chromatographia* 69 (1/2) (2009) 157–161.
- [94] A.V. Gaikwad, V.G. Rajurkar, T. Shivakumar, G.Y. Dama, H.L. Tara, Simultaneous estimation of ramipril & valsartan in tablets by HPTLC, *Indo-Global J. Pharm. Sci.* 1 (1) (2011) 99–112.
- [95] T. Inglot, A. Gumieniczek, P. Maczka, Thin layer chromatography (TLC) and video-scanning in a new method for the analysis of angiotensin AT1 receptor antagonists, telmisartan and valsartan, in pharmaceuticals, *Ann. Univ. Mariae Curie-Skłodowska Lublin-Polonia XXIII* (4) (2010) 43–50.
- [96] S.P. Lokhande, S.P. Gupta, K. Sureshkumar, J. Dharumar, G. Garg, N. Upmanyu, Development and validation of a HPTLC method for the simultaneous estimation of ramipril and valsartan, *Int. J. Pharm. Teach. Pract.* 3 (1) (2012) 225–227.
- [97] D.D. Tsvetkova, D.P. Obreshkova, Application of validated TLC-densitometric method for simultaneous identification and determination of losartan potassium, telmisartan, and valsartan in tablets, *J. Planar Chromatogr.* 24 (4) (2012) 326–330.
- [98] A.R. Shrivastava, C.R. Barhate, C.J. Kapadia, Stress degradation studies on valsartan using validated stability-indicating high-performance thin-layer chromatography, *J. Planar Chromatogr.* 22 (6) (2009) 411–416.
- [99] N.J. Shah, B.N. Suhagia, R.R. Shah, N.M. Patel, HPTLC method for the simultaneous estimation of valsartan and hydrochlorothiazide in tablet dosage form, *Indian J. Pharm. Sci.* 71 (1) (2009) 72–74. <http://www.ijpsonline.com/article.asp?issn=0250-474X;year=2009;volume=71;issue=1;page=72;epage=74;aurlast=Shah> (accessed 30.03.14).

- [100] K. Czerwińska, A.P. Mazurek, Identification and determination of angiotensin II receptor antagonists with densitometric method, *Acta Pol. Pharm.* 63 (5) (2006) 402–404.
- [101] CAMAG Camag Bibliography Service CBS 51-112, Camag Muttenz (2014).
- [102] T.W. Inglot, K. Dabrowska, G. Misztal, The normal-phase retention behaviour of some angiotensin-II receptor antagonists, *J. Planar Chromatogr.* 20 (4) (2007) 293–301.
- [103] T.W. Inglot, K. Dabrowska, A. Gumieniczek, The reversed-phase retention behaviour of some angiotensin-II receptor antagonists, *J. Planar Chromatogr.* 22 (2) (2009) 145–155.
- [104] G. Iriarte, N. Ferreiros, I. Ibarrondo, R.M. Alonso, M.I. Maguregi, R.M. Jiménez, L. Gonzalez, Optimization via experimental design of an SPE-HPLC-UV-fluorescence method for the determination of valsartan and its metabolite in human plasma samples, *J. Sep. Sci.* 29 (2006) 2265–2283.
- [105] S. Kailasam, Simultaneous Determination of Therapeutic Drug Analytes in Human Plasma Using LC/MS/MS. Application Note, Agilent Technologies India, Bangalore, India, 2011.
- [106] H. Li, Y. Wang, Y. Jiang, Y. Tang, J. Wang, L. Zhao, J. Gu, A liquid chromatography/tandem mass spectrophotometry method for the simultaneous quantification of valsartan and hydrochlorothiazide in human plasma, *J. Chromatogr. B* 852 (2007) 436–442.
- [107] M. Levi, G. Wuerzner, E. Ezan, A. Pruvost, Direct analysis of valsartan and candesartan in human plasma and urines by on-line solid phase extraction coupled to electrospray tandem mass spectrometry, *J. Chromatogr. B* 877 (2009) 919–926.
- [108] N. Daneshtalab, R.Z. Lewanczuk, F. Jamali, High-performance liquid chromatographic analysis of angiotensin II receptor antagonist valsartan using a liquid extraction method, *J. Chromatogr. B* 766 (2002) 345–349.
- [109] M. Celebier, M.S. Kaynak, S. Altinöz, S. Şahin, Validated HPLC method development: the simultaneous analysis of amlodipine and valsartan in samples for liver perfusion studies, *Hacettepe Univ. J. Fac. Pharm.* 28 (1) (2008) 15–30.
- [110] J. Macek, J. Klíma, P. Ptáček, Rapid determination of valsartan in human plasma by protein precipitation and high-performance liquid chromatography, *J. Chromatogr. B* 832 (2006) 169–172.
- [111] A. Zarghi, A. Shafaati, S.M. Foroutan, H. Movahed, Rapid quantification of valsartan in human plasma by liquid chromatography using a monolithic column and a fluorescence detection: application in pharmacokinetics studies, *Sci. Pharm.* 76 (2008) 439–450.
- [112] M.R. Brunetto, Y. Contreras, S. Clavijo, D. Torres, Y. Delgado, F. Ovalles, C. Ayala, M. Gallignani, J.M. Estela, V.C. Martín, Determination of losartan, telmisartan, and valsartan by direct injection of human urine into a column-switching liquid chromatographic system with fluorescence detection, *J. Pharm. Biomed. Anal.* 50 (2009) 194–199.
- [113] A.C.F. Spínola, S. Almeida, A. Filipe, R. Neves, F. Trabelsi, A. Farré, Results of single-center, single-dose, randomized-sequence, open-label, two-way crossover bioequivalence study of two formulations of valsartan 160-mg tablets in healthy volunteers under fasting conditions, *Clin. Ther.* 31 (9) (2009) 1992–2001.
- [114] P.S. Selvan, K.V. Gowda, U. Mandal, W.D. Sam Solomon, T.K. Pal, Simultaneous determination of fixed dose combination of nebivolol and valsartan in human plasma by liquid chromatographic-tandem mass spectrometry and its application to pharmacokinetic study, *J. Chromatogr. B* 858 (2007) 143–150.
- [115] G. Iriarte, N. Ferreiros, I. Ibarrondo, R.M. Alonso, M.I. Maguregi, R.M. Jiménez, Biovalidation of an SPE-HPLC-UV-fluorescence method for the determination of valsartan and its metabolite valeryl-4-hydroxy-valsartan in human plasma, *J. Sep. Sci.* 30 (2007) 2231–2240.

- [116] R.N. Sharma, S.S. Pancholi, Simple RP-HPLC method for the determination of triple drug combination of valsartan, amlodipine and hydrochlorothiazide in human plasma, *Acta Pharm.* 62 (2012) 45–58.
- [117] A.R. Dixit, S.R. Rajput, S.G. Patel, Preparation and bioavailability assessment of SMEDDS containing valsartan, *AAPS PharmSciTech* 11 (1) (2010) 314–321.
- [118] L. Kristoffersen, E.L. Øistad, M.S. Opdal, M. Krogh, E. Lundanes, A.S. Christophersen, Simultaneous determination of 6 beta-blockers, 3 calcium-channel antagonists, 4 angiotensin-II antagonists and 1 antiarrhythmic drug in post-mortem whole blood by automated solid phase extraction and liquid chromatography mass spectrometry method development and robustness testing by experimental design, *J. Chromatogr. B* 850 (2007) 147–160.
- [119] G.S. Ganesh, P. Deme, K. Madhusudana, R. Sista, Simultaneous determination of amlodipine, valsartan and hydrochlorothiazide by LC-ESI-MS/MS and its application to pharmacokinetics in rats, *J. Pharm. Anal.* 4 (2014) 399–406. <http://dx.doi.org/10.1016/j.jpha.2013.12.003>.
- [120] L. González, R.M. Alonso, R.M. Jimenez, A high-performance liquid chromatographic method for screening angiotensin II receptor antagonists in human urine, *Chromatographia* 52 (2000) 735–740.

Valsartan

ORIGINALITY REPORT

20%

SIMILARITY INDEX

19%

INTERNET SOURCES

11%

PUBLICATIONS

0%

STUDENT PAPERS

PRIMARY SOURCES

1	www.tandfonline.com Internet Source	2%
2	www.hplc-eshop.ch Internet Source	1%
3	www.yumpu.com Internet Source	1%
4	gnu.inflibnet.ac.in Internet Source	1%
5	www.scribd.com Internet Source	1%
6	www.science.gov Internet Source	1%

7	documents.mx Internet Source	1%
8	www.jove.com Internet Source	1%
9	www.hplc.sk Internet Source	1%
10	www.omicsonline.org Internet Source	<1%
11	www.agilent.com Internet Source	<1%
12	cn.agilent.com Internet Source	<1%
13	www.rroj.com Internet Source	<1%
14	www.ukessays.com Internet Source	<1%
15	www.usp.org Internet Source	<1%

16	pandelisperakakis.info Internet Source	<1 %
17	Muszalska, Izabela, Agnieszka Sobczak, Agnieszka Dołhań, and Anna Jelińska. "Analysis of Sartans: A Review", Journal of Pharmaceutical Sciences, 2014. Publication	<1 %
18	baadalsg.inflibnet.ac.in Internet Source	<1 %
19	link.springer.com Internet Source	<1 %
20	Miroslava Lachová, Jozef Lehotay, Gabriela Karasová, Ivan Skačáni, D. W. Armstrong. " Isolation of -Theanine from Plant Material using a Molecularly Imprinted Polymer ", Journal of Liquid Chromatography & Related Technologies, 2007 Publication	<1 %
21	www.eurasianjournals.com Internet Source	<1 %
22	www.jmess.org Internet Source	<1 %
23	www.guidechem.com	

Internet Source

<1 %

24

www.chemspider.com

Internet Source

<1 %

25

www.mdpi.com

Internet Source

<1 %

26

shodhganga.inflibnet.ac.in

Internet Source

<1 %

27

Benzi, Micaela, Marco Bobba, Marco Demartini, Valentina Gianotti, Fabio Gosetti, Emilio Marengo, Eleonora Mazzucco, Elisa Robotti, and Davide Zampieri. "Organic Bases", Food Analysis by HPLC Third Edition, 2012.

Publication

<1 %

28

sci-pharm.at

Internet Source

<1 %

29

V. F. Samanidou, S. A. Nisyriou, I. N. Papadoyannis. "Residue Analysis of Penicillins in Food Products of Animal Origin by HPLC: A Review", Journal of Liquid Chromatography & Related Technologies, 2007

Publication

<1 %

- 30 Izabela Muszalska, Agnieszka Sobczak, Agnieszka Dołhań, Anna Jelińska. "Analysis of Sartans: A Review", Journal of Pharmaceutical Sciences, 2014 <1 %
Publication
-
- 31 teses.usp.br <1 %
Internet Source
-
- 32 Mikaël Levi, Grégoire Wuerzner, Eric Ezan, Alain Pruvost. "Direct analysis of valsartan or candesartan in human plasma and urines by on-line solid phase extraction coupled to electrospray tandem mass spectrometry", Journal of Chromatography B, 2009 <1 %
Publication
-
- 33 www.scilit.net <1 %
Internet Source
-
- 34 docplayer.net <1 %
Internet Source
-
- 35 K. S. Kokilambigai, R. Seetharaman, K. S. Lakshmi. "Critical Review on the Analytical Techniques for the Determination of the Oldest Statin—Atorvastatin—in Bulk, Pharmaceutical Formulations and Biological Fluids", Critical Reviews in Analytical Chemistry, 2017 <1 %
Publication

36

Internet Source

<1%

37

impactfactor.org

Internet Source

<1%

38

ijpba.info

Internet Source

<1%

39

knauer.net

Internet Source

<1%

40

www.nicpbp.org.cn

Internet Source

<1%

41

Romina M. Bianchini, Patricia M. Castellano, Teodoro S. Kaufman. "Characterization of two new potential impurities of Valsartan obtained under photodegradation stress condition", Journal of Pharmaceutical and Biomedical Analysis, 2011

Publication

<1%

42

www.asianpharmtech.com

Internet Source

<1%

43

gmp.kiev.ua

Internet Source

<1%

44	www.plantphysiol.org Internet Source	<1 %
45	Chi-Yu Lu, Yung-Ming Chang, Wei-Lung Tseng, Chia-Hsien Feng, Chin-Yuan Lu. "Analysis of angiotensin II receptor antagonist and protein markers at microliter level plasma by LC–MS/MS", Journal of Pharmaceutical and Biomedical Analysis, 2009 Publication	<1 %
46	www.antec.it Internet Source	<1 %
47	www.chiraltech.com Internet Source	<1 %
48	derpharmachemica.com Internet Source	<1 %
49	www.bcmst.ru Internet Source	<1 %
50	lib.bioinfo.pl Internet Source	<1 %
51	www.e-journals.in Internet Source	<1 %

52	www.portanna.org Internet Source	<1%
53	www.hplc.eu Internet Source	<1%
54	www.mysciencework.com Internet Source	<1%
55	www.drogallega.es Internet Source	<1%
56	www.crawfordscientific.com Internet Source	<1%
57	www.ijpsonline.com Internet Source	<1%
58	academicjournals.org Internet Source	<1%
59	bmccomplementaltemed.biomedcentral.com Internet Source	<1%
60	www.google.com.tr Internet Source	<1%

61	www.niper.gov.in Internet Source	<1 %
62	vercopak.com.tw Internet Source	<1 %
63	doras.dcu.ie Internet Source	<1 %
64	www.esrnexus.com Internet Source	<1 %
65	Niki Vergi-Athanasiou, Julia Atta-Politou, Michael Koupparis, John Spyropoulos. "Development and Validation of an HPLC Method, with Fluorescence Detection, for Simultaneous Determination of Paroxetine and its Metabolites in Plasma", Journal of Liquid Chromatography & Related Technologies, 2007 Publication	<1 %
66	www.mz-at.de Internet Source	<1 %
67	Liberato Brum Junior, Danieli Cátia Ceni, Marcio Fronza, Paulo Renato de Oliveira, Sérgio Luiz Dalmora. "Validation of an LC-Tandem MS/MS Method for the Determination of Etoricoxib in Human Plasma and Pharmaceutical Formulations", Journal of Liquid	<1 %

Chromatography & Related Technologies, 2007

Publication

68	www.teknokroma.es Internet Source	<1 %
69	www.ijrps.pharmascope.org Internet Source	<1 %
70	addi.ehu.es Internet Source	<1 %
71	recipe.de Internet Source	<1 %
72	nrl.northumbria.ac.uk Internet Source	<1 %
73	C. Kowalski, M. Pomorska, T. Sławik. "Development of HPLC with UV-VIS Detection for the Determination of the Level of Oxytetracycline in the Biological Matrix", Journal of Liquid Chromatography & Related Technologies, 2006 Publication	<1 %
74	es.scribd.com Internet Source	<1 %

75	Peyami Sari, Jianguo Sun, Majid Razzak, Ian G. Tucker. "HPLC Assay of Levamisole and Abamectin in Sheep Plasma for Application to Pharmacokinetic Studies", Journal of Liquid Chromatography & Related Technologies, 2007 Publication	<1%
76	Su-Geun Yang, So-Ra Park, Dae-Duk Kim, Suk-Jae Chung, Chang-Koo Shim. "A Simple HPLC Method for the Determination of Cyclosporin A in Human Whole Blood", Journal of Liquid Chromatography & Related Technologies, 2007 Publication	<1%
77	sure.sunderland.ac.uk Internet Source	<1%
78	kprime.net Internet Source	<1%
79	eng.hi138.com Internet Source	<1%
80	Ponnayyan Sulochana, Suresh, Kuldeep Sharma, Ramesh Mullangi, and Sathesh Kumar Sukumaran. "Review of the validated HPLC and LC-MS/MS methods for determination of drugs used in clinical practice for Alzheimer's disease : Bioanalytical methods for quantitation of drugs to treat Alzheimer's", Biomedical Chromatography, 2014. Publication	<1%

- 81 Selvan, P.S.. "Simultaneous determination of fixed dose combination of nebivolol and valsartan in human plasma by liquid chromatographic-tandem mass spectrometry and its application to pharmacokinetic study", Journal of Chromatography B, 20071015 <1 %
Publication
-
- 82 Agbaba, Danica, and Katarina Nikolić. "TLC of Antihypertensive and Antihypotensive Drugs", Chromatographic Science Series, 2013. <1 %
Publication
-
- 83 Darko Ivanović, Anđelija Malenović, Biljana Jančić, Mirjana Medenica, Marija Mašković. "Monitoring of Impurity Level of Valsartan and Hydrochlorothiazide Employing an RP–HPLC Gradient Mode", Journal of Liquid Chromatography & Related Technologies, 2007 <1 %
Publication
-
- 84 Cafer Saka, Ömer Şahin. "Determination of Serotonin-Norepinephrine Reuptake Inhibitor Antidepressants in Pharmaceuticals and Biological Material", Critical Reviews in Analytical Chemistry, 2013 <1 %
Publication
-
- 85 Irena Baranowska, Sylwia Magiera, Jacek Baranowski. "Clinical applications of fast liquid chromatography: A review on the analysis of cardiovascular drugs and their metabolites", Journal of Chromatography B, 2013 <1 %
Publication

Sampath, A.. "Identification and characterization of potential impurities of valsartan, AT1

Exclude quotes

Off

Exclude matches

Off

Exclude bibliography

On